
Cleo

Kyle Johnsen, Nathan Cruzado

Oct 13, 2023

CONTENTS

1	Closed Loop processing	3
2	Electrode recording	5
3	Optogenetic stimulation	7
4	Getting started	9
5	Related resources	11
5.1	Publications	11
6	Documentation contents	13
6.1	Overview	13
6.2	Tutorials	18
6.3	Reference	74
7	Indices and tables	117
	Python Module Index	119
	Index	121

Hello there! Cleo has the goal of bridging theory and experiment for mesoscale neuroscience, facilitating electrode recording, optogenetic stimulation, and closed-loop experiments (e.g., real-time input and output processing) with the [Brian 2](#) spiking neural network simulator. We hope users will find these components useful for prototyping experiments, innovating methods, and testing observations about a hypotheses *in silico*, incorporating into spiking neural network models laboratory techniques ranging from passive observation to complex model-based feedback control. Cleo also serves as an extensible, modular base for developing additional recording and stimulation modules for Brian simulations.

This package was developed by [Kyle Johnsen](#) and Nathan Cruzado under the direction of [Chris Rozell](#) at Georgia Institute of Technology.

CLOSED LOOP PROCESSING

Cleo allows for flexible I/O processing in real time, enabling the simulation of closed-loop experiments such as event-triggered or feedback control. The user can also add latency to closed-loop stimulation to study the effects of computation delays.

ELECTRODE RECORDING

Cleo provides functions for configuring electrode arrays and placing them in arbitrary locations in the simulation. The user can then specify parameters for probabilistic spike detection or a spike-based LFP approximation developed by [Teleńczuk et al., 2020](#).

OPTOGENETIC STIMULATION

By providing an optic fiber-light propagation model, Cleo enables users to flexibly add photostimulation to their model. Both a four-state Markov state model of opsin dynamics is available, as well as a minimal proportional current option for compatibility with simple neuron models. Parameters are provided for the common blue light/ChR2 setup.

GETTING STARTED

Just use pip to install—the name on PyPI is `cleosim`:

```
pip install cleosim
```

Then head to the [overview section of the documentation](#) for a more detailed discussion of motivation, structure, and basic usage.

RELATED RESOURCES

Those using Cleo to simulate closed-loop control experiments may be interested in software developed for the execution of real-time, *in-vivo* experiments. Developed by members of [Chris Rozell](#)'s and [Garrett Stanley](#)'s labs at Georgia Tech, the [CLOCTools repository](#) can serve these users in two ways:

1. By providing utilities and interfaces with experimental platforms for moving from simulation to reality.
2. By providing performant control and estimation algorithms for feedback control. Although Cleo enables closed-loop manipulation of network simulations, it does not include any advanced control algorithms itself. The `ldsCtrlEst` library implements adaptive linear dynamical system-based control while the `hmm` library can generate and decode systems with discrete latent states and observations.

5.1 Publications

CLOC Tools: A Library of Tools for Closed-Loop Neuroscience A.A. Willats, M.F. Bolus, K.A. Johnsen, G.B. Stanley, and C.J. Rozell. *In prep*, 2022.

State-Aware Control of Switching Neural Dynamics A.A. Willats, M.F. Bolus, C.J. Whitmire, G.B. Stanley, and C.J. Rozell. *In prep*, 2022.

Closed-Loop Identifiability in Neural Circuits A. Willats, M. O'Shaughnessy, and C. Rozell. *In prep*, 2022.

State-space optimal feedback control of optogenetically driven neural activity M.F. Bolus, A.A. Willats, C.J. Rozell and G.B. Stanley. *Journal of Neural Engineering*, 18(3), pp. 036006, March 2021.

Design strategies for dynamic closed-loop optogenetic neurocontrol in vivo M.F. Bolus, A.A. Willats, C.J. Whitmire, C.J. Rozell and G.B. Stanley. *Journal of Neural Engineering*, 15(2), pp. 026011, January 2018.

DOCUMENTATION CONTENTS

6.1 Overview

6.1.1 Introduction

Who is this package for?

Cleo (Closed Loop, Electrophysiology, and Optogenetics Simulator) is a Python package developed to bridge theory and experiment for mesoscale neuroscience. We envision two primary uses cases:

1. For prototyping closed-loop control of neural activity *in silico*. Animal experiments are costly to set up and debug, especially with the added complexity of real-time intervention—our aim is to enable researchers, given a decent spiking model of the system of interest, to assess whether the type of control they desire is feasible and/or what configuration(s) would be most conducive to their goals.
2. The complexity of experimental interfaces means it's not always clear what a model would look like in a real experiment. Cleo can help anyone interested in observing or manipulating a model while taking into account the constraints present in real experiments. Because Cleo is built around the [Brian simulator](#), we especially hope this is helpful for existing Brian users who for whatever reason would like a convenient way to inject recorders (e.g., electrodes) or stimulators (e.g., optogenetics) into the core network simulation.

What is closed-loop control?

In short, determining the inputs to deliver to a system from its outputs. In neuroscience terms, making the stimulation parameters a function of the data recorded in real time.

Structure and design

Cleo wraps a spiking network simulator and allows for the injection of stimulators and/or recorders. The models used to emulate these devices are often non-trivial to implement or use in a flexible manner, so Cleo aims to make device injection and configuration as painless as possible, requiring minimal modification to the original network.

Cleo also orchestrates communication between the simulator and a user-configured [IOProcessor](#) object, modeling how experiment hardware takes samples, processes signals, and controls stimulation devices in real time.

For an explanation of why we choose to prioritize spiking network models and how we chose Brian as the underlying simulator, see [Design rationale](#).

Why closed-loop control in neuroscience?

Fast, real-time, closed-loop control of neural activity enables intervention in processes that are too fast or unpredictable to control manually or with pre-defined stimulation, such as sensory information processing, motor planning, and oscillatory activity. Closed-loop control in a *reactive* sense enables the experimenter to respond to discrete events of interest, such as the arrival of a traveling wave or sharp wave ripple, whereas *feedback* control deals with driving the system towards a desired point or along a desired state trajectory. The latter has the effect of rejecting noise and disturbances, reducing variability across time and across trials, allowing the researcher to perform inference with less data and on a finer scale. Additionally, closed-loop control can compensate for model mismatch, allowing it to reach more complex targets where open-loop control based on imperfect models is bound to fail.

6.1.2 Installation

Make sure you have Python ≥ 3.7 , then use pip: `pip install cleosim`.

Note: The name on PyPI is `cleosim` since `cleo` was already taken, but in code it is still used as `import cleo`. The other Cleo appears to actually be a fairly well developed package, so I'm sorry if you need to use it along with this Cleo in the same environment. In that case, [there are workarounds](#).

Or, if you're a developer, [install poetry](#) and run `poetry install` from the repository root.

6.1.3 Usage

Brian network model

The starting point for using Cleo is a Brian spiking neural network model of the system of interest. For those new to Brian, the [docs](#) are a great resource. If you have a model built with another simulator or modeling language, you may be able to [import it to Brian via NeuroML](#).

Perhaps the biggest change you may have to make to an existing model to make it compatible with Cleo's optogenetics and electrode recording is to give the neurons of interest coordinates in space. See the [Tutorials](#) or the [cleo.coords](#) module for more info.

You'll need your model in a Brian [Network](#) object before you move on. E.g.,:

```
net = brian2.Network(...)
```

CLSimulator

Once you have a network model, you can construct a [CLSimulator](#) object:

```
sim = cleo.CLSimulator(net)
```

The simulator object wraps the Brian network and coordinates device injection, processing input and output, and running the simulation.

Recording

Recording devices take measurements of the Brian network. Some extremely simple implementations (which do little more than wrap Brian monitors) are available in the `cleo.recorders` module.

To use a *Recorder*, you must inject it into the simulator via `inject()`:

```
rec = MyRecorder('recorder_name', ...) # note that all devices need a unique name
sim.inject(rec, neuron_group1, neuron_group2, ...) # can pass in additional arguments
```

The recorder will only record from the neuron groups specified on injection, allowing for such scenarios as singling out a cell type to record from.

Electrodes

Electrode recording is the main recording modality currently implemented in Cleo. See the *Electrode recording* tutorial for more detail, but in brief, usages consists of:

1. Constructing a *Probe* object with coordinates at the desired contact locations
 - Convenience functions for generating shank probe coordinates exist. See *Specifying electrode coordinates*.
2. Specifying the signals to be recorded. Currently there are three implemented. See *Specifying signals to record*.
 - Multi-unit activity
 - Sorted spikes
 - TKLFP: Teleńczuk kernel approximation of LFP
3. Injection into the simulator

Stimulation

Stimulator devices manipulate the Brian network. Usage is similar to recorders:

```
stim = MyStimulator('stimulator_name', ...) # again, all devices need a unique name
# again, specify neuron groups device will affect and any additional arguments needed
sim.inject(stim, neuron_group1, neuron_group2, ...)
```

As with recorders, you can inject stimulators per neuron group to produce a targeted effect.

Optogenetics

Optogenetics is the main stimulator device currently implemented by Cleo. This takes the form of an *OptogeneticIntervention*, which, on injection, adds a light source at the specified location and transfects the neurons (via Brian “synapses” that deliver current according to an opsin model, leaving the neuron model equations untouched).

Out of the box you can access a four-state Markov model of channelrhodopsin-2 (ChR2) and parameters for a 473-nm blue optic fiber light source.:

```
from cleo.opto import *
opto = OptogeneticIntervention(
    name="...",
```

(continues on next page)

(continued from previous page)

```
opsin_model=FourStateModel(params=chr2_4s),
light_model_params=default_blue,
location=(0, 0, 0.5) * mm,
)
```

Note, however, that Markov opsin dynamics models require target neurons to have membrane potentials in realistic ranges and an *Iopto* term defined in amperes. If you need to interface with a model without these features, you may want to use the simplified `ProportionalCurrentModel`. You can find more details, including a comparison between the two model types, in the [optogenetics tutorial](#).

These model and parameter settings were designed to be flexible enough that an interested user should be able to imitate and replace them with other opsins, light sources, etc. See the [Optogenetic stimulation](#) tutorial for more detail.

IO Processor

Just as in a real experiment where the experiment hardware must be connected to signal processing equipment and/or computers for recording and control, the `CLSimulator` must be connected to an `IOProcessor`:

```
sim.set_io_processor(...)
```

If you are only recording, you may want to use the `RecordOnlyProcessor`. Otherwise you will want to implement the `LatencyIOProcessor`, which not only takes samples at the specified rate, but processes the data and delivers input to the network after a user-defined delay, emulating the latency inherent in real experiments. You define your processor by creating a subclass and defining the `process()` function:

```
class MyProcessor(LatencyIOProcessor):

    def process(self, state_dict, sample_time_ms):
        # state_dict contains a {'recorder_name': value} dict of network
        foo = state_dict['foo_recorder']
        out = ... # do something with sampled spikes
        delay_ms = 3
        t_out_ms = sample_time_ms + delay_ms
        # output must be a {'stimulator_name': value} dict setting stimulator values
        return {'stim': out}, t_out_ms

my_proc = MyProcessor(sample_period_ms=1)
sim.set_io_processor(my_proc)
```

See [On-off control](#) for a minimal working example or [PI control](#) for more advanced features, including decomposing the processing into blocks with accompanying stochastic delay objects.

Running experiments

Use `CLSimulator`'s `run()` function with the desired duration:

```
sim.run(500*ms, ...) # kwargs are passed to Brian's run function
```

Use `CLSimulator`'s `reset()` function to restore the default state (right after initialization/injection) for the network and all devices. This could be useful for running a simulation multiple times under different conditions.

To facilitate access to data after the simulation, many classes offer a `save_history` option on construction. If true, that object will store relevant variables as attributes. For example,:

```
sorted_spikes = cleo.ephys.SortedSpiking(...)
...
sim.run(...)

plt.plot(sorted_spikes.t_ms, sorted_spikes.i)
```

6.1.4 Design rationale

Why not prototype with more abstract models?

Cleo aims to be practical, and as such provides models at the level of abstraction corresponding to the variables the experimenter has available to manipulate. This means models of spatially defined, spiking neural networks.

Of course, neuroscience is studied at many spatial and temporal scales. While other projects may be better suited for larger segments of the brain and/or longer timescales (such as [HNN](#) or BMTK’s [PopNet](#) or [FilterNet](#)), this project caters to finer-grained models because they can directly simulate the effects of alternate experimental configurations. For example, how would the model change when swapping one opsin for another, using multiple opsins simultaneously, or with heterogeneous expression? How does recording or stimulating one cell type vs. another affect the experiment? Would using a more sophisticated control algorithm be worth the extra compute time, and thus later stimulus delivery, compared to a simpler controller?

Questions like these could be answered using an abstract dynamical system model of a neural circuit, but they would require the extra step of mapping the afore-mentioned details to a suitable abstraction—e.g., estimating a transfer function to model optogenetic stimulation for a given opsin and light configuration. Thus, we haven’t emphasized these sorts of models so far in our development of Cleo, though they should be possible to implement in Brian if you are interested. For example, one could develop a Poisson linear dynamical system (PLDS), record spiking output, and configure stimulation to act directly on the system’s latent state.

And just as experiment prototyping could be done on a more abstract level, it could also be done on an even more realistic level, which we did not deem necessary. That brings us to the next point...

Why Brian?

Brian is a relatively new spiking neural network simulator written in Python. Here are some of its advantages:

- Flexibility: allowing (and requiring!) the user to define models mathematically rather than selecting from a pre-defined library of cell types and features. This enables us to define arbitrary models for recorders and stimulators and easily interface with the simulation
- Ease of use: it’s all just Python
- Speed

[NEST](#) is a popular alternative to Brian also strong in point neuron simulations. However, it appears to be less flexible, and thus harder to extend. [NEURON](#) is another popular alternative to Brian. Its main advantage is its first-class support of detailed, morphological, multi-compartment neurons. In fact, strong alternatives to Brian for this project were BioNet ([docs](#), [paper](#)) and NetPyNE ([docs](#), [paper](#)), which already offer a high-level interface to NEURON with extracellular potential recording. Optogenetics could be incorporated with [pre-existing .hoc code](#), though the light model would need to be implemented. From brief examination of the [source code of BioNet](#), it appears that closed-loop stimulation would not be too difficult to add. It is unclear for NetPyNE.

In the end, we chose Brian since our priority was to model circuit/population-level dynamics over molecular/intra-neuron dynamics. Also, Brian does have support for multi-compartment neurons, albeit less fully featured, if that is needed.

6.1.5 Future development

Here are some features which are missing but could be useful to add:

- Better support for multiple opsins simultaneously. At present the user would have to include a separate variable for each new opsin current, which makes changing the number of different opsins inconvenient
- Support for multiple light sources affecting a single opsin transfection—whether the light sources have the same or different wavelengths
- Electrode microstimulation
- A more accurate LFP signal (only usable for morphological neurons) based on the volume conductor forward model as in [LFPy](#) or [Vertex](#)
- The [Mazzoni-Lindén LFP approximation](#) for LIF point-neuron networks
- Imaging as a recording modality

6.2 Tutorials

6.2.1 Electrode recording

How to insert electrodes to measure different spiking and extracellular signals from a Brian network simulation.

Preamble:

```
from brian2 import * # includes numpy
import cleo
from cleo import *
# the default cython compilation target isn't worth it for
# this trivial example
prefs.codegen.target = "numpy"
np.random.seed(1919)

cleo.utilities.style_plots_for_docs()

# colors
c = {
    'light': '#df87e1',
    'main': '#C500CC',
    'dark': '#8000B4',
    'exc': '#d6755e',
    'inh': '#056eee',
    'accent': '#36827F',
}
```

Network setup

First we create a toy E-I network with Poisson firing rates and assign coordinates:


```
N = 500
n_e = int(N * 0.8)
n_i = int(N * 0.2)

exc = PoissonGroup(n_e, 10 * Hz, name="exc")
inh = PoissonGroup(n_i, 30 * Hz, name="inh")

net = Network([exc, inh])
sim = CLSimulator(net)

cleo.coords.assign_coords_rand_rect_prism(
    exc, xlim=(-0.2, 0.2), ylim=(-0.2, 0.2), zlim=(0.7, 0.9)
)
cleo.coords.assign_coords_rand_rect_prism(
    inh, xlim=(-0.2, 0.2), ylim=(-0.2, 0.2), zlim=(0.7, 0.9)
)
cleo.viz.plot(exc, inh, colors=[c['exc'], c['inh']], scatterargs={'alpha': .6})
```

```
(<Figure size 640x480 with 1 Axes>,
 <Axes3D: xlabel='x (mm)', ylabel='y (mm)', zlabel='z (mm)'>)
```

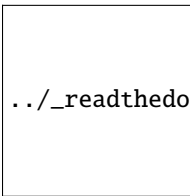

../_readthedocs/jupyter_execute/electrodes_3_1.png

Specifying electrode coordinates

Now we insert an electrode shank probe in the center of the population by injecting an Probe device. Note that Probe takes arbitrary coordinates as arguments, so you can place contacts wherever you wish. However, the `cleo.ephys` module provides convenience functions to easily generate coordinates [common in NeuroNexus probes](#). Here are some examples:

```
from cleo import ephys
from mpl_toolkits.mplot3d import Axes3D

array_length = 0.4 * mm # length of the array itself, not the shank
tetr_coords = ephys.tetrode_shank_coords(array_length, tetrode_count=3)
poly2_coords = ephys.poly2_shank_coords(
    array_length, channel_count=32, intercol_space=50 * umeter
)
poly3_coords = ephys.poly3_shank_coords(
    array_length, channel_count=32, intercol_space=30 * umeter
)
# by default start_location (location of first contact) is at (0, 0, 0)
single_shank = ephys.linear_shank_coords(
```

(continues on next page)


(continued from previous page)

```

    array_length, channel_count=8, start_location=(-0.2, 0, 0) * mm
)
# tile vector determines length and direction of tiling (repeating)
multishank = ephys.tile_coords(single_shank, num_tiles=3, tile_vector=(0.4, 0, 0) * mm)

fig = plt.figure(figsize=(8, 8))
fig.suptitle("Example array configurations")
for i, (coords, title) in enumerate(
    [
        (tetr_coords, "3-tetrode shank"),
        (poly2_coords, "32-channel Poly2 shank"),
        (poly3_coords, "32-channel Poly3 shank"),
        (multishank, "Multi-shank"),
    ],
    start=1,
):
    ax = fig.add_subplot(2, 2, i, projection="3d")
    x, y, z = coords.T / umeter
    ax.scatter(x, y, z, marker="x", c="black")
    ax.set(
        title=title,
        xlabel="x (m)",
        ylabel="y (m)",
        zlabel="z (m)",
        xlim=(-200, 200),
        ylim=(-200, 200),
        zlim=(400, 0),
    )
)

```

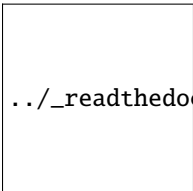

../_readthedocs/jupyter_execute/electrodes_5_0.png

As seen above, the `tile_coords` function can be used to repeat a single shank to produce coordinates for a multi-shank probe. Likewise it can be used to repeat multi-shank coordinates to achieve a 3D recording array (what NeuroNexus calls a [MatrixArray](#)).

For our example we will use a simple linear array. We configure the probe so it has 32 contacts ranging from 0.2 to 1.2 mm in depth. We could specify the orientation, but by default shank coordinates extend downwards (in the positive z direction).

We can add the electrode to the plotting function to visualize it along with the neurons:

```

coords = ephys.linear_shank_coords(1 * mm, 32, start_location=(0, 0, 0.2) * mm)
probe = ephys.Probe(coords, save_history=True)
cleo.viz.plot(
    exc, inh, colors=[c['exc'], c['inh']], zlim=(0, 1.2), devices=[probe], scatterargs={
        ↪ 'alpha': .3}
)

```



```
(<Figure size 640x480 with 1 Axes>,
<Axes3D: xlabel='x (mm)', ylabel='y (mm)', zlabel='z (mm)')>)
```

```
../_readthedocs/jupyter_execute/electrodes_7_1.png
```

Specifying signals to record

This looks right, but we need to specify what signals we want to pick up with our electrode. Let's try the two basic spiking signals and an LFP approximation for point neurons.

The two spiking signals (sorted and multi-unit) take the same parameters, mainly `r_perfect_detection`, within which all spikes will be detected, and `r_half_detection`, at which distance a spike has only a 50% chance of being detected. My choice to set these parameters at 50 and 100 μ m is arbitrary, though from [at least some published data](#) that seems reasonable.

We use default parameters for the Teleńczuk kernel LFP approximation method (TKLFP), but will need to specify cell type (excitatory or inhibitory) and sampling period (if unavailable from a connected IO processor) upon injection.

```
mua = ephys.MultiUnitSpiking(
    r_perfect_detection=0.05 * mm,
    r_half_detection=0.1 * mm,
)
ss = ephys.SortedSpiking(0.05 * mm, 0.1 * mm)
tklfp = ephys.TKLFPSignal()

probe.add_signals(mua, ss, tklfp)

from cleo.ioproc import RecordOnlyProcessor

sim.set_io_processor(RecordOnlyProcessor(sample_period_ms=1))
sim.inject(probe, exc, tklfp_type="exc")
sim.inject(probe, inh, tklfp_type="inh")
```

```
CLSimulator(io_processor=<cleo.ioproc.base.RecordOnlyProcessor object at 0x7f87cac0f790>,
→ devices={Probe(brian_objects={<SpikeMonitor, recording from 'spikemonitor_5'>,
→ <SpikeMonitor, recording from 'spikemonitor_1'>, <SpikeMonitor, recording from
→ 'spikemonitor'>, <SpikeMonitor, recording from 'spikemonitor_4'>, <SpikeMonitor,
→ recording from 'spikemonitor_2'>, <SpikeMonitor, recording from 'spikemonitor_3'>},
→ sim=..., name='Probe', signals=[MultiUnitSpiking(name='MultiUnitSpiking', brian_
→ objects={<SpikeMonitor, recording from 'spikemonitor_3'>, <SpikeMonitor, recording
→ from 'spikemonitor'>}, probe=..., r_perfect_detection=50. * umetre, r_half_
→ detection=100. * umetre, cutoff_probability=0.01, save_history=True),
→ SortedSpiking(name='SortedSpiking', brian_objects={<SpikeMonitor, recording from
→ 'spikemonitor_1'>, <SpikeMonitor, recording from 'spikemonitor_4'>}, probe=..., r_
→ perfect_detection=50. * umetre, r_half_detection=100. * umetre, cutoff_probability=0.
→ 01, save_history=True), TKLFPSignal(name='TKLFPSignal', brian_objects={<SpikeMonitor,
→ recording from 'spikemonitor_2'>, <SpikeMonitor, recording from 'spikemonitor_5'>},
→ probe=..., uLFP_threshold_uV=0.001, save_history=True)], probe=NOTHING)})
```

Simulation and results

Now we'll run the simulation:

```
sim.run(150*ms)
```

And plot the output of the three signals we've recorded:


```
from matplotlib.colors import ListedColormap, LinearSegmentedColormap
fig, axs = plt.subplots(3, 1, figsize=(8, 9), sharex=True)

# assuming all neurons are detectable for c=ss.i >= n_e to work
# in practice this will often not be the case and we'd have to map
# from probe index to neuron group index using ss.i_probe_by_ing.inverse
exc_inh_cmap = ListedColormap([c['exc'], c['inh']])
axs[0].scatter(ss.t_ms, ss.i, marker=".", c=ss.i >= n_e, cmap=exc_inh_cmap)
axs[0].set(title="sorted spikes", ylabel="neuron index")

axs[1].scatter(mua.t_ms, mua.i, marker=".", s=2, c='white')
axs[1].set(title="multi-unit spikes", ylabel="channel index")

lines = axs[2].plot(tklfp.lfp_uV)
axs[2].set(
    title="Teleńczuk kernel LFP approximation", xlabel="t (ms)", ylabel="LFP (V)"
)

# color-code channel depth
from mpl_toolkits.axes_grid1.inset_locator import inset_axes
depth_cmap = LinearSegmentedColormap.from_list('cleo', ['white', c['dark']])
axins = inset_axes(axs[2], width='95%', height='3%', loc="upper center")
for i in range(32):
    l = lines[i]
    l.set_color(depth_cmap(i / 31))
from matplotlib.colors import Normalize
channel_mappable = plt.cm.ScalarMappable(Normalize(0, 1.2), depth_cmap)
fig.colorbar(
    channel_mappable,
    cax=axins,
    orientation="horizontal",
    ticks=[],
    label="channel depth (0.2 to 1.2 mm)",
);
```

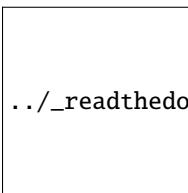

../_readthedocs/jupyter_execute/electrodes_13_0.png

Or, to see the LFP as a function of depth better:

```
fig, axs = plt.subplots(1, 2, figsize=(8, 9))
channel_offsets = -12 * np.arange(32)
```

(continues on next page)


(continued from previous page)

```

lfp_to_plot = tkllfp.lfp_uV + channel_offsets
axs[0].plot(lfp_to_plot, color="w")
axs[0].set(
    yticks=channel_offsets,
    yticklabels=range(1, 33),
    xlabel="t (ms)",
    ylabel="channel #",
)
cmap = LinearSegmentedColormap.from_list('lfp', ['#131416', c['main']])
im = axs[1].imshow(tkllfp.lfp_uV.T, aspect="auto", cmap=cmap)
axs[1].set(xlabel="t (ms)")
fig.colorbar(im, aspect=40, label="LFP (V)")

```

```
<matplotlib.colorbar.Colorbar at 0x7f87ca70dc60>
```

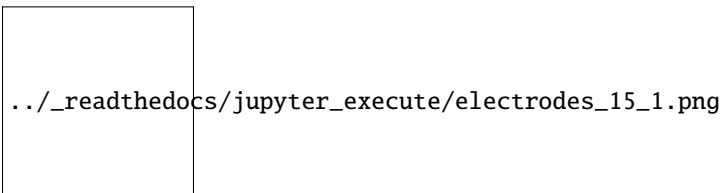

6.2.2 Optogenetic stimulation

How to inject an optogenetic intervention (opsin and optic fiber) into a simulation.

Preamble:

```

%load_ext autoreload
%autoreload 2
from brian2 import *
import matplotlib.pyplot as plt

import cleo
from cleo import *

cleo.utilities.style_plots_for_docs()

# numpy faster than cython for lightweight example
prefs.codegen.target = 'numpy'
# for reproducibility
np.random.seed(1866)

```

The autoreload extension is already loaded. To reload it, use:

```
%reload_ext autoreload
```

Create a Markov opsin-compatible network

Cleo enables two basic approaches to modeling opsin currents. One is a fairly accurate Markov state model and the other is a simple proportional current model. We will look at the Markov model here.

The established Markov opsin models (as presented in [Evans et al., 2016](#)), are conductance-based and so depend on somewhat realistic membrane voltages. Hence, we'll use the [exponential integrate-and-fire \(EIF\) model](#). Note that we follow the conventions used in neuron modeling, where current is positive, rather than the conventions in opsin modeling, where the photocurrent is negative.

We'll use a small neuron group, biased by Poisson input spikes.

```
n = 10
ng = NeuronGroup(
    n,
    """
    dv/dt = -(v - E_L) + Delta_T*exp((v-theta)/Delta_T) + Rm*I) / tau_m : volt
    I : amp
    """,
    threshold="v>30*mV",
    reset="v=-55*mV",
    namespace={
        "tau_m": 20 * ms,
        "Rm": 500 * Mohm,
        "theta": -50 * mV,
        "Delta_T": 2 * mV,
        "E_L": -70*mV,
    },
)
ng.v = -70 * mV

input_group = PoissonInput(ng, "v", n, 100 * Hz, 1 * mV)

mon = SpikeMonitor(ng)

net = Network(ng, input_group, mon)
ng.equations
```

$$\frac{dv}{dt} = \frac{\Delta_T e^{\frac{-\theta+v}{\Delta_T}} + E_L + IRm - v}{I} \quad \begin{array}{l} \text{(unit of } v: \text{ V)} \\ \text{(unit: A)} \end{array}$$

Assign coordinates and configure optogenetic model

The `OptogeneticIntervention` class implements the chosen opsin kinetics model with specified parameters. A standard four-state Markov model as well as channelrhodopsin-2 (ChR2) parameters are included with cleo and are accessible in the `cleo.opto` module. For extending to other models (such as three-state or six-state), see the [source code](#)—the state equations, opsin-specific parameters, and light wavelength-specific parameters (if not using 473-nm blue) would be needed.

For reference, cleo draws heavily on [Foutz et al., 2012](#) for the light propagation model and on [Evans et al., 2016](#) for the opsin kinetics model.

```

from cleo.coords import assign_coords_rand_rect_prism
from cleo.opto import *
from cleo.light import *

assign_coords_rand_rect_prism(ng, xlim=(-0.1, 0.1), ylim=(-0.1, 0.1), zlim=(0.4, 0.6))

opsin = chr2_4s()
fiber = Light(
    coords=(0, 0, 0.2) * mm,
    light_model=fiber473nm(),
    name="fiber",
)

cleo.viz.plot(
    ng,
    colors=["xkcd:fuchsia"],
    xlim=(-0.2, 0.2),
    ylim=(-0.2, 0.2),
    zlim=(0, 0.8),
    devices=[fiber],
)

```

```

(<Figure size 640x480 with 1 Axes>,
 <Axes3D: xlabel='x (mm)', ylabel='y (mm)', zlabel='z (mm)'>)

```

../_readthedocs/jupyter_execute/optogenetics_5_1.png

Open-loop optogenetic stimulation

We need to inject our optogenetic intervention into the simulator. `cleo` handles all the object creation and equations needed to interact with the existing Brian model without the need to alter it, with the possible exception of adding a variable to represent the opsin current. This needs to be specified upon injection with `Iopto_var_name=...` if not the default `Iopto`. The membrane potential variable name also needs to be specified (with `v_var_name=...`) if not the default `v`.

```

sim = CLSimulator(net)
sim.inject(opsin, ng, Iopto_var_name='I')
sim.inject(fiber, ng)

```

```

CLSimulator(io_processor=None, devices={FourStateOpsin(brian_objects=
→ {NeuronGroup(clock=Clock(dt=100. * usecond, name='defaultclock'), when=start, order=0,
→ name='light_agg_ChR2_neurongroup'), Synapses(clock=Clock(dt=100. * usecond, name=
→ 'defaultclock'), when=start, order=0, name='syn_ChR2_neurongroup')}, sim=..., name=
→ 'ChR2', save_history=True, spectrum=[(400, 0.34), (422, 0.65), (460, 0.96), (470, 1),
→ (473, 1), (500, 0.57), (520, 0.22), (540, 0.06), (560, 0.01)], required_vars=[('Iopto',
→ amp), ('v', volt)], g0=114. * nsiemens, gamma=0.00742, phim=2.33e+23 * (second ** -1),
→ / (meter ** 2), k1=4.15 * khertz, k2=0.868 * khertz, p=0.833, Gf0=37.3 * hertz, kf=58.
→ 1 * hertz, Gb0=16.1 * hertz, kb=63. * hertz, q=1.94, Gd1=105. * hertz, Gd2=13.8
→ hertz, Gr0=0.33 * hertz, E=0. * volt, v0=43. * mvolt, v1=17.1 * mvolt, model="\n
→ dG1/dt = Gd1*O1 + Gr0*C2 - Gd1*C1 : 1 (clock-driven)\n          dO1/dt = Gd1*C1 + Gb*O2
→ - (Gd1+Gf)*O1 : 1 (clock-driven)\n          dO2/dt = Ga2*C2 + Gf*O1 - (Gd2+Gb)*O2 : 1
→ (clock-driven)\n          C2 = 1 - C1 - O1 - O2 : 1\n\n          Theta = int(phi_pre >
→ 0*phi_pre) : 1\n          Hp = Theta * phi_pre**p/(phi_pre**p + phim**p) : 1\n
→ Ga1 = k1*Hp : hertz\n          Ga2 = k2*Hp : hertz\n          Hq = Theta * phi_pre**q/(phi

```

6.2. Tutorials 25

```

→ dG1/dt = Gd1*O1 + Gr0*C2 - Gd1*C1 : 1 (clock-driven)\n          dO1/dt = Gd1*C1 + Gb*O2
→ - (Gd1+Gf)*O1 : 1 (clock-driven)\n          dO2/dt = Ga2*C2 + Gf*O1 - (Gd2+Gb)*O2 : 1
→ (clock-driven)\n          C2 = 1 - C1 - O1 - O2 : 1\n\n          Theta = int(phi_pre >
→ 0*phi_pre) : 1\n          Hp = Theta * phi_pre**p/(phi_pre**p + phim**p) : 1\n
→ Ga1 = k1*Hp : hertz\n          Ga2 = k2*Hp : hertz\n          Hq = Theta * phi_pre**q/(phi

```

IO processor setup

Here we design an IO processor that ignores measurements and simply sets the light intensity according to the `stimulus(t)` function:

```
from cleo.ioproc import LatencyIOProcessor

def stimulus(time_ms):
    f = 30
    return 1 * (1 + np.sin(2*np.pi*f * time_ms/1000))

class OpenLoopOpto(LatencyIOProcessor):
    def __init__(self):
        super().__init__(sample_period_ms=1)

    # since this is open-loop, we don't use state_dict
    def process(self, state_dict, time_ms):
        opto_intensity = stimulus(time_ms)
        # return output dict and time
        return {"fiber": opto_intensity, time_ms}

sim.set_io_processor(OpenLoopOpto())
```

```
CLSSimulator(io_processor=<__main__.OpenLoopOpto object at 0x7fa986024280>, devices=
→ {FourStateOpsin(brian_objects={NeuronGroup(clock=Clock(dt=100. * usecond, name=
→ 'defaultclock'), when=start, order=0, name='light_agg_ChR2_neurongroup'),
→ Synapses(clock=Clock(dt=100. * usecond, name='defaultclock'), when=start, order=0,
→ name='syn_ChR2_neurongroup')}, sim=..., name='ChR2', save_history=True, spectrum=[(400,
→ 0.34), (422, 0.65), (460, 0.96), (470, 1), (473, 1), (500, 0.57), (520, 0.22), (540,
→ 0.06), (560, 0.01)], required_vars=[('Iopto', amp), ('v', volt)], g0=114. * nsiemens,
→ gamma=0.00742, phim=2.33e+23 * (second ** -1) / (meter ** 2), k1=4.15 * khertz, k2=0.
→ 868 * khertz, p=0.833, Gf0=37.3 * hertz, kf=58.1 * hertz, Gb0=16.1 * hertz, kb=63. *
→ hertz, q=1.94, Gd1=105. * hertz, Gd2=13.8 * hertz, Gr0=0.33 * hertz, E=0. * volt,
→ v0=43. * mvolt, v1=17.1 * mvolt, model="\n          dC1/dt = Gd1*O1 + Gr0*C2 - Ga1*C1 :
→ 1 (clock-driven)\n          dO1/dt = Ga1*C1 + Gb*O2 - (Gd1+Gf)*O1 : 1 (clock-driven)\n
→          dO2/dt = Ga2*C2 + Gf*O1 - (Gd2+Gb)*O2 : 1 (clock-driven)\n          C2 = 1 - C1 -
→ O1 - O2 : 1\n\n          Theta = int(phi_pre > 0*phi_pre) : 1\n          Hp = Theta * phi_
→ pre**p/(phi_pre**p + phim**p) : 1\n          Ga1 = k1*Hp : hertz\n          Ga2 = k2*Hp :
→ hertz\n          Hq = Theta * phi_pre**q/(phi_pre**q + phim**q) : 1\n          Gf = kf*Hq
→ + Gf0 : hertz\n          Gb = kb*Hq + Gb0 : hertz\n\n          fphi = O1 + gamma*O2 : 1\n
→          # v1/v0 when v-E == 0 via l'Hopital's rule\n          fv = f_unless_x0(\n
→          (1 - exp(-(V_VAR_NAME_post-E)/v0)) / ((V_VAR_NAME_post-E)/v1),\n          V_VAR_
→ NAME_post - E,\n          v1/v0\n          ) : 1\n\n          IOPTO_VAR_NAME_post = -
→ g0*fphi*f_v*(V_VAR_NAME_post-E)*rho_rel : ampere (summed)\n          rho_rel : 1", extra_
→ namespace={'f_unless_x0': <brian2.core.functions.Function object at 0x7fa988672c80>}),
→ Light(brian_objects=set(), sim=..., name='fiber', save_history=True, value=array([0.]),
→ light_model=OpticFiber(R0=100. * umetre, NAfib=0.37, K=125. * metre ** -1, S=7370. *
→ metre ** -1, ntis=1.36), coords=array([[ 0.,  0., 200.]]) * umetre, wavelength=0.473
→ * umetre, direction=array([0., 0., 1.]), max_Irr0_mW_per_mm2=None, max_Irr0_mW_per_mm2_
→ viz=None, default_value=array([0.])))
```

(continues on next page)


(continued from previous page)

Run simulation and plot results

```
sim.run(100*ms)

fig, (ax1, ax2) = plt.subplots(2, 1, sharex=True)
ax1.plot(mon.t / ms, mon.i[:, '|'])
ax1.set(ylabel='neuron index', title='spiking')
t_sim=np.linspace(0, 100, 1000)
ax2.plot(t_sim, stimulus(t_sim))
ax2.set(ylabel=r'$I_{rr_0}$ (mm/mW$^2$)', title='optogenetic stimulus', xlabel='time (ms)');
```

```
INFO      No numerical integration method specified for group 'neurongroup', using
↳method 'euler' (took 0.01s, trying other methods took 0.04s). [brian2.stateupdaters.
↳base.method_choice]
INFO      No numerical integration method specified for group 'syn_ChR2_neurongroup',
↳using method 'euler' (took 0.02s, trying other methods took 0.05s). [brian2.
↳stateupdaters.base.method_choice]
```

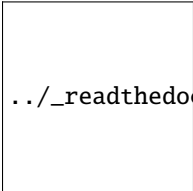

../_readthedocs/jupyter_execute/optogenetics_11_1.png

```
opsin.synapses[ng.name].equations
```

$$\begin{aligned}
C2 &= -C_1 - O_1 - O_2 + 1 & (\text{unit of } C_2: \text{ rad}) \\
Theta &= \text{int}(\phi_{pre} > 0) & (\text{unit of } \Theta: \text{ rad}) \\
fphi &= O_1 + O_2\gamma & (\text{unit of } fphi: \text{ rad}) \\
fv &= f_{unlessx0} \left(\frac{v_1 \cdot \left(1 - e^{\frac{E - v_{post}}{v_0}} \right)}{-E + v_{post}}, -E + v_{post}, \frac{v_1}{v_0} \right) & (\text{unit of } fv: \text{ rad}) \\
Hp &= \frac{\Theta \phi_{pre}^p}{\phi_{pre}^p + phimp} & (\text{unit of } Hp: \text{ rad}) \\
Hq &= \frac{\Theta \phi_{pre}^q}{\phi_{pre}^q + phimq} & (\text{unit of } Hq: \text{ rad}) \\
Ga1 &= Hpk_1 & (\text{unit of } Ga_1: \text{ Hz}) \\
Ga2 &= Hpk_2 & (\text{unit of } Ga_2: \text{ Hz}) \\
Gf &= Gf_0 + Hqkf & (\text{unit of } Gf: \text{ Hz}) \\
Gb &= Gb_0 + Hqkb & (\text{unit of } Gb: \text{ Hz}) \\
\frac{dC_1}{dt} &= -C_1 Ga_1 + C_2 Gr_0 + Gd_1 O_1 & (\text{unit of } C_1: \text{ rad, flags: clock-driven}) \\
\frac{dO_1}{dt} &= C_1 Ga_1 + Gb O_2 - O_1 (Gd_1 + Gf) & (\text{unit of } O_1: \text{ rad, flags: clock-driven}) \\
\frac{dO_2}{dt} &= C_2 Ga_2 + Gf O_1 - O_2 (Gb + Gd_2) & (\text{unit of } O_2: \text{ rad, flags: clock-driven}) \\
\rho_{rel} & & (\text{unit: rad})
\end{aligned}$$

Conclusion

We can see clearly that firing rate correlates with light intensity as expected.

As a recap, in this tutorial we've seen how to:

- configure an `OptogeneticIntervention`,
- inject it into the simulation,
- and control its light intensity in an open-loop fashion.

Appendix: alternative opsin and neuron models

Because it would be a pain and an obstacle to reproducibility to have to replace all pre-existing simple neuron models with more sophisticated ones with proper voltage ranges and units, we provide an approximation that is much more flexible, requiring only a current term, of any unit, in the target neurons.

The Markov models of opsin dynamics we've used so far produce a rise, peak, and fall to a steady-state plateau current when subjected to sustained light. Since they are conductance-based, the current also varies with membrane voltage, including during spikes. The `ProportionalCurrentModel`, on the other hand, simply delivers current proportional to light intensity. This should be adequate for a wide range of use cases where the exact opsin current dynamics on short timescales don't matter so much and a sort of average current-light relationship will suffice.

Speaking of realistic membrane voltages, does the Markov model's voltage-dependent current render it unsuitable for the most basic leaky integrate-and-fire (LIF) neuron model? LIF neurons reset on reaching their rheobase threshold,

staying perpetually in a subthreshold region producing exaggerated opsin currents. How much does this affect the output? We will explore this question by comparing a variety of opsin/neuron model combinations.

First, we introduce exponential integrate-and-fire neurons, which maintain simplicity while modeling an upward membrane potential swing during a spike. For more info, see the [related section in the Neuronal Dynamics online textbook](#) and their [example parameters table](#).

```
neuron_params = {
    "tau_m": 20 * ms,
    "Rm": 500 * Mohm,
    "theta": -50 * mV,
    "Delta_T": 2 * mV,
    "E_L": -70*mV,
}

def prep_ng(ng, neuron_type, markov_opsin):
    ng.v = neuron_params['E_L']
    assign_coords_rand_rect_prism(ng, xlim=(0, 0), ylim=(0, 0), zlim=(0, 0))
    state_mon = StateMonitor(ng, ("Iopto", "v"), record=True)
    spike_mon = SpikeMonitor(ng)
    return neuron_type, ng, state_mon, spike_mon, markov_opsin

experiments = []

eif = NeuronGroup(
    1,
    """
    dv/dt = -(v - E_L) + Delta_T*exp((v-theta)/Delta_T) + Rm*Iopto) / tau_m : volt
    Iopto : amp
    """,
    threshold="v > -10*mV",
    reset="v = E_L - 0*mV",
    namespace=neuron_params,
)

experiments.append(prepare_ng(eif, 'EIF', True))
```

Configure LIF models

Here we define LIF neurons with biological parameters for the sake of comparison, but the ProportionalCurrentModel is compatible with models of any voltage range and units, so long as it has an Iopto term.

```
def prep_lif(markov_opsin):
    ng = NeuronGroup(
        1,
        """dv/dt = -(v - E_L) + Rm*Iopto) / tau_m : volt
        Iopto : amp""",
        threshold="v > theta + 4*mV",
        reset="v = E_L - 0*mV",
        namespace=neuron_params,
    )
    return prep_ng(ng, "LIF", markov_opsin)
```

(continues on next page)

(continued from previous page)

```
experiments.append(prepare_lif(True))
experiments.append(prepare_lif(False))
```

Comparing to more realistic models

To see how well simplified neuron and opsin models do, we'll also compare to the more complex `AdEx neuron` (with “tonic” firing pattern parameters) and a Hodgkin-Huxley model (code from `Neuronal Dynamics`).

```
adex = NeuronGroup(
    1,
    """dv/dt = -(v - E_L) + 2*mV*exp((v-theta)/Delta_T) + Rm*(Iopto-w) / tau_m : volt
    dw/dt = (0*nsiemens*(v-E_L) - w) / (100*ms) : amp
    Iopto : amp""",
    threshold="v>=-10*mV",
    reset="v=-55*mV; w+=5*pamp",
    namespace=neuron_params,
)
experiments.append(prepare_ng(adex, "AdEx", True))

# Parameters
# Cm = 1*ufarad*cm**-2 * area
Cm = neuron_params["tau_m"] / neuron_params["Rm"]
# area = 5000*umetre**2
area = Cm / (1*ufarad*cm**-2)
gl = 0.3*msiemens*cm**-2 * area
El = -65*mV
EK = -90*mV
ENa = 50*mV
g_na = 40*msiemens*cm**-2 * area
g_kd = 35*msiemens*cm**-2 * area
VT = -63*mV

# The model
eqs = Equations('''
dv/dt = (gl*(El-v) - g_na*(m*m*m)*h*(v-ENa) - g_kd*(n*n*n*n)*(v-EK) + Iopto)/Cm : volt
dm/dt = 0.32*(mV**-1)*4*mV/exprel((13.*mV-v+VT)/(4*mV))/ms*(1-m)-0.28*(mV**-1)*5*mV/
    exprel((v-VT-40.*mV)/(5*mV))/ms*m : 1
dn/dt = 0.032*(mV**-1)*5*mV/exprel((15.*mV-v+VT)/(5*mV))/ms*(1.-n)-.5*exp((10.*mV-v+VT)/
    (40.*mV))/ms*n : 1
dh/dt = 0.128*exp((17.*mV-v+VT)/(18.*mV))/ms*(1.-h)-4./(1+exp((40.*mV-v+VT)/(5.*mV)))/
    ms*h : 1
Iopto : amp
''')
# Threshold and refractoriness are only used for spike counting
hh = NeuronGroup(1, eqs,
    threshold='v > -40*mV',
    reset='',
    method='exponential_euler')
experiments.append(prepare_ng(hh, "HH", True))
```

Opsin configuration

Note that the only parameter we need to set for the simple opsin model is the gain on light intensity, `I_per_Irr`. This term defines what the neuron receives for every 1 mW/mm2 of light intensity. Here that term is defined in amperes, but it could have been unitless for a simpler model.

The gain is tuned somewhat by hand (in relation to the membrane resistance and the 20 mV gap between rest and threshold potential) to achieve similar outputs to the Markov model.

```
light = Light(light_model=fiber473nm())
simple_opsin = ProportionalCurrentOpsin(
    name="simple_opsin",
    # handpicked gain to make firing rate roughly comparable to EIF
    I_per_Irr=140/neuron_params['Rm']*20*mV,
)
markov_opsin = chr2_4s()
markov_opsin.name = "markov_opsin"
```

Simulation

And we set up the simulator:

```
net = Network()
sim = CLSimulator(net)
for ng_type, ng, state_mon, spike_mon, use_markov_opsin in experiments:
    net.add(ng, state_mon, spike_mon)
    sim.inject(light, ng)
    if use_markov_opsin:
        sim.inject(markov_opsin, ng)
    else:
        sim.inject(simple_opsin, ng)
```

We'll now run the simulation with light pulses of increasing amplitudes to observe the effect on the current.

```
# hand-picked range of amplitudes to show 0 to moderate firing rates
for Irr0_mW_per_mm2 in np.linspace(0.015, 0.05, 5):
    light.update(Irr0_mW_per_mm2)
    sim.run(60 * ms)
    light.update(0)
    sim.run(60 * ms)
```

```
INFO      No numerical integration method specified for group 'neurongroup_1', using
↳method 'euler' (took 0.01s, trying other methods took 0.02s). [brian2.stateupdaters.
↳base.method_choice]
INFO      No numerical integration method specified for group 'neurongroup_2', using
↳method 'exact' (took 0.04s). [brian2.stateupdaters.base.method_choice]
INFO      No numerical integration method specified for group 'neurongroup_3', using
↳method 'exact' (took 0.02s). [brian2.stateupdaters.base.method_choice]
INFO      No numerical integration method specified for group 'neurongroup_4', using
↳method 'euler' (took 0.01s, trying other methods took 0.05s). [brian2.stateupdaters.
↳base.method_choice]
```

```

WARNING      'n' is an internal variable of group 'neurongroup_5', but also exists in the_
↳run namespace with the value 10. The internal variable will be used. [brian2.groups.
↳group.Group.resolve.resolution_conflict]
INFO         No numerical integration method specified for group 'syn_markov_opsin_
↳neurongroup_1', using method 'euler' (took 0.01s, trying other methods took 0.01s).↳
↳[brian2.stateupdaters.base.method_choice]
INFO         No numerical integration method specified for group 'syn_markov_opsin_
↳neurongroup_2', using method 'euler' (took 0.01s, trying other methods took 0.01s).↳
↳[brian2.stateupdaters.base.method_choice]
INFO         No numerical integration method specified for group 'syn_markov_opsin_
↳neurongroup_4', using method 'euler' (took 0.01s, trying other methods took 0.01s).↳
↳[brian2.stateupdaters.base.method_choice]
INFO         No numerical integration method specified for group 'syn_markov_opsin_
↳neurongroup_5', using method 'euler' (took 0.01s, trying other methods took 0.01s).↳
↳[brian2.stateupdaters.base.method_choice]

```

Results

```

c1 = '#8000b4'
c2 = '#df87e1'

fig, axs = plt.subplots(
    len(experiments), 1, figsize=(8, 2*len(experiments)), sharex=True
)

for ax, (ng_type, _, state_mon, spike_mon, markov_opsin) in zip(axs, experiments):
    ax.plot(state_mon.t / ms, state_mon.Iopto[0] / pamp, c=c1, label="$I_{opto}$ (pA)")
    ax.plot(state_mon.t / ms, state_mon.v[0] / mV, c=c2, label="v (mV)")
    opsin_name = "4-state Markov" if markov_opsin else "proportional current"
    ax.set(title=f"{ng_type} neuron, {opsin_name} opsin model")

axs[-1].set_xlabel('t (ms)')
axs[0].legend();

max_ylim = max([ax.get_ylim()[1] for ax in axs])
for ax in axs:
    ax.set_ylim([-75, 200])

fig.tight_layout()

```

../_readthedocs/jupyter_execute/optogenetics_28_0.png

Qualitatively we can see that the proportional current model doesn't capture the rise, peak, plateau, and fall dynamics that a Markov model can produce, but is a reasonable approximation if all you need is a roughly linear light intensity-firing rate relationship. We also see that a variety of neuron/opsin model combinations all produce similar firing responses to light.

6.2.3 Multi-channel, bidirectional optogenetics

Cleo supports the simultaneous use of multiple Light devices, multiple channels per device, and multiple opsins per neuron group. Here we'll see how to use all these features.

```
# boilerplate
%load_ext autoreload
%autoreload 2
from brian2 import *
import matplotlib.pyplot as plt

import cleo
from cleo import *

cleo.utilities.style_plots_for_docs()

# numpy faster than cython for lightweight example
prefs.codegen.target = 'numpy'
# for reproducibility
np.random.seed(1866)

# colors
c = {
    'main': '#C500CC',
    '473nm': '#72b5f2',
    '590nm': (1, .875, 0),
}
```

The autoreload extension is already loaded. To reload it, use:

```
%reload_ext autoreload
```

Network setup

We'll use excitatory and inhibitory populations of exponential integrate-and-fire neurons.

```
n = 500
ng = NeuronGroup(
    n,
    """
    dv/dt = (
        -(v - E_L)
        + Delta_T*exp((v-theta)/Delta_T)
        + Rm*(I_exc + I_inh + I_bg)
    ) / tau_m + sigma*sqrt(2/tau_m)*xi: volt
    I_exc : amp
    I_inh : amp
    """,
    threshold="v > -10*mV",
    reset="v=E_L",
    namespace={
        "tau_m": 20 * ms,
        "Rm": 500 * Mohm,
```

(continues on next page)

(continued from previous page)

```

        "theta": -50 * mV,
        "Delta_T": 2 * mV,
        "E_L": -70*mV,
        "sigma": 5 * mV,
        "I_bg": 60 * pamp,
    },
)
ng.v = np.random.uniform(-70, -50, n) * mV

W = 1 * mV
p_S = 0.3
n_neighbors = 40
S_ee = Synapses(ng, model='w: 1', on_pre="v_post+=W*w/sqrt(N)")
S_ee.connect(condition='abs(i-j)<=n_neighbors and i!=j')
S_ee.w = np.exp(np.random.randn(int(S_ee.N - 0))) * np.random.choice([-1, 1], int(S_ee.N - 0))
    ↪ - 0))

spike_mon = SpikeMonitor(ng)
# for visualizing currents
Imon = StateMonitor(ng, ('I_exc', 'I_inh'), record=range(50))
net = Network(ng, S_ee, spike_mon, Imon)
sim = cleo.CLSimulator(net)

```

```

from cleo.coords import assign_coords_grid_rect_prism
xmax_mm = 2
assign_coords_grid_rect_prism(ng, xlim=(0, xmax_mm), ylim=(-.15, .15), zlim=(0, .3),
    ↪ shape=(20, 5, 5))
cleo.viz.plot(ng, colors=[c['main']])

```

```

(<Figure size 640x480 with 1 Axes>,
 <Axes3D: xlabel='x (um)', ylabel='y (um)', zlabel='z (um)'>)

```

```

../_readthedocs/jupyter_execute/multi_opto_4_1.png

```

Injecting a multi-channel Light

A `Light` device can have multiple channels; all the user needs is to specify the coordinates (and optionally direction) of each light source (channel). A `LightModel` (e.g., that of an optical fiber) defines how light propagates from each source.

Here we inject 590 nm light for activating Vf-Chrimson. Lacking a more rigorous quantification, we assume absorption and scattering coefficients of 590 nm light in the brain are roughly 0.8 times that of 470 nm light (see [Jacques 2013](#) for some justification).

```

from cleo.light import Light, OpticFiber
from cleo.opto import vfchrimson_4s

```

(continues on next page)

(continued from previous page)

```

n_fibers = 4
coords = np.zeros((n_fibers, 3))
end_space = 1 / (2 * n_fibers)
coords[:, 0] = np.linspace(end_space, 1 - end_space, n_fibers) * xmax_mm
coords[:, 1] = -0.025
amber_fibers = Light(
    name="amber fibers",
    coords=coords * mm,
    light_model=OpticFiber(K=0.125 * 0.8 / mm, S=7.37 * 0.8 / mm),
    wavelength=590 * nmeter,
)
sim.inject(amber_fibers, ng)
vfc = vfchromson_4s()
sim.inject(vfc, ng, Iopto_var_name="I_exc")
cleo.viz.plot(ng, colors=[c["main"]], sim=sim)

```

```

(<Figure size 640x480 with 1 Axes>,
 <Axes3D: xlabel='x (um)', ylabel='y (um)', zlabel='z (um)'>)

```

../_readthedocs/jupyter_execute/multi_opto_6_1.png

Bidirectional control via a second opsin

Here we will demonstrate increasing and decreasing activity in the same experiment by injecting an inhibitory opsin with a minimally overlapping activation spectrum. Alternatively, we could achieve bidirectional control with excitatory opsins on excitatory and inhibitory neurons.

We will use the anion channel GtACR2 which is maximally activated by 470 nm light.

Let's first visualize how the action spectra for the two opsins overlap:

```

gtacr2 = cleo.opto.gtacr2_4s()
cleo.light.plot_spectra(vfc, gtacr2)

```

```

(<Figure size 640x480 with 1 Axes>,
 <Axes: title={'center': 'Action/excitation spectra'}, xlabel=' (nm)', ylabel=''>)

```

../_readthedocs/jupyter_execute/multi_opto_8_1.png

We see that GtACR2 will be totally unaffected by amber (590 nm) light, while Vf-Chrimson will be slightly activated by blue (470 nm) light.

```

coords[:, 1] = 0.025
blue_fibers = Light(
    name="blue fibers",
    coords=coords * mm,
    light_model=cleo.light.fiber473nm(),
)
sim.inject(blue_fibers, ng)
sim.inject(gtacr2, ng, Iopto_var_name="I_inh")
cleo.viz.plot(ng, colors=[c["main"]], sim=sim)

```

```

WARNING /home/kyle/Dropbox (GaTech)/projects/cleo/cleo/light/light_dependence.py:107:
↳UserWarning: = 590.0 nm is outside the range of the action spectrum data for GtACR2.
↳Assuming = 0.
    warnings.warn(
[py.warnings]

```

```

(<Figure size 640x480 with 1 Axes>,
 <Axes3D: xlabel='x (um)', ylabel='y (um)', zlabel='z (um)'>)

```

../_readthedocs/jupyter_execute/multi_opto_10_2.png

Open-loop stimulation

We will now design a stimulus pattern to demonstrate bidirectional control segregated by channel.

```

from cleo.ioproc import LatencyIOProcessor

class OpenLoopOpto(LatencyIOProcessor):
    def __init__(self):
        super().__init__(sample_period_ms=1)

    # since this is open-loop, we don't use state_dict
    def process(self, state_dict, time_ms):
        amplitude_mW_mm2 = 1
        time_offsets = np.array([0, -20, -40, -60])
        t = time_ms + time_offsets
        amber = ((t >= 20) & (t < 40)) * amplitude_mW_mm2
        blue = ((t >= 60) & (t < 63)) * amplitude_mW_mm2

        # return output dict and time
        return ({"amber fibers": amber, "blue fibers": blue}, time_ms)

sim.set_io_processor(OpenLoopOpto())

```

```

CLSSimulator(io_processor=<__main__.OpenLoopOpto object at 0x7fe758a361d0>, devices=
↳{BansalFourStateOpsin(sim=..., name='GtACR2', save_history=True, on_pre='',
↳spectrum=[(400, 0.4), (410, 0.49), (420, 0.56), (430, 0.65), (440, 0.82), (450, 0.88),
↳(460, 0.88), (470, 1.0), (480, 0.91), (490, 0.67), (500, 0.41), (510, 0.21), (520, 0.
↳12), (530, 0.06), (540, 0.02), (550, 0.0), (560, 0.0)], required_vars=[('Iopto', amp),
36('v', volt)], Gd1=17. * hertz, Gd2=10. * hertz, Gr0=0.5 * hertz, Gr1=0.5 * hertz,
↳phim=2.e+23 * (second ** -1) / (meter ** 2), k1=40. * khertz, k2=20. * khertz, Gf0=1.
↳* hertz, Gb0=3. * hertz, kf=1. * hertz, kb=5. * hertz, gamma=0.05, p=1, q=0.1, E=-69.5
↳* mvolt, model='n
dC1/dt = Gd1*C01 + Gr0*C2 - Ga1*C1 : 1 (clock-driven)\n
dC2/dt = Gd2*C02 + Gr1*C1 - Ga2*C2 : 1 (clock-driven)\n
dC3/dt = Gd3*C03 + Gr2*C2 - Ga3*C3 : 1 (clock-driven)\n
dC4/dt = Gd4*C04 + Gr3*C3 - Ga4*C4 : 1 (clock-driven)\n
dC5/dt = Gd5*C05 + Gr4*C4 - Ga5*C5 : 1 (clock-driven)\n
dC6/dt = Gd6*C06 + Gr5*C5 - Ga6*C6 : 1 (clock-driven)\n
dC7/dt = Gd7*C07 + Gr6*C6 - Ga7*C7 : 1 (clock-driven)\n
dC8/dt = Gd8*C08 + Gr7*C7 - Ga8*C8 : 1 (clock-driven)\n
dC9/dt = Gd9*C09 + Gr8*C8 - Ga9*C9 : 1 (clock-driven)\n
dC10/dt = Gd10*C10 + Gr9*C9 - Ga10*C10 : 1 (clock-driven)\n
dC11/dt = Gd11*C11 + Gr10*C10 - Ga11*C11 : 1 (clock-driven)\n
dC12/dt = Gd12*C12 + Gr11*C11 - Ga12*C12 : 1 (clock-driven)\n
dC13/dt = Gd13*C13 + Gr12*C12 - Ga13*C13 : 1 (clock-driven)\n
dC14/dt = Gd14*C14 + Gr13*C13 - Ga14*C14 : 1 (clock-driven)\n
dC15/dt = Gd15*C15 + Gr14*C14 - Ga15*C15 : 1 (clock-driven)\n
dC16/dt = Gd16*C16 + Gr15*C15 - Ga16*C16 : 1 (clock-driven)\n
dC17/dt = Gd17*C17 + Gr16*C16 - Ga17*C17 : 1 (clock-driven)\n
dC18/dt = Gd18*C18 + Gr17*C17 - Ga18*C18 : 1 (clock-driven)\n
dC19/dt = Gd19*C19 + Gr18*C18 - Ga19*C19 : 1 (clock-driven)\n
dC20/dt = Gd20*C20 + Gr19*C19 - Ga20*C20 : 1 (clock-driven)\n
dC21/dt = Gd21*C21 + Gr20*C20 - Ga21*C21 : 1 (clock-driven)\n
dC22/dt = Gd22*C22 + Gr21*C21 - Ga22*C22 : 1 (clock-driven)\n
dC23/dt = Gd23*C23 + Gr22*C22 - Ga23*C23 : 1 (clock-driven)\n
dC24/dt = Gd24*C24 + Gr23*C23 - Ga24*C24 : 1 (clock-driven)\n
dC25/dt = Gd25*C25 + Gr24*C24 - Ga25*C25 : 1 (clock-driven)\n
dC26/dt = Gd26*C26 + Gr25*C25 - Ga26*C26 : 1 (clock-driven)\n
dC27/dt = Gd27*C27 + Gr26*C26 - Ga27*C27 : 1 (clock-driven)\n
dC28/dt = Gd28*C28 + Gr27*C27 - Ga28*C28 : 1 (clock-driven)\n
dC29/dt = Gd29*C29 + Gr28*C28 - Ga29*C29 : 1 (clock-driven)\n
dC30/dt = Gd30*C30 + Gr29*C29 - Ga30*C30 : 1 (clock-driven)\n
dC31/dt = Gd31*C31 + Gr30*C30 - Ga31*C31 : 1 (clock-driven)\n
dC32/dt = Gd32*C32 + Gr31*C31 - Ga32*C32 : 1 (clock-driven)\n
dC33/dt = Gd33*C33 + Gr32*C32 - Ga33*C33 : 1 (clock-driven)\n
dC34/dt = Gd34*C34 + Gr33*C33 - Ga34*C34 : 1 (clock-driven)\n
dC35/dt = Gd35*C35 + Gr34*C34 - Ga35*C35 : 1 (clock-driven)\n
dC36/dt = Gd36*C36 + Gr35*C35 - Ga36*C36 : 1 (clock-driven)\n
dC37/dt = Gd37*C37 + Gr36*C36 - Ga37*C37 : 1 (clock-driven)\n
dC38/dt = Gd38*C38 + Gr37*C37 - Ga38*C38 : 1 (clock-driven)\n
dC39/dt = Gd39*C39 + Gr38*C38 - Ga39*C39 : 1 (clock-driven)\n
dC40/dt = Gd40*C40 + Gr39*C39 - Ga40*C40 : 1 (clock-driven)\n
dC41/dt = Gd41*C41 + Gr40*C40 - Ga41*C41 : 1 (clock-driven)\n
dC42/dt = Gd42*C42 + Gr41*C41 - Ga42*C42 : 1 (clock-driven)\n
dC43/dt = Gd43*C43 + Gr42*C42 - Ga43*C43 : 1 (clock-driven)\n
dC44/dt = Gd44*C44 + Gr43*C43 - Ga44*C44 : 1 (clock-driven)\n
dC45/dt = Gd45*C45 + Gr44*C44 - Ga45*C45 : 1 (clock-driven)\n
dC46/dt = Gd46*C46 + Gr45*C45 - Ga46*C46 : 1 (clock-driven)\n
dC47/dt = Gd47*C47 + Gr46*C46 - Ga47*C47 : 1 (clock-driven)\n
dC48/dt = Gd48*C48 + Gr47*C47 - Ga48*C48 : 1 (clock-driven)\n
dC49/dt = Gd49*C49 + Gr48*C48 - Ga49*C49 : 1 (clock-driven)\n
dC50/dt = Gd50*C50 + Gr49*C49 - Ga50*C50 : 1 (clock-driven)\n
dC51/dt = Gd51*C51 + Gr50*C50 - Ga51*C51 : 1 (clock-driven)\n
dC52/dt = Gd52*C52 + Gr51*C51 - Ga52*C52 : 1 (clock-driven)\n
dC53/dt = Gd53*C53 + Gr52*C52 - Ga53*C53 : 1 (clock-driven)\n
dC54/dt = Gd54*C54 + Gr53*C53 - Ga54*C54 : 1 (clock-driven)\n
dC55/dt = Gd55*C55 + Gr54*C54 - Ga55*C55 : 1 (clock-driven)\n
dC56/dt = Gd56*C56 + Gr55*C55 - Ga56*C56 : 1 (clock-driven)\n
dC57/dt = Gd57*C57 + Gr56*C56 - Ga57*C57 : 1 (clock-driven)\n
dC58/dt = Gd58*C58 + Gr57*C57 - Ga58*C58 : 1 (clock-driven)\n
dC59/dt = Gd59*C59 + Gr58*C58 - Ga59*C59 : 1 (clock-driven)\n
dC60/dt = Gd60*C60 + Gr59*C59 - Ga60*C60 : 1 (clock-driven)\n
dC61/dt = Gd61*C61 + Gr60*C60 - Ga61*C61 : 1 (clock-driven)\n
dC62/dt = Gd62*C62 + Gr61*C61 - Ga62*C62 : 1 (clock-driven)\n
dC63/dt = Gd63*C63 + Gr62*C62 - Ga63*C63 : 1 (clock-driven)\n
dC64/dt = Gd64*C64 + Gr63*C63 - Ga64*C64 : 1 (clock-driven)\n
dC65/dt = Gd65*C65 + Gr64*C64 - Ga65*C65 : 1 (clock-driven)\n
dC66/dt = Gd66*C66 + Gr65*C65 - Ga66*C66 : 1 (clock-driven)\n
dC67/dt = Gd67*C67 + Gr66*C66 - Ga67*C67 : 1 (clock-driven)\n
dC68/dt = Gd68*C68 + Gr67*C67 - Ga68*C68 : 1 (clock-driven)\n
dC69/dt = Gd69*C69 + Gr68*C68 - Ga69*C69 : 1 (clock-driven)\n
dC70/dt = Gd70*C70 + Gr69*C69 - Ga70*C70 : 1 (clock-driven)\n
dC71/dt = Gd71*C71 + Gr70*C70 - Ga71*C71 : 1 (clock-driven)\n
dC72/dt = Gd72*C72 + Gr71*C71 - Ga72*C72 : 1 (clock-driven)\n
dC73/dt = Gd73*C73 + Gr72*C72 - Ga73*C73 : 1 (clock-driven)\n
dC74/dt = Gd74*C74 + Gr73*C73 - Ga74*C74 : 1 (clock-driven)\n
dC75/dt = Gd75*C75 + Gr74*C74 - Ga75*C75 : 1 (clock-driven)\n
dC76/dt = Gd76*C76 + Gr75*C75 - Ga76*C76 : 1 (clock-driven)\n
dC77/dt = Gd77*C77 + Gr76*C76 - Ga77*C77 : 1 (clock-driven)\n
dC78/dt = Gd78*C78 + Gr77*C77 - Ga78*C78 : 1 (clock-driven)\n
dC79/dt = Gd79*C79 + Gr78*C78 - Ga79*C79 : 1 (clock-driven)\n
dC80/dt = Gd80*C80 + Gr79*C79 - Ga80*C80 : 1 (clock-driven)\n
dC81/dt = Gd81*C81 + Gr80*C80 - Ga81*C81 : 1 (clock-driven)\n
dC82/dt = Gd82*C82 + Gr81*C81 - Ga82*C82 : 1 (clock-driven)\n
dC83/dt = Gd83*C83 + Gr82*C82 - Ga83*C83 : 1 (clock-driven)\n
dC84/dt = Gd84*C84 + Gr83*C83 - Ga84*C84 : 1 (clock-driven)\n
dC85/dt = Gd85*C85 + Gr84*C84 - Ga85*C85 : 1 (clock-driven)\n
dC86/dt = Gd86*C86 + Gr85*C85 - Ga86*C86 : 1 (clock-driven)\n
dC87/dt = Gd87*C87 + Gr86*C86 - Ga87*C87 : 1 (clock-driven)\n
dC88/dt = Gd88*C88 + Gr87*C87 - Ga88*C88 : 1 (clock-driven)\n
dC89/dt = Gd89*C89 + Gr88*C88 - Ga89*C89 : 1 (clock-driven)\n
dC90/dt = Gd90*C90 + Gr89*C89 - Ga90*C90 : 1 (clock-driven)\n
dC91/dt = Gd91*C91 + Gr90*C90 - Ga91*C91 : 1 (clock-driven)\n
dC92/dt = Gd92*C92 + Gr91*C91 - Ga92*C92 : 1 (clock-driven)\n
dC93/dt = Gd93*C93 + Gr92*C92 - Ga93*C93 : 1 (clock-driven)\n
dC94/dt = Gd94*C94 + Gr93*C93 - Ga94*C94 : 1 (clock-driven)\n
dC95/dt = Gd95*C95 + Gr94*C94 - Ga95*C95 : 1 (clock-driven)\n
dC96/dt = Gd96*C96 + Gr95*C95 - Ga96*C96 : 1 (clock-driven)\n
dC97/dt = Gd97*C97 + Gr96*C96 - Ga97*C97 : 1 (clock-driven)\n
dC98/dt = Gd98*C98 + Gr97*C97 - Ga98*C98 : 1 (clock-driven)\n
dC99/dt = Gd99*C99 + Gr98*C98 - Ga99*C99 : 1 (clock-driven)\n
dC100/dt = Gd100*C100 + Gr99*C99 - Ga100*C100 : 1 (clock-driven)\n
dC101/dt = Gd101*C101 + Gr100*C100 - Ga101*C101 : 1 (clock-driven)\n
dC102/dt = Gd102*C102 + Gr101*C101 - Ga102*C102 : 1 (clock-driven)\n
dC103/dt = Gd103*C103 + Gr102*C102 - Ga103*C103 : 1 (clock-driven)\n
dC104/dt = Gd104*C104 + Gr103*C103 - Ga104*C104 : 1 (clock-driven)\n
dC105/dt = Gd105*C105 + Gr104*C104 - Ga105*C105 : 1 (clock-driven)\n
dC106/dt = Gd106*C106 + Gr105*C105 - Ga106*C106 : 1 (clock-driven)\n
dC107/dt = Gd107*C107 + Gr106*C106 - Ga107*C107 : 1 (clock-driven)\n
dC108/dt = Gd108*C108 + Gr107*C107 - Ga108*C108 : 1 (clock-driven)\n
dC109/dt = Gd109*C109 + Gr108*C108 - Ga109*C109 : 1 (clock-driven)\n
dC110/dt = Gd110*C110 + Gr109*C109 - Ga110*C110 : 1 (clock-driven)\n
dC111/dt = Gd111*C111 + Gr110*C110 - Ga111*C111 : 1 (clock-driven)\n
dC112/dt = Gd112*C112 + Gr111*C111 - Ga112*C112 : 1 (clock-driven)\n
dC113/dt = Gd113*C113 + Gr112*C112 - Ga113*C113 : 1 (clock-driven)\n
dC114/dt = Gd114*C114 + Gr113*C113 - Ga114*C114 : 1 (clock-driven)\n
dC115/dt = Gd115*C115 + Gr114*C114 - Ga115*C115 : 1 (clock-driven)\n
dC116/dt = Gd116*C116 + Gr115*C115 - Ga116*C116 : 1 (clock-driven)\n
dC117/dt = Gd117*C117 + Gr116*C116 - Ga117*C117 : 1 (clock-driven)\n
dC118/dt = Gd118*C118 + Gr117*C117 - Ga118*C118 : 1 (clock-driven)\n
dC119/dt = Gd119*C119 + Gr118*C118 - Ga119*C119 : 1 (clock-driven)\n
dC120/dt = Gd120*C120 + Gr119*C119 - Ga120*C120 : 1 (clock-driven)\n
dC121/dt = Gd121*C121 + Gr120*C120 - Ga121*C121 : 1 (clock-driven)\n
dC122/dt = Gd122*C122 + Gr121*C121 - Ga122*C122 : 1 (clock-driven)\n
dC123/dt = Gd123*C123 + Gr122*C122 - Ga123*C123 : 1 (clock-driven)\n
dC124/dt = Gd124*C124 + Gr123*C123 - Ga124*C124 : 1 (clock-driven)\n
dC125/dt = Gd125*C125 + Gr124*C124 - Ga125*C125 : 1 (clock-driven)\n
dC126/dt = Gd126*C126 + Gr125*C125 - Ga126*C126 : 1 (clock-driven)\n
dC127/dt = Gd127*C127 + Gr126*C126 - Ga127*C127 : 1 (clock-driven)\n
dC128/dt = Gd128*C128 + Gr127*C127 - Ga128*C128 : 1 (clock-driven)\n
dC129/dt = Gd129*C129 + Gr128*C128 - Ga129*C129 : 1 (clock-driven)\n
dC130/dt = Gd130*C130 + Gr129*C129 - Ga130*C130 : 1 (clock-driven)\n
dC131/dt = Gd131*C131 + Gr130*C130 - Ga131*C131 : 1 (clock-driven)\n
dC132/dt = Gd132*C132 + Gr131*C131 - Ga132*C132 : 1 (clock-driven)\n
dC133/dt = Gd133*C133 + Gr132*C132 - Ga133*C133 : 1 (clock-driven)\n
dC134/dt = Gd134*C134 + Gr133*C133 - Ga134*C134 : 1 (clock-driven)\n
dC135/dt = Gd135*C135 + Gr134*C134 - Ga135*C135 : 1 (clock-driven)\n
dC136/dt = Gd136*C136 + Gr135*C135 - Ga136*C136 : 1 (clock-driven)\n
dC137/dt = Gd137*C137 + Gr136*C136 - Ga137*C137 : 1 (clock-driven)\n
dC138/dt = Gd138*C138 + Gr137*C137 - Ga138*C138 : 1 (clock-driven)\n
dC139/dt = Gd139*C139 + Gr138*C138 - Ga139*C139 : 1 (clock-driven)\n
dC140/dt = Gd140*C140 + Gr139*C139 - Ga140*C140 : 1 (clock-driven)\n
dC141/dt = Gd141*C141 + Gr140*C140 - Ga141*C141 : 1 (clock-driven)\n
dC142/dt = Gd142*C142 + Gr141*C141 - Ga142*C142 : 1 (clock-driven)\n
dC143/dt = Gd143*C143 + Gr142*C142 - Ga143*C143 : 1 (clock-driven)\n
dC144/dt = Gd144*C144 + Gr143*C143 - Ga144*C144 : 1 (clock-driven)\n
dC145/dt = Gd145*C145 + Gr144*C144 - Ga145*C145 : 1 (clock-driven)\n
dC146/dt = Gd146*C146 + Gr145*C145 - Ga146*C146 : 1 (clock-driven)\n
dC147/dt = Gd147*C147 + Gr146*C146 - Ga147*C147 : 1 (clock-driven)\n
dC148/dt = Gd148*C148 + Gr147*C147 - Ga148*C148 : 1 (clock-driven)\n
dC149/dt = Gd149*C149 + Gr148*C148 - Ga149*C149 : 1 (clock-driven)\n
dC150/dt = Gd150*C150 + Gr149*C149 - Ga150*C150 : 1 (clock-driven)\n
dC151/dt = Gd151*C151 + Gr150*C150 - Ga151*C151 : 1 (clock-driven)\n
dC152/dt = Gd152*C152 + Gr151*C151 - Ga152*C152 : 1 (clock-driven)\n
dC153/dt = Gd153*C153 + Gr152*C152 - Ga153*C153 : 1 (clock-driven)\n
dC154/dt = Gd154*C154 + Gr153*C153 - Ga154*C154 : 1 (clock-driven)\n
dC155/dt = Gd155*C155 + Gr154*C154 - Ga155*C155 : 1 (clock-driven)\n
dC156/dt = Gd156*C156 + Gr155*C155 - Ga156*C156 : 1 (clock-driven)\n
dC157/dt = Gd157*C157 + Gr156*C156 - Ga157*C157 : 1 (clock-driven)\n
dC158/dt = Gd158*C158 + Gr157*C157 - Ga158*C158 : 1 (clock-driven)\n
dC159/dt = Gd159*C159 + Gr158*C158 - Ga159*C159 : 1 (clock-driven)\n
dC160/dt = Gd160*C160 + Gr159*C159 - Ga160*C160 : 1 (clock-driven)\n
dC161/dt = Gd161*C161 + Gr160*C160 - Ga161*C161 : 1 (clock-driven)\n
dC162/dt = Gd162*C162 + Gr161*C161 - Ga162*C162 : 1 (clock-driven)\n
dC163/dt = Gd163*C163 + Gr162*C162 - Ga163*C163 : 1 (clock-driven)\n
dC164/dt = Gd164*C164 + Gr163*C163 - Ga164*C164 : 1 (clock-driven)\n
dC165/dt = Gd165*C165 + Gr164*C164 - Ga165*C165 : 1 (clock-driven)\n
dC166/dt = Gd166*C166 + Gr165*C165 - Ga166*C166 : 1 (clock-driven)\n
dC167/dt = Gd167*C167 + Gr166*C166 - Ga167*C167 : 1 (clock-driven)\n
dC168/dt = Gd168*C168 + Gr167*C167 - Ga168*C168 : 1 (clock-driven)\n
dC169/dt = Gd169*C169 + Gr168*C168 - Ga169*C169 : 1 (clock-driven)\n
dC170/dt = Gd170*C170 + Gr169*C169 - Ga170*C170 : 1 (clock-driven)\n
dC171/dt = Gd171*C171 + Gr170*C170 - Ga171*C171 : 1 (clock-driven)\n
dC172/dt = Gd172*C172 + Gr171*C171 - Ga172*C172 : 1 (clock-driven)\n
dC173/dt = Gd173*C173 + Gr172*C172 - Ga173*C173 : 1 (clock-driven)\n
dC174/dt = Gd174*C174 + Gr173*C173 - Ga174*C174 : 1 (clock-driven)\n
dC175/dt = Gd175*C175 + Gr174*C174 - Ga175*C175 : 1 (clock-driven)\n
dC176/dt = Gd176*C176 + Gr175*C175 - Ga176*C176 : 1 (clock-driven)\n
dC177/dt = Gd177*C177 + Gr176*C176 - Ga177*C177 : 1 (clock-driven)\n
dC178/dt = Gd178*C178 + Gr177*C177 - Ga178*C178 : 1 (clock-driven)\n
dC179/dt = Gd179*C179 + Gr178*C178 - Ga179*C179 : 1 (clock-driven)\n
dC180/dt = Gd180*C180 + Gr179*C179 - Ga180*C180 : 1 (clock-driven)\n
dC181/dt = Gd181*C181 + Gr180*C180 - Ga181*C181 : 1 (clock-driven)\n
dC182/dt = Gd182*C182 + Gr181*C181 - Ga182*C182 : 1 (clock-driven)\n
dC183/dt = Gd183*C183 + Gr182*C182 - Ga183*C183 : 1 (clock-driven)\n
dC184/dt = Gd184*C184 + Gr183*C183 - Ga184*C184 : 1 (clock-driven)\n
dC185/dt = Gd185*C185 + Gr184*C184 - Ga185*C185 : 1 (clock-driven)\n
dC186/dt = Gd186*C186 + Gr185*C185 - Ga186*C186 : 1 (clock-driven)\n
dC187/dt = Gd187*C187 + Gr186*C186 - Ga187*C187 : 1 (clock-driven)\n
dC188/dt = Gd188*C188 + Gr187*C187 - Ga188*C188 : 1 (clock-driven)\n
dC189/dt = Gd189*C189 + Gr188*C188 - Ga189*C189 : 1 (clock-driven)\n
dC190/dt = Gd190*C190 + Gr189*C189 - Ga190*C190 : 1 (clock-driven)\n
dC191/dt = Gd191*C191 + Gr190*C190 - Ga191*C191 : 1 (clock-driven)\n
dC192/dt = Gd192*C192 + Gr191*C191 - Ga192*C192 : 1 (clock-driven)\n
dC193/dt = Gd193*C193 + Gr192*C192 - Ga193*C193 : 1 (clock-driven)\n
dC194/dt = Gd194*C194 + Gr193*C193 - Ga194*C194 : 1 (clock-driven)\n
dC195/dt = Gd195*C195 + Gr194*C194 - Ga195*C195 : 1 (clock-driven)\n
dC196/dt = Gd196*C196 + Gr195*C195 - Ga196*C196 : 1 (clock-driven)\n
dC197/dt = Gd197*C197 + Gr196*C196 - Ga197*C197 : 1 (clock-driven)\n
dC198/dt = Gd198*C198 + Gr197*C197 - Ga198*C198 : 1 (clock-driven)\n
dC199/dt = Gd199*C199 + Gr198*C198 - Ga199*C199 : 1 (clock-driven)\n
dC200/dt = Gd200*C200 + Gr199*C199 - Ga200*C200 : 1 (clock-driven)\n
dC201/dt = Gd201*C201 + Gr200*C200 - Ga201*C201 : 1 (clock-driven)\n
dC202/dt = Gd202*C202 + Gr201*C201 - Ga202*C202 : 1 (clock-driven)\n
dC203/dt = Gd203*C203 + Gr202*C202 - Ga203*C203 : 1 (clock-driven)\n
dC204/dt = Gd204*C204 + Gr203*C203 - Ga204*C204 : 1 (clock-driven)\n
dC205/dt = Gd205*C205 + Gr204*C204 - Ga205*C205 : 1 (clock-driven)\n
dC206/dt = Gd206*C206 + Gr205*C205 - Ga206*C206 : 1 (clock-driven)\n
dC207/dt = Gd207*C207 + Gr206*C206 - Ga207*C207 : 1 (clock-driven)\n
dC208/dt = Gd208*C208 + Gr207*C207 - Ga208*C208 : 1 (clock-driven)\n
dC209/dt = Gd209*C209 + Gr208*C208 - Ga209*C209 : 1 (clock-driven)\n
dC210/dt = Gd210*C210 + Gr209*C209 - Ga210*C210 : 1 (clock-driven)\n
dC211/dt = Gd211*C211 + Gr210*C210 - Ga211*C211 : 1 (clock-driven)\n
dC212/dt = Gd212*C212 + Gr211*C211 - Ga212*C212 : 1 (clock-driven)\n
dC213/dt = Gd213*C213 + Gr212*C212 - Ga213*C213 : 1 (clock-driven)\n
dC214/dt = Gd214*C214 + Gr213*C213 - Ga214*C214 : 1 (clock-driven)\n
dC215/dt = Gd215*C215 + Gr214*C214 - Ga215*C215 : 1 (clock-driven)\n
dC216/dt = Gd216*C216 + Gr215*C215 - Ga216*C216 : 1 (clock-driven)\n
dC217/dt = Gd217*C217 + Gr216*C216 - Ga217*C217 : 1 (clock-driven)\n
dC218/dt = Gd218*C218 + Gr217*C217 - Ga218*C218 : 1 (clock-driven)\n
dC219/dt = Gd219*C219 + Gr218*C218 - Ga219*C219 : 1 (clock-driven)\n
dC220/dt = Gd220*C220 + Gr219*C219 - Ga220*C220 : 1 (clock-driven)\n
dC221/dt = Gd221*C221 + Gr220*C220 - Ga221*C221 : 1 (clock-driven)\n
dC222/dt = Gd222*C222 + Gr221*C221 - Ga222*C222 : 1 (clock-driven)\n
dC223/dt = Gd223*C223 + Gr222*C222 - Ga223*C223 : 1 (clock-driven)\n
dC224/dt = Gd224*C224 + Gr223*C223 - Ga224*C224 : 1 (clock-driven)\n
dC225/dt = Gd225*C225 + Gr224*C224 - Ga225*C225 : 1 (clock-driven)\n
dC226/dt = Gd226*C226 + Gr225*C225 - Ga226*C226 : 1 (clock-driven)\n
dC227/dt = Gd227*C227 + Gr226*C226 - Ga227*C227 : 1 (clock-driven)\n
dC228/dt = Gd228*C228 + Gr227*C227 - Ga228*C228 : 1 (clock-driven)\n
dC229/dt = Gd229*C229 + Gr228*C228 - Ga229*C229 : 1 (clock-driven)\n
dC230/dt = Gd230*C230 + Gr229*C229 - Ga230*C230 : 1 (clock-driven)\n
dC231/dt = Gd231*C231 + Gr230*C230 - Ga231*C231 : 1 (clock-driven)\n
dC232/dt = Gd232*C232 + Gr231*C231 - Ga232*C232 : 1 (clock-driven)\n
dC233/dt = Gd233*C233 + Gr232*C232 - Ga233*C233 : 1 (clock-driven)\n
dC234/dt = Gd234*C234 + Gr233*C233 - Ga234*C234 : 1 (clock-driven)\n
dC235/dt = Gd235*C235 + Gr234*C234 - Ga235*C235 : 1 (clock-driven)\n
dC236/dt = Gd236*C236 + Gr235*C235 - Ga236*C236 : 1 (clock-driven)\n
dC237/dt = Gd237*C237 + Gr236*C236 - Ga237*C237 : 1 (clock-driven)\n
dC238/dt = Gd238*C238 + Gr237*C237 - Ga238*C238 : 1 (clock-driven)\n
dC239/dt = Gd239*C239 + Gr238*C238 - Ga239*C239 : 1 (clock-driven)\n
dC240/dt = Gd240*C240 + Gr239*C239 - Ga240*C240 : 1 (clock-driven)\n
dC241/dt = Gd241*C241 + Gr240*C240 - Ga241*C241 : 1 (clock-driven)\n
dC242/dt = Gd242*C242 + Gr241*C241 - Ga242*C242 : 1 (clock-driven)\n
dC243/dt = Gd243*C243 + Gr242*C242 - Ga243*C243 : 1 (clock-driven)\n
dC244/dt = Gd244*C244 + Gr243*C243 - Ga244*C244 : 1 (clock-driven)\n
dC245/dt = Gd245*C245 + Gr244*C244 - Ga245*C245 : 1 (clock-driven)\n
dC246/dt = Gd246*C246 + Gr245*C245 - Ga246*C246 : 1 (clock-driven)\n
dC247/dt = Gd247*C247 + Gr246*C246 - Ga247*C247 : 1 (clock-driven)\n
dC248/dt = Gd248*C248 + Gr247*C247 - Ga248*C248 : 1 (clock-driven)\n
dC249/dt = Gd249*C249 + Gr248*C248 - Ga249*C249 : 1 (clock-driven)\n
dC250/dt = Gd250*C250 + Gr249*C249 - Ga250*C250 : 1 (clock-driven)\n
dC251/dt = Gd251*C251 + Gr250*C250 - Ga251*C251 : 1 (clock-driven)\n
dC252/dt = Gd252*C252 + Gr251*C251 - Ga252*C252 : 1 (clock-driven)\n
dC253/dt = Gd253*C253 + Gr252*C252 - Ga253*C253 : 1 (clock-driven)\n
dC254/dt = Gd254*C254 + Gr253*C253 - Ga254*C254 : 1 (clock-driven)\n
dC255/dt = Gd255*C255 + Gr254*C254 - Ga255*C255 : 1 (clock-driven)\n
dC256/dt = Gd256*C256 + Gr255*C255 - Ga256*C256 : 1 (clock-driven)\n
dC257/dt = Gd257*C257 + Gr256*C256 - Ga257*C257 : 1 (clock-driven)\n
dC258/dt = Gd258*C258 + Gr257*C257 - Ga258*C258 : 1 (clock-driven)\n
dC259/dt = Gd259*C259 + Gr258*C258 - Ga259*C259 : 1 (clock-driven)\n
dC260/dt = Gd260*C260 + Gr259*C259 - Ga260*C260 : 1 (clock-driven)\n
dC261/dt = Gd261*C261 + Gr260*C260 - Ga261*C261 : 1 (clock-driven)\n
dC262/dt = Gd262*C262 + Gr261*C261 - Ga262*C262 : 1 (clock-driven)\n
dC263/dt = Gd263*C263 + Gr262*C262 - Ga263*C263 : 1 (clock-driven)\n
dC264/dt = Gd264*C264 + Gr263*C263 - Ga264*C264 : 1 (clock-driven)\n
dC265/dt = Gd265*C265 + Gr264*C264 - Ga265*C265 : 1 (clock-driven)\n
dC266/dt = Gd266*C266 + Gr265*C265 - Ga266*C266 : 1 (clock-driven)\n
dC267/dt = Gd267*C267 + Gr266*C266 - Ga267*C267 : 1 (clock-driven)\n
dC268/dt = Gd268*C268 + Gr267*C267 - Ga268*C268 : 1 (clock-driven)\n
dC269/dt = Gd269*C269 + Gr268*C268 - Ga269*C269 : 1 (clock-driven)\n
dC270/dt = Gd270*C270 + Gr269*C269 - Ga270*C270 : 1 (clock-driven)\n
dC271/dt = Gd271*C271 + Gr2
```


(continued from previous page)

Run simulation and plot results

```
sim.reset()
sim.run(200*ms)
```

```
INFO      No numerical integration method specified for group 'neurongroup', using
↳method 'euler' (took 0.15s, trying other methods took 0.00s). [brian2.stateupdaters.
↳base.method_choice]
INFO      No numerical integration method specified for group 'syn_GtACR2_neurongroup',
↳using method 'euler' (took 0.02s, trying other methods took 0.06s). [brian2.
↳stateupdaters.base.method_choice]
INFO      No numerical integration method specified for group 'syn_VfChrimson_
↳neurongroup', using method 'euler' (took 0.01s, trying other methods took 0.01s).
↳[brian2.stateupdaters.base.method_choice]
```

```
fig, (ax1, ax2) = plt.subplots(2, 1, sharex=True)

ax1.plot(spike_mon.t / ms, spike_mon.i[:, :], ",")
ax1.set(ylabel="neuron index", title="spiking")

ax2.step(
    blue_fibers.t_ms, blue_fibers.values + np.arange(n_fibers) * 1.3 + 0.1, c=c["473nm"]
)
ax2.step(amber_fibers.t_ms, amber_fibers.values + np.arange(n_fibers) * 1.3, c=c["590nm
↳"])
ax2.set(
    yticks=[],
    ylabel="intensity per channel",
    title="photostimulation",
    xlabel="time (ms)",
)
```


```
[[],
 Text(0, 0.5, 'intensity per channel'),
 Text(0.5, 1.0, 'photostimulation'),
 Text(0.5, 0, 'time (ms)')]
```

```
../_readthedocs/jupyter_execute/multi_opto_15_1.png
```

As you can see, Vf-Chrimson has fast dynamics, enabling high-frequency control. GtACR2, on the other hand, continues acting long after the light is removed due to its slow deactivation kinetics. We can confirm this by plotting the current from each of the opsins in the first segment of neurons:

```
fig, (ax1, ax2) = plt.subplots(2, 1, sharex=True)
ax1.plot(Imon.t / ms, Imon.I_exc.T / namp, lw=0.2, c=c["590nm"])
ax1.set(title="Vf-Chrimson current", ylabel="$I_{exc}$ (nA)")
ax2.plot(Imon.t / ms, Imon.I_inh.T / namp, lw=0.2, c=c["473nm"])
ax2.set(title="GtACR2 current", ylabel="$I_{inh}$ (nA)")
```

```
[Text(0.5, 1.0, 'GtACR2 current'), Text(0, 0.5, '$I_{inh}$ (nA)')]
```

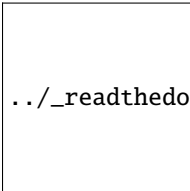

../_readthedocs/jupyter_execute/multi_opto_17_1.png

We can also confirm that blue light has a non-negligible effect on Vf-Chrimson.

Conclusion

As a recap, in this tutorial we've seen how to:

- configure a multi-channel Light device,
- use more than one of them simultaneously, and
- inject multiple opsins of overlapping action spectra into the same network.

6.2.4 On-off control

Here we will see how to set up a minimum, working closed loop with a very simple threshold-triggered control scheme.

Preamble:

```
from brian2 import *
from cleo import *

import matplotlib.pyplot as plt

utilities.style_plots_for_docs()

# the default cython compilation target isn't worth it for
# this trivial example
prefs.codegen.target = "numpy"
```

Set up network

We will use a simple leaky integrate-and-fire network with Poisson spike train input. We use Brian's standard SpikeMonitor to view resulting spikes here for simplicity, but see the electrodes tutorial for a more realistic electrode recording scheme.

```
n = 10
population = NeuronGroup(n, '''
    dv/dt = (-v - 70*mV + Rm*I) / tau : volt
    tau: second
    Rm: ohm
    I: amp''',
    threshold='v>-50*mV',
    reset='v=-70*mV'
)
population.tau = 10*ms
population.Rm = 100*Mohm
population.I = 0*mA
population.v = -70*mV

input_group = PoissonGroup(n, np.linspace(0, 100, n)*Hz + 10*Hz)

S = Synapses(input_group, population, on_pre='v+=5*mV')
S.connect(condition='abs(i-j)<=3')

pop_mon = SpikeMonitor(population)

net = Network([population, input_group, S, pop_mon])

print("Recorded population's equations:")
population.user_equations
```

Recorded population's equations:

$$\frac{dv}{dt} = \frac{IRm - 70mV - v}{\tau} \quad (\text{unit of } v: \text{ V})$$

$$\tau \quad (\text{unit: s})$$

$$Rm \quad (\text{unit: ohm})$$

$$I \quad (\text{unit: A})$$

Run simulation

```
net.run(200*ms)
```


```
INFO      No numerical integration method specified for group 'neurongroup', using
method 'exact' (took 0.15s). [brian2.stateupdaters.base.method_choice]
```

```
sptrains = pop_mon.spike_trains()
fig, ax = plt.subplots()
ax.eventplot([t / ms for t in sptrains.values()], lineoffsets=list(sptrains.keys()))
ax.set(title='population spiking', ylabel='neuron index', xlabel='time (ms)')
```

```

[Text(0.5, 1.0, 'population spiking'),
 Text(0, 0.5, 'neuron index'),
 Text(0.5, 0, 'time (ms)')]

```

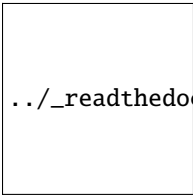

```

../_readthedocs/jupyter_execute/on_off_ctrl_6_1.png

```

Because lower neuron indices receive very little input, we see no spikes for neuron 0. Let's change that with closed-loop control.

IO processor setup

We use the `IOProcessor` class to define interactions with the network. To achieve our goal of making neuron 0 fire, we'll use a contrived, simplistic setup where

1. the recorder reports the voltage of a given neuron (of index 5 in our case),
2. the controller outputs a pulse whenever that voltage is below a certain threshold, and
3. the stimulator applies that pulse to the specified neuron.

So if everything is wired correctly, we'll see bursts of activity in just the first neuron.

```

from cleo.recorders import RateRecorder, VoltageRecorder
from cleo.stimulators import StateVariableSetter

i_rec = int(n / 2)
i_ctrl = 0
sim = CLSimulator(net)
v_rec = VoltageRecorder(name="rec")
sim.inject(v_rec, population[i_rec])
sim.inject(
    StateVariableSetter(name="stim", variable_to_ctrl="I", unit=nA), population[i_ctrl]
)

```

```

CLSimulator(io_processor=None, devices={VoltageRecorder(brian_objects={<StateMonitor,
→ recording ['v'] from 'neurongroup_subgroup'>}, sim=..., name='rec', voltage_var_name='v
→ ', mon=<StateMonitor, recording ['v'] from 'neurongroup_subgroup'>),
→ StateVariableSetter(brian_objects=set(), sim=..., name='stim', value=0, default_
→ value=0, save_history=True, variable_to_ctrl='I', unit=nA, neuron_groups=[<Subgroup
→ 'neurongroup_subgroup_1' of 'neurongroup' from 0 to 1>])})

```

We need to implement the `LatencyIOProcessor` object. For a more sophisticated case we'd use `ProcessingBlock` objects to decompose the computation in the `process` function.

```

from cleo.ioproc import LatencyIOProcessor

trigger_threshold = -60*mV
class ReactivePulseIOProcessor(LatencyIOProcessor):
    def __init__(self, pulse_current=1):

```

(continues on next page)

(continued from previous page)

```

super().__init__(sample_period_ms=1)
self.pulse_current = pulse_current
self.out = {}

def process(self, state_dict, time_ms):
    v = state_dict['rec']
    if v is not None and v < trigger_threshold:
        self.out['stim'] = self.pulse_current
    else:
        self.out['stim'] = 0

    return (self.out, time_ms)

sim.set_io_processor(ReactivePulseIOProcessor(pulse_current=1))

```

```

CLSimulator(io_processor=<__main__.ReactivePulseIOProcessor object at 0x7f57975a1600>,
→ devices={VoltageRecorder(brian_objects={<StateMonitor, recording ['v'] from
→ 'neurongroup_subgroup'>}, sim=..., name='rec', voltage_var_name='v', mon=<StateMonitor,
→ recording ['v'] from 'neurongroup_subgroup'>), StateVariableSetter(brian_
→ objects=set(), sim=..., name='stim', value=0, default_value=0, save_history=True,
→ variable_to_ctrl='I', unit=namp, neuron_groups=[<Subgroup 'neurongroup_subgroup_1' of
→ 'neurongroup' from 0 to 1>])})

```


And run the simulation:

```
sim.run(200*ms)
```

```

fig, (ax1, ax2) = plt.subplots(2, 1, sharex=True)
ax1.plot(pop_mon.t / ms, pop_mon.i[:, "|"])
ax1.plot(
    pop_mon.t[pop_mon.i == i_ctrl] / ms,
    pop_mon.i[pop_mon.i == i_ctrl],
    "|",
    c="#C500CC",
)
ax1.set(title="population spiking", ylabel="neuron index", xlabel="time (ms)")
ax2.fill_between(
    v_rec.mon.t / ms, (v_rec.mon.v.T < trigger_threshold)[:], 0, color=(0.0, 0.72, 1.0)
)
ax2.set(title="pulses", xlabel="time (ms)", ylabel="pulse on/off (1/0)", yticks=[0, 1])
plt.tight_layout()


```

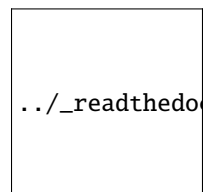

```
../_readthedocs/jupyter_execute/on_off_ctrl_14_0.png
```

Yes, we see the IO processor triggering pulses as expected. And here's a plot of neuron 5's voltage to confirm that those pulses are indeed where we expect them to be, whenever the voltage is below -60 mV.

```
fig, ax = plt.subplots()
ax.set(title=f"Voltage for neuron {i_rec}", ylabel="v (mV)", xlabel='time (ms)')
ax.plot(v_rec.mon.t/ms, v_rec.mon.v.T / mV);
ax.hlines(-60, 0, 400, color='#c500cc');
ax.legend(['v', 'threshold'], loc='upper right');
```

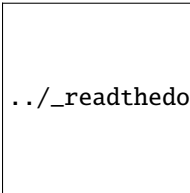

../_readthedocs/jupyter_execute/on_off_ctrl_16_0.png

Conclusion

In this tutorial we've seen the basics of configuring an IOProcessor to implement a closed-loop intervention on a Brian network simulation.

6.2.5 PI control

In this tutorial we'll introduce

1. PI control, a commonly used model-free control method,
2. the concept of decomposing the IOProcessor's computation into ProcessingBlocks, and
3. modeling computation delays on those blocks to reflect hardware and algorithmic speed limitations present in a real experiment.

Preamble:

```
from brian2 import *
import matplotlib.pyplot as plt
from cleo import *

utilities.style_plots_for_docs()

np.random.seed(7000)

# the default cython compilation target isn't worth it for
# this trivial example
prefs.codegen.target = "numpy"
```

Create the Brian network

We'll create a population of 10 LIF neurons mainly driven by feedforward input but with some recurrent connections as well.

```
n = 10
population = NeuronGroup(n, '''
    dv/dt = (-v - 70*mV + Rm*I) / tau : volt
    tau: second
```

(continues on next page)

(continued from previous page)

```

        Rm: ohm
        I: amp'',
        threshold='v>-50*mV',
        reset='v=-70*mV'
    )
    population.tau = 10*ms
    population.Rm = 100*Mohm
    population.I = 0*mA
    population.v = -70*mV

    input_group = PoissonGroup(n, np.linspace(20, 200, n)*Hz)

    S = Synapses(input_group, population, on_pre='v+=5*mV')
    S.connect(condition=f'abs(i-j)<={3}')
    S2 = Synapses(population, population, on_pre='v+=2*mV')
    S2.connect(p=0.2)

    pop_mon = SpikeMonitor(population)

    net = Network(population, input_group, S, S2, pop_mon)
    population.equations

```

$$\frac{dv}{dt} = \frac{IRm - 70mV - v}{\tau} \quad \begin{array}{l} \text{(unit of } v: \text{ V)} \\ \text{(unit: A)} \\ \text{(unit: ohm)} \\ \text{(unit: s)} \end{array}$$

Run simulation without control:

```
net.run(100*ms)
```

```
INFO      No numerical integration method specified for group 'neurongroup', using
↳method 'exact' (took 0.06s). [brian2.stateupdaters.base.method_choice]
```

```

fig, ax = plt.subplots()
ax.scatter(pop_mon.t / ms, pop_mon.i, marker='|', s=200);
ax.set(title='population spiking', ylabel='neuron index', xlabel='time (ms)');

```

```
../_readthedocs/jupyter_execute/PI_ctrl_6_0.png
```

Constructing a closed-loop simulation

We will use the popular model-free PI control to control a single neuron's firing rate. PI stands for proportional-integral, referring to a feedback gain *proportional* to the instantaneous error as well as the *integrated* error over time.

First we construct a CLSimulator from the network:

```
from cleo import CLSimulator
sim = CLSimulator(net)
```

Then, to control neuron i , we need to:

1. capture spiking using a GroundTruthSpikeRecorder

```
from cleo.recorders import GroundTruthSpikeRecorder
i = 0 # neuron to control
rec = GroundTruthSpikeRecorder(name='spike_rec')
sim.inject(rec, population[i])
```

```
CLSimulator(io_processor=None, devices={GroundTruthSpikeRecorder(brian_objects={
    <SpikeMonitor, recording from 'spikemonitor_1'>}, sim=..., name='spike_rec', _mon=
    <SpikeMonitor, recording from 'spikemonitor_1'>, _num_spikes_seen=0, neuron_group=
    <Subgroup 'neurongroup_subgroup' of 'neurongroup' from 0 to 1>))})
```

1. define the firing rate trajectory we want our target neuron to follow

```
# the target firing rate trajectory, as a function of time
def target_Hz(t_ms):
    if t_ms < 250: # constant target at first
        return 400
    else: # sinusoidal afterwards
        a = 200
        t_s = t_ms / 1000
        return a + a * np.sin(2 * np.pi * 20 * t_s)
```

1. estimating its firing rate from incoming spikes using a FiringRateEstimator
2. compute the stimulus intensity with a PIController
3. output that value for a StateVariableSetter stimulator to use

Here we initialize blocks when the IOProcessor is created and define how to process network output and set the control signal in the process function.

```
from cleo.ioproc import (
    LatencyIOProcessor,
    FiringRateEstimator,
    ConstantDelay,
    PIController,
)

class PIRateIOProcessor(LatencyIOProcessor):
    delta = 1 # ms

    def __init__(self):
        super().__init__(sample_period_ms=self.delta, processing="parallel")
```

(continues on next page)

(continued from previous page)

```

self.rate_estimator = FiringRateEstimator(
    tau_ms=15,
    sample_period_ms=self.delta,
    delay=ConstantDelay(4.1), # latency in ms
    save_history=True, # lets us plot later
)

# using hand-tuned gains that seem reasonable
self.pi_controller = PIController(
    target_Hz,
    Kp=0.005,
    Ki=0.04,
    sample_period_ms=self.delta,
    delay=ConstantDelay(2.87), # latency in ms
    save_history=True, # lets us plot later
)

def process(self, state_dict, sample_time_ms):
    spikes = state_dict["spike_rec"]
    # feed output and out_time through each block
    out, time_ms = self.rate_estimator.process(
        spikes, sample_time_ms, sample_time_ms=sample_time_ms
    )
    out, time_ms = self.pi_controller.process(
        out, time_ms, sample_time_ms=sample_time_ms
    )
    # this dictionary output format allows for the flexibility
    # of controlling multiple stimulators
    if out < 0: # limit to positive current
        out = 0
    out_dict = {"I_stim": out}
    # time_ms at the end reflects the delays added by each block
    return out_dict, time_ms

io_processor = PIRateIOProcessor()
sim.set_io_processor(io_processor)

```

```

CLSimulator(io_processor=<__main__.PIRateIOProcessor object at 0x7f6be05fdb70>, devices=
↳ {GroundTruthSpikeRecorder(brian_objects={<SpikeMonitor, recording from 'spikemonitor_1
↳ '>}, sim=..., name='spike_rec', _mon=<SpikeMonitor, recording from 'spikemonitor_1'>, _
↳ num_spikes_seen=0, neuron_group=<Subgroup 'neurongroup_subgroup' of 'neurongroup' from
↳ 0 to 1>}))

```

Note that we can set delays for individual ProcessingBlocks in the IO processor to better approximate the experiment. We use simple constant delays here, but a GaussianDelay class is also available and others could be easily implemented.

Now we inject the stimulator:

```

from cleo.stimulators import StateVariableSetter
sim.inject(
    StateVariableSetter(

```

(continues on next page)

(continued from previous page)

```

        name='I_stim', variable_to_ctrl='I', unit=nA),
        population[i]
    )

```

```

CLSimulator(io_processor=<__main__.PIRateIOProcessor object at 0x7f6be05fdb70>, devices=
    ↳ {StateVariableSetter(brian_objects=set(), sim=..., name='I_stim', value=0, default_
    ↳ value=0, save_history=True, variable_to_ctrl='I', unit=nA, neuron_groups=[<Subgroup
    ↳ 'neurongroup_subgroup_1' of 'neurongroup' from 0 to 1>]),
    ↳ GroundTruthSpikeRecorder(brian_objects={<SpikeMonitor, recording from 'spikemonitor_1'>
    ↳ }, sim=..., name='spike_rec', _mon=<SpikeMonitor, recording from 'spikemonitor_1'>, _
    ↳ num_spikes_seen=0, neuron_group=<Subgroup 'neurongroup_subgroup' of 'neurongroup' from
    ↳ 0 to 1>}))

```

Run the simulation

```
sim.run(300*ms)
```

```

fig, (ax1, ax2, ax3) = plt.subplots(3, 1, sharex=True);
ax1.plot(pop_mon.t / ms, pop_mon.i[:], '|');
ax1.plot(pop_mon.t[pop_mon.i == i]/ms, pop_mon.i[pop_mon.i==i], '|', c='xkcd:hot pink')
ax1.set(title='population spiking', ylabel='neuron index')

ax2.plot(io_processor.rate_estimator.t_in_ms, io_processor.rate_estimator.values, c=
    ↳ 'xkcd:hot pink');
ax2.plot(io_processor.rate_estimator.t_in_ms, [target_Hz(t) for t in io_processor.rate_
    ↳ estimator.t_in_ms],\
        c='xkcd:green');
ax2.set(ylabel='firing rate (Hz)', title=f'neuron {i} activity');
ax2.legend(['estimated firing rate', 'target firing rate']);

ax3.plot(io_processor.pi_controller.t_out_ms, io_processor.pi_controller.values, c=
    ↳ 'xkcd:cerulean')
ax3.set(title='control input', ylabel='$I_{stim}$ (nA)', xlabel='time (ms)')

fig.tight_layout()
fig.show()

```

```

WARNING /tmp/ipykernel_28526/1609733203.py:16: UserWarning: Matplotlib is currently
    ↳ using module://matplotlib_inline.backend_inline, which is a non-GUI backend, so cannot
    ↳ show the figure.
    fig.show()
    [py.warnings]

```

```
../_readthedocs/jupyter_execute/PI_ctrl_19_1.png
```

Note the lag in keeping up with the target firing rate, which can be directly attributed to the ~7 ms delay we coded in to the IO processor.

Conclusion

In this tutorial, we've learned how to

- use PI control to interact with a Brian simulation,
- decompose processing steps into blocks, and
- assign delays to processing blocks to model real-life latency.

6.2.6 LQR optimal control using `ldsctrllest`

This tutorial will be more comprehensive than the others, bringing together all of `cleo`'s main capabilities—electrode recording, optogenetics, and latency modeling—as well as introducing more sophisticated model-based feedback control. To achieve the latter, we will use the `ldsctrllest` Python bindings to the `ldsCtrlEst` C++ library.

Preamble:

```
from brian2 import *
import matplotlib.pyplot as plt
import cleo

cleo.utilities.style_plots_for_docs()

# numpy faster than cython for lightweight example
prefs.codegen.target = 'numpy'
np.random.seed(1856)
```

Network setup

As in the optogenetics tutorial, we'll use a trivial network of a small neuron group biased by Poisson input spikes. We'll use the exponential integrate-and-fire neuron model, which maintains simplicity while modeling an upward membrane potential swing when spiking.

```
n = 2
ng = NeuronGroup(
    n,
    """
    dv/dt = -(v - E_L) + Delta_T*exp((v-theta)/Delta_T) + Rm*I) / tau_m : volt
    I : amp
    """,
    threshold="v>30*mV",
    reset="v=-55*mV",
    namespace={
        "tau_m": 20 * ms,
        "Rm": 500 * Mohm,
        "theta": -50 * mV,
        "Delta_T": 2 * mV,
        "E_L": -70 * mV,
    },
)
```

(continues on next page)

(continued from previous page)

```

)
ng.v = -70 * mV

input_group = PoissonInput(ng, "v", 10, 100 * Hz, 2.5 * mV)

net = Network(ng, input_group)

```

Coordinates, stimulation, and recording

Here we assign coordinates to the neurons and configure the optogenetic intervention and recording setup:

```

from cleo.coords import assign_coords_rand_rect_prism
from cleo.opto import *
from cleo.light import Light, fiber473nm
from cleo.ephys import Probe, SortedSpiking

hor_lim = 0.05
assign_coords_rand_rect_prism(
    ng, xlim=(-hor_lim, hor_lim), ylim=(-hor_lim, hor_lim), zlim=(0.4, 0.6)
)

fiber = Light(
    name="fiber",
    light_model=fiber473nm(),
    coords=(0, 0, 0.3) * mm,
)

opsin = chr2_4s()

spikes = SortedSpiking(
    name="spikes",
    r_perfect_detection=40 * umeter,
    r_half_detection=80 * umeter,
)

probe = Probe(
    coords=[0, 0, 0.5] * mm,
    signals=[spikes],
    save_history=True,
)

cleo.viz.plot(
    ng,
    colors=["xkcd:fuchsia"],
    xlim=(-0.2, 0.2),
    ylim=(-0.2, 0.2),
    zlim=(0.3, 0.8),
    devices=[probe, fiber],
    scatterargs={"alpha": 1},
    axis_scale_unit=mm,
)

```

```
(<Figure size 640x480 with 1 Axes>,
<Axes3D: xlabel='x (mm)', ylabel='y (mm)', zlabel='z (mm)'>)
```

```
../_readthedocs/jupyter_execute/lqr_ctrl_ldsctrllest_5_1.png
```

Looks right. Let's set up the simulation and inject the devices:

```
sim = cleo.CLSimulator(net)
sim.inject(fiber, ng)
sim.inject(opsin, ng, Iopto_var_name='I')
sim.inject(probe, ng)
```

```
CLSimulator(io_processor=None, devices={Probe(sim=..., name='Probe', save_history=True,
→ signals=[SortedSpiking(name='spikes', brian_objects={<SpikeMonitor, recording from
→ 'spikemonitor'>}, probe=..., r_perfect_detection=40. * umetre, r_half_detection=80. *
→ umetre, cutoff_probability=0.01)], probe=NOTHING), FourStateOpsin(sim=..., name='ChR2',
→ save_history=True, on_pre='', spectrum=[(400, 0.34), (422, 0.65), (460, 0.96), (470,
→ 1), (473, 1), (500, 0.57), (520, 0.22), (540, 0.06), (560, 0.01)], required_vars=[(
→ 'Iopto', amp), ('v', volt)], g0=114. * nsiemens, gamma=0.00742, phim=2.33e+23 *
→ (second ** -1) / (meter ** 2), k1=4.15 * khertz, k2=0.868 * khertz, p=0.833, Gf0=37.3
→ * hertz, kf=58.1 * hertz, Gb0=16.1 * hertz, kb=63. * hertz, q=1.94, Gd1=105. * hertz,
→ Gd2=13.8 * hertz, Gr0=0.33 * hertz, E=0. * volt, v0=43. * mvolt, v1=17.1 * mvolt,
→ model="\n          dC1/dt = Gd1*O1 + Gr0*C2 - Ga1*C1 : 1 (clock-driven)\n          dO1/dt
→ = Ga1*C1 + Gb*O2 - (Gd1+Gf)*O1 : 1 (clock-driven)\n          dO2/dt = Ga2*C2 + Gf*O1 -
→ (Gd2+Gb)*O2 : 1 (clock-driven)\n          C2 = 1 - C1 - O1 - O2 : 1\n\n          Theta =
→ int(phi_pre > 0*phi_pre) : 1\n          Hp = Theta * phi_pre**p/(phi_pre**p + phim**p) :
→ 1\n          Ga1 = k1*Hp : hertz\n          Ga2 = k2*Hp : hertz\n          Hq = Theta * phi_
→ pre**q/(phi_pre**q + phim**q) : 1\n          Gf = kf*Hq + Gf0 : hertz\n          Gb =
→ kb*Hq + Gb0 : hertz\n\n          fphi = O1 + gamma*O2 : 1\n          # v1/v0 when v-E == 0
→ via l'Hopital's rule\n          fv = f_unless_x0(\n          (1 - exp(-(V_VAR_NAME_
→ post-E)/v0)) / ((V_VAR_NAME_post-E)/v1),\n          V_VAR_NAME_post - E,\n
→ v1/v0\n          ) : 1\n\n          IOPTO_VAR_NAME_post = -g0*fphi*f_v*(V_VAR_NAME_post-
→ E)*rho_rel : ampere (summed)\n          rho_rel : 1", extra_namespace={'f_unless_x0':
→ <brian2.core.functions.Function object at 0x7f5971fba6e0>}), Light(sim=..., name='fiber
→ ', save_history=True, value=array([0.]), light_model=OpticFiber(R0=100. * umetre,
→ NAfib=0.37, K=125. * metre ** -1, S=7370. * metre ** -1, ntis=1.36), wavelength=0.473
→ * umetre, direction=array([0., 0., 1.]), max_Irr0_mW_per_mm2=None, max_Irr0_mW_per_mm2_
→ viz=None))
```

Prepare controller

Our goal will be to control two neuron's firing rates simultaneously. To do this, we will use the LQR technique explained in Bolus et al., 2021 ("State-space optimal feedback control of optogenetically driven neural activity").

Fit model

Our controller needs a model of the system's dynamics, which we can obtain by fitting to training data. We will generate training data using Gaussian random walk inputs. `ldsCtrlEst` is designed for data coming from an experiment, organized into trials, so we will run the simulation repeatedly, resetting after each run. Here u represents the input and z the spike output.

We will intentionally use very little training data so the importance of adaptive control will become apparent later on.

```
n_trials = 5
n_samp = 100
u = []
z = []
n_u = 1 # 1-dimensional input (just one optogenetic actuator)
n_z = 2 # we'll be controlling two neurons
for trial in range(n_trials):
    # one-sided normally distributed training data, stdev of 10 mW/mm2
    u_trial = 10*np.abs(np.random.randn(n_u, n_samp))
    u.append(u_trial)
    z.append(np.zeros((n_z, n_samp)))
```

The IO processor is simple enough here that we won't bother separating steps using `:class:~cleo.ioproc.ProcessingBlock` objects, which is recommended for more complex scenarios where modularity is more important.

```
from cleo.ioproc import LatencyIOProcessor

class TrainingStimIOP(LatencyIOProcessor):
    i_samp = 0
    i_trial = 0

    # here we just feed in the training inputs and record the outputs
    def process(self, state_dict, sample_time_ms):
        i, t, z_t = state_dict['Probe']['spikes']
        z[self.i_trial][:, self.i_samp] = z_t[:n_z] # just first two neurons
        out = {'fiber': u[self.i_trial][:, self.i_samp]}
        self.i_samp += 1
        return out, sample_time_ms

training_stim_iop = TrainingStimIOP(sample_period_ms=1)
sim.set_io_processor(training_stim_iop)

for i_trial in range(n_trials):
    training_stim_iop.i_trial = i_trial
    training_stim_iop.i_samp = 0
    sim.run(n_samp*ms)
    sim.reset()
```

```
INFO      No numerical integration method specified for group 'neurongroup', using
↳method 'euler' (took 0.02s, trying other methods took 0.07s). [brian2.stateupdaters.
↳base.method_choice]
```

```
INFO      No numerical integration method specified for group 'syn_ChR2_neurongroup',
↳using method 'euler' (took 0.02s, trying other methods took 0.09s). [brian2.
↳stateupdaters.base.method_choice]
```

Now we have u and z in the form we need for `ldsctrlest`'s fitting functions: n_{trial} -length lists of n by n_{samp} arrays. We will now fit Gaussian linear dynamical systems using the SSID algorithm. See [the documentation](#) for more detailed explanations.

```
import ldsctrlest as lds
import ldsctrlest.gaussian as glDs
n_x_fit = 2 # latent dimensionality of system
n_h = 50 # size of block Hankel data matrix
dt = 0.001 # timestep (in seconds)
u_train = lds.UniformMatrixList(u, free_dim=2)
z_train = lds.UniformMatrixList(z, free_dim=2)
ssid = glDs.FitSSID(n_x_fit, n_h, dt, u_train, z_train)
fit, sing_vals = ssid.Run(lds.SSIDWt.kMOESP)
```

Design controller

LQR optimal control

We now use the fit parameters to create the controller system and set additional parameters. The feedback gain, K_c , is especially important, determining how the controller responds to the current “error”—the difference between where the system is (estimated to be) now and where we want it to be. The field of optimal control deals with how to design the controller so as to minimize a cost function reflecting what we care about.

With a linear system (obtained from the fitting procedure above) and quadratic per-timestep cost function L penalizing distance from the reference x^* and the input u

$$L = \frac{1}{2}(x - x^*)^T Q (x - x^*) + \frac{1}{2}u^T R u$$

we can use the closed-form optimal solution called the Linear Quadratic Regulator (LQR).

$$K = (R + B^T P B)^{-1} (B^T P A) \quad u = -Kx$$

The P matrix is obtained by numerically solving the discrete algebraic Riccati equation:

$$P = A^T P A - (A^T P B) (R + B^T P B)^{-1} (B^T P A) + Q$$

```
fit_sys = glDs.System(fit)
# upper and lower bounds on control signal (optic fiber light intensity)
u_lb = 0 # mW/mm2
u_ub = 30 # mW/mm2
controller = glDs.Controller(fit_sys, u_lb, u_ub)
# careful not to use this anymore since controller made a copy
```

(continues on next page)

(continued from previous page)

```

del fit_sys

from scipy.linalg import solve_discrete_are
# cost matrices
# Q reflects how much we care about state error
# we use C'C since we really care about output error, not latent state
Q_cost = controller.sys.C.T @ controller.sys.C
R_cost = 1e-4 * np.eye(n_u) # reflects how much we care about minimizing the stimulus
A, B = controller.sys.A, controller.sys.B
P = solve_discrete_are(A, B, Q_cost, R_cost)
controller.Kc = np.linalg.inv(R_cost + B.T @ P @ B) @ (B.T @ P @ A)
controller.Print()

```

```

***** SYSTEM *****
x:
    0
    0

P:
    1.0000e-06      0
      0    1.0000e-06

A:
    0.5877    0.7122
   -0.7679    0.4800

B:
   -0.0152
   -0.0124

g:
    1.0000

m:
    0
    0

Q:
    0.0783    0.0510
    0.0510    0.1465

Q_m:
    1.0000e-06      0
      0    1.0000e-06

d:
    0
    0

C:
   -0.1074    0.0030
   -0.0810    0.1761

```

(continues on next page)

(continued from previous page)

```

y:
    0
    0

R:
    0.0943    0.0164
    0.0164    0.1574

g_design :    1.0000

u_lb : 0
u_ub : 30

```

We now configure the IOProcessor to use our controller:

```

class CtrlLoop(LatencyIOProcessor):
    def __init__(self, samp_period_ms, controller, y_ref: callable):
        super().__init__(samp_period_ms)
        self.controller = controller
        self.sys = controller.sys
        self.y_ref = y_ref
        self.do_control = False # allows us to turn on and off control

        # for post hoc visualization/analysis:
        self.u = np.empty((n_u, 0))
        self.x_hat = np.empty((n_x_fit, 0))
        self.y_hat = np.empty((n_z, 0))
        self.z = np.empty((n_z, 0))

    def process(self, state_dict, sample_time_ms):
        i, t, z_t = state_dict["Probe"]["spikes"]
        z_t = z_t[:n_z].reshape((-1, 1)) # just first n_z neurons
        self.controller.y_ref = self.y_ref(sample_time_ms)

        u_t = self.controller.ControlOutputReference(z_t, do_control=self.do_control)
        out = {fiber.name: u_t.squeeze()}

        # record variables from this timestep
        self.u = np.hstack([self.u, u_t])
        self.y_hat = np.hstack([self.y_hat, self.sys.y])
        self.x_hat = np.hstack([self.x_hat, self.sys.x])
        self.z = np.hstack([self.z, z_t])

        return out, sample_time_ms + 3 # 3 ms delay

y_ref = 200 * dt # target rate in Hz
ctrl_loop = CtrlLoop(
    samp_period_ms=1, controller=controller, y_ref=lambda t: np.ones((n_z, 1)) * y_ref
)

```

Run the experiment

We'll now run the simulation with and without control to compare.


```
sim.set_io_processor(ctrl_loop)
T0 = 100
sim.run(T0*ms)

ctrl_loop.do_control = True
T1 = 350
sim.run(T1*ms)
```

```
WARNING      'dt' is an internal variable of group 'syn_ChR2_neurongroup', but also exists_
↳ in the run namespace with the value 0.001. The internal variable will be used. [brian2.
↳ groups.group.Group.resolve.resolution_conflict]
```

Now we plot the results to see how well the controller was able to match the desired firing rate:

```
fig, (ax1, ax2, ax3) = plt.subplots(3, 1, sharex=True, figsize=(8,8))
c1 = "#C500CC"
c2 = "#df87e1"
spikes1 = spikes.t_ms[spikes.i == 0]
spikes2 = spikes.t_ms[spikes.i == 1]
ax1.eventplot([spikes1, spikes2], lineoffsets=[1, 2], colors=[c1, c2], lw=0.5)
ax1.set(ylabel='neuron index', ylim=(0.5, 2.5), title='spikes')
ax2.set(ylabel='spikes/s', title='real-time estimated firing rates')
ax2.plot(ctrl_loop.y_hat[0]/dt, c=c1, alpha=0.7, label='$\hat{y}_1$')
ax2.plot(ctrl_loop.y_hat[1]/dt, c=c2, alpha=0.7, label='$\hat{y}_2$')
ax2.hlines(y_ref/dt, 100, T0+T1, color='green', label='reference')
ax2.hlines(0, 0, 100, color='green')
ax2.axvline(T0, c='xkcd:red', linestyle=':', label='ctrl on')
ax2.legend(loc="right")
ax3.plot(range(T0+T1), ctrl_loop.u.T)
ax3.set(xlabel='t (ms)', ylabel='u (mW/mm$^2$)', title='input signal');
```

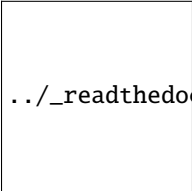

```
../_readthedocs/jupyter_execute/lqr_ctrl_ldsctrllest_21_0.png
```

Looks all right, but in addition to the system's estimated firing rate let's count the spikes over the control period to see how well we hit the target on average:

```
print("Results (spikes/second):")
print('baseline =', np.sum(ctrl_loop.z[:, :T0], axis=1)/(T0/1000))
print('target =', [y_ref*1000, y_ref*1000])
print('lqr achieved =', (np.sum(ctrl_loop.z[:, T0:T0+T1], axis=1)/(T1/1000)).round(1))
```

```
Results (spikes/second):
baseline = [100. 120.]
target = [200.0, 200.0]
lqr achieved = [928.6 748.6]
```

We can see that the system consistently underestimates the true firing rate. And as we could expect, we weren't able to maintain the target firing rate with both neurons simultaneously since one was exposed to more light than the other. However, the controller was able to achieve something. See the appendix for how we can avoid overshooting with both neurons, which should be avoidable.

Conclusion

As a recap, in this tutorial we've seen how to:

- inject optogenetic stimulation into an existing Brian network
- inject an electrode into an existing Brian network to record spikes
- generate training data and fit a Gaussian linear dynamical system to the spiking output using `ldsctrllest`
- configure an `ldsctrllest` LQR controller based on that linear system and design optimal gains
- use that controller in running a complete simulated feedback control experiment

Appendix


Adaptive control

`ldsCtrlEst` also provides an *adaptive* variation on LQR, capable of inferring state beyond our static, linear model and thus able to account for unmodeled disturbances and noise. Let's see how it compares:

```
controller.sys.do_adapt_m = True # enable adaptive disturbance estimation
# set covariance for the disturbance state
# larger values mean the system more readily ascribes changes to unmodeled disturbance
controller.sys.Q_m = 1e-2 * np.eye(n_x_fit)
controller.control_type = lds.kControlTypeAdaptM # enable adaptive control
```

```
ctrl_loop.sys.do_adapt_m = True
T2 = 350
sim.run(T2*ms)
T = T0 + T1 + T2
```

```
fig, (ax1, ax2, ax3) = plt.subplots(3, 1, sharex=True, figsize=(8,8))
spikes1 = spikes.t_ms[spikes.i == 0]
spikes2 = spikes.t_ms[spikes.i == 1]
ax1.eventplot([spikes1, spikes2], lineoffsets=[1, 2], colors=[c1, c2], lw=0.5)
ax1.set(ylabel='neuron index', ylim=(.5, 2.5), title='spikes')
ax2.set(ylabel='spikes/s', title='real-time estimated firing rates')
ax2.plot(ctrl_loop.y_hat[0]/dt, c=c1, alpha=0.7, label='$\hat{y}_1$')
ax2.plot(ctrl_loop.y_hat[1]/dt, c=c2, alpha=0.7, label='$\hat{y}_2$')
ax2.hlines(y_ref/dt, 100, T, color='green', label='reference')
ax2.hlines(0, 0, 100, color='green')
ax2.axvline(T0, c='xkcd:red', linestyle=':', label='ctrl on')
ax2.axvline(T0+T1, c='xkcd:red', linestyle='--', label='adapt on')
ax2.legend(loc="upper left")
ax3.plot(range(T), ctrl_loop.u.T)
ax3.set(xlabel='t (ms)', ylabel='u (mW/mm$^2$)', title='input signal');
```

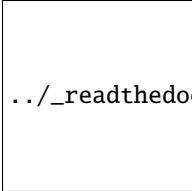

```
../_readthedocs/jupyter_execute/lqr_ctrl_ldsctrllest_29_0.png
```

We can see the effect most easily in the input signal, which has much more variation now. Let's confirm that the firing rates were better balanced around the target:

```
print("Results (spikes/second):")
print('baseline =', np.sum(ctrl_loop.z[:, :T0], axis=1)/(T0/1000))
print('target =', [y_ref*1000, y_ref*1000])
print('static achieved =', (np.sum(ctrl_loop.z[:, T0:T0+T1], axis=1)/(T1/1000)).round(1))
print('adaptive achieved =', (np.sum(ctrl_loop.z[:, T0+T1:T], axis=1)/(T2/1000)).
      ↪round(1))
```


```
Results (spikes/second):
baseline = [100. 120.]
target = [200.0, 200.0]
static achieved = [928.6 748.6]
adaptive achieved = [354.3 317.1]
```

That looks better. Adaptive control achieves a balance between the two neurons, as we would expect.

Post-hoc firing rate estimate

To see if the system's online estimation of firing rates is reasonable, we compute a Gaussian-smoothed version with a 25-ms standard deviation:

```
from scipy.stats import norm
kernel = norm.pdf(np.linspace(-75, 75, 151), scale=25) # 25-ms Gaussian window
smoothed1 = np.convolve(ctrl_loop.z[0, :], kernel, mode='same')
smoothed2 = np.convolve(ctrl_loop.z[1, :], kernel, mode='same')
plt.axhline(y_ref, c='g')
plt.axvline(T0, c='r', ls=':')
plt.axvline(T0+T1, c='r', ls='--')
plt.xlabel("t (ms)")
plt.ylabel("spikes/s")
plt.title("Gaussian-smoothed firing rate")
plt.plot(smoothed1, c=c1)
plt.plot(smoothed2, c=c2);
```

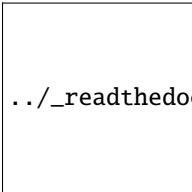

```
../_readthedocs/jupyter_execute/lqr_ctrl_ldsctrllest_34_0.png
```

6.2.7 All-optical control

In this tutorial we'll see how to

- configure and inject a microscope for two-photon imaging and stimulation
- inject a calcium indicator for neurons visible to the microscope
- set up two-photon optogenetics by
 - injecting an opsin with an extended action spectrum and
 - injecting a 2P light source targeting neurons imaged by the microscope

Preamble

```
import brian2 as b2
from brian2 import np
import cleo
from cleo import opto, imaging, light
import matplotlib.pyplot as plt

# for reproducibility
rng = np.random.default_rng(92)
np.random.seed(92)

cleo.utilities.style_plots_for_docs()
```

Brian network setup

All we need are LIF neurons that will spike in response to photostimulation for us to see the resulting fluorescence traces.

```
ng = b2.NeuronGroup(
    200,
    """dv/dt = -(v - E_L) + Rm*Iopto) / tau_m : volt
    Iopto : amp""",
    threshold="v > -50*mV",
    reset="v=E_L",
    namespace={
        "tau_m": 20 * b2.ms,
        "Rm": 500 * b2.Mohm,
        "E_L": -70 * b2.mV,
    },
)
ng.v = -70 * b2.mV
cleo.coords.assign_coords_rand_rect_prism(ng, [-0.2, 0.2], [-0.2, 0.2], [0, 0.4])
sim = cleo.CLSimulator(b2.Network(ng))
```

Microscope configuration

By default a scope selects neurons based on focus depth and assigns them a signal-to-noise ratio (SNR) based on soma size. (Larger cells' SNR decays more slowly with distance from the focal plane.) `focus_depth` and `soma_radius` are taken from scope parameters but can be overridden on injection.

We'll use `GCaMP6f` as our indicator and leave the default parameters, though those can be changed in the function call `gcamp6f()`.

```
scope = imaging.Scope(
    focus_depth=100 * b2.um,
    img_width=500 * b2.um,
    sensor=imaging.gcamp6f(),
)
```

Heterogeneous expression is simulated by providing a `rho_rel_generator` function that takes the number of neurons as an argument and returns a vector of relative expression levels. This is done on scope injection (as opposed to sensor injection afterward) because the scope needs this information to select ROIs with sufficiently high SNR.

```
expr_level_gen = lambda n: rng.lognormal(0, 0.2, n)
sim.inject(scope, ng, rho_rel_generator=expr_level_gen) # uses scope's parameters
# optional overrides (we'll use a small soma to intentionally get few targets)
sim.inject(scope, ng, focus_depth=200 * b2.um, soma_radius=5 * b2.um, rho_rel_
generator=expr_level_gen)
```

```
CLSimulator(io_processor=None, devices={Scope(sim=..., name='Scope', save_history=True,
sensor=GECI(sim=None, name='gcamp6f', save_history=True, model='\n          dCa/dt = -
gamma * (Ca - Ca_rest) / (1 + kappa_S + kappa_B) : mmolar (clock-driven)\n
kappa_B = B_T * K_d / (Ca + K_d)**2 : 1\n\n          CaB_active = Ca_rest + b :
mmolar # add tiny bit to avoid /0\n          db/dt = beta : mmolar (clock-driven)\n
lam = 1/tau_off + 1/tau_on : 1/second\n          kap = 1/tau_off : 1/
second\n          dbeta/dt = (          # should be M/s/s\n
A * (lam - kap) * (Ca - Ca_rest) # M/s/s\n          - (kap + lam) * beta
# M/s/s\n          - kap * lam * b # M/s/s\n          ) : mmolar/second
(clock-driven)\n          \nexc_factor = 1 : 1\n\n          dFF_baseline = 1 / (1
+ (K_d / Ca_rest) ** n_H) : 1\n          dFF = exc_factor * rho_rel * dFF_max * (\n
1 / (1 + (K_d / CaB_active) ** n_H)\n          - dFF_baseline\n
) : 1\n          rho_rel : 1\n          ', on_pre='Ca += dCa_T / (1 + kappa_S
+ kappa_B)', sigma_noise=0.03748181818181818, dFF_1AP=0.097755000000000001, location=
'cytoplasm', cal_model=DynamicCalcium(on_pre='Ca += dCa_T / (1 + kappa_S + kappa_B)',
model='\n          dCa/dt = -gamma * (Ca - Ca_rest) / (1 + kappa_S + kappa_B) :
mmolar (clock-driven)\n          kappa_B = B_T * K_d / (Ca + K_d)**2 : 1', Ca_
rest=50. * nmolar, gamma=292.3 * hertz, B_T=200. * umolar, kappa_S=110, dCa_T=7.6 *
umolar), bind_act_model=DoubExpCalBindingActivation(model='\n          CaB_active =
Ca_rest + b : mmolar # add tiny bit to avoid /0\n          db/dt = beta : mmolar
(clock-driven)\n          lam = 1/tau_off + 1/tau_on : 1/second\n          kap = 1/
tau_off : 1/second\n          dbeta/dt = (          # should be M/s/s\n
A * (lam - kap) * (Ca - Ca_rest) # M/s/s\n          - (kap + lam) *
beta # M/s/s\n          - kap * lam * b # M/s/s\n          ) :
mmolar/second (clock-driven)\n          ', A=7.61251 * khertz, tau_on=1.17164616 *
second, tau_off=10.14020867 * msecond, Ca_rest=50. * nmolar), exc_
model=NullExcitation(model='exc_factor = 1 : 1'), fluor_model='\n          dFF_
baseline = 1 / (1 + (K_d / Ca_rest) ** n_H) : 1\n          dFF = exc_factor * rho_
rel * dFF_max * (\n          1 / (1 + (K_d / CaB_active) ** n_H)\n
- dFF_baseline\n
) : 1\n          rho_rel : 1\n          ', (K_d=290 * nm, n_H=2.7, dFF_max=25.2), img_width=0.5 * mmetre, focus_depth=100. * umetre,
location=array([0., 0., 0.]) * metre, direction=array([0., 0., 1.]), soma_radius=10. *
umetre, snr_cutoff=1)})
```

(continued from previous page)

Signal strength (F/F for 1 spike) per ROI is thus defined as $dFF_{1AP} * \rho_{rel}$.

While Cleo adjusts noise levels according to distance from the focal plane automatically, other adjustments can be made manually by calling `target_neurons_in_plane()` to reflect factors such as cell heterogeneity or indicator expression or increased noise with depth.

Neurons with SNR (defined as $dFF_{1AP} / \sigma_{noise}$) below the scope's `snr_cutoff` parameter are discarded.

```
i_targets, noise_focus_factor, focus_coords = scope.target_neurons_in_plane(ng, focus_
↳ depth=300 * b2.um, soma_radius=5 * b2.um)
# scale noise std randomly to simulate biological variability
std_noise = scope.sensor.sigma_noise * noise_focus_factor * rng.normal(1, .2, len(i_
↳ targets))
sim.inject(scope, ng, focus_depth=None, i_targets=i_targets, sigma_noise=std_noise, rho_
↳ rel_generator=expr_level_gen)

cleo.viz.plot(ng, colors=['#c500cc'], sim=sim)
```

(<Figure size 640x480 with 1 Axes>,
<Axes3D: xlabel='x (um)', ylabel='y (um)', zlabel='z (um)'>)

../_readthedocs/jupyter_execute/all_optical_12_1.png

We can see targets off the visualized plane, resulting from our injections at 200 and 300 m depths, resembling a multi-plane imaging experiment.

When `focus_depth` is set to `None`, corresponding to [sculpted holographic imaging](#), the scope will select all neurons in the volume, or the user can specify a list of neurons to select via `i_targets` on injection.

After all targets are specified, the sensor protein must also be injected. Doing it this way allows for efficiency: the sensor protein is only simulated for the imaged neurons. Expression levels are as sampled previously on scope injection.

```
scope.inject_sensor_for_targets()
```

2p stimulation configuration

Future work will address how to better quantify the relationship between 2P and the 1P action spectra, but for now we'll pretend that Vf-Chrimson activation at 1060 nm is 1/100th that of peak 1P activation for the same power density.

```
opsin = opto.vfchrimson_4s()
opsin.spectrum.append((1060, .01))
sim.inject(opsin, ng)
```


```
CLSSimulator(io_processor=None, devices={GECI(sim=..., name='gcamp6f', save_history=True,
↳ model='n
      dCa/dt = -gamma * (Ca - Ca_rest) / (1 + kappa_S + kappa_B) :
↳ mmolar (clock-driven)\n      kappa_B = B_T * K_d / (Ca + K_d)**2 : 1\n\n
↳ CaB_active = Ca_rest + b : mmolar # add tiny bit to avoid /0\n      db/dt =
↳ beta : mmolar (clock-driven)\n      lam = 1/tau_off + 1/tau_on : 1/second\n
↳ kap = 1/tau_off : 1/second\n      dbeta/dt = (
↳ be M/s/s\n      A * (lam - kap) * (Ca - Ca_rest) # M/s/s\n
6.2 Tutorials (kap + lam) * beta # M/s/s\n      - kap * lam * b # M/s/s\n
↳ ) : mmolar/second (clock-driven)\n      \nexc_factor = 1 : 1\n\n
↳ dFF_baseline = 1 / (1 + (K_d / Ca_rest) ** n_H) : 1\n      dFF = exc_factor *
↳ rho_rel * dFF_max * (\n      1 / (1 + (K_d / CaB_active) ** n_H)\n
```

(continued from previous page)

Now, to target these neurons with 2P laser power, we use `tp_light_from_scope()` to create a `Light` object with a 2P laser *GaussianEllipsoid* light profile centered on each of the scope's targets:

```
laser = light.tp_light_from_scope(scope, wavelength=1060*b2.nmeter, name='laser')
sim.inject(laser, ng)
cleo.viz.plot(ng, colors=['#c500cc'], sim=sim)
```

```
(<Figure size 640x480 with 1 Axes>,
 <Axes3D: xlabel='x (um)', ylabel='y (um)', zlabel='z (um)'>)
```

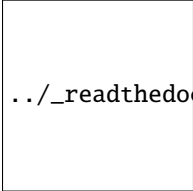

```
../_readthedocs/jupyter_execute/all_optical_19_1.png
```

Simulating an all-optical experiment

We'll stimulate each neuron with a different number of laser pulses and see the resulting calcium traces.

```
from cleo.ioproc import LatencyIOProcessor

# for seeing ground-truth spikes
i_all_targets = scope.i_targets_for_neuron_group(ng)
smon = b2.SpikeMonitor(ng, record=i_all_targets)
sim.network.add(smon)

amplitude_mW = 2.5
pulse_width_ms = 2
interpulse_ms = 20

num_pulses = np.arange(1, scope.n + 1)
train_ends_ms = interpulse_ms * num_pulses

class IOProc(LatencyIOProcessor):
    def process(self, state_dict, time_ms):
        t_in_cycle = time_ms % interpulse_ms
        on = (0 <= t_in_cycle) & (t_in_cycle < pulse_width_ms) & (time_ms < (interpulse_
ms * num_pulses))

        return {"laser": on * amplitude_mW}, time_ms

sim.set_io_processor(IOProc(1))

sim.reset()
tmax = 1000
sim.run(tmax * b2.ms)
```



```

INFO      No numerical integration method specified for group 'neurongroup', using_
↳method 'exact' (took 0.18s). [brian2.stateupdaters.base.method_choice]
INFO      No numerical integration method specified for group 'syn_VfChrimson_
↳neurongroup', using method 'euler' (took 0.02s, trying other methods took 0.06s)._
↳[brian2.stateupdaters.base.method_choice]
INFO      No numerical integration method specified for group 'syn_gcamp6f_neurongroup',
↳ using method 'euler' (took 0.02s, trying other methods took 0.05s). [brian2.
↳stateupdaters.base.method_choice]

```

Now we'll plot the result:

```

import matplotlib.pyplot as plt
sep = .5
n2plot = len(i_all_targets)
fig, (ax1, ax2) = plt.subplots(2, 1, figsize=(8, 6), sharex=True)

spt = np.empty((0,))
spi = np.empty((0,))
for i_target, i_ng in enumerate(i_all_targets):
    single_spt = smon.t[smon.i == i_ng]
    spt = np.concatenate((spt, single_spt))
    spi = np.concatenate((spi, np.full_like(single_spt, i_target)))
ax1.scatter(spt / b2.ms, spi, marker='|', label='spikes')
t_pulse_ms = np.arange(0, num_pulses.max() * interpulse_ms, interpulse_ms)
ax1.scatter(t_pulse_ms, np.full_like(t_pulse_ms, -1), marker='|', color='r', label=
↳'pulses')
ax1.set(ylabel='neuron index', title='ground-truth spikes', xlim=(0, tmax))
ax1.legend()

ax2.plot(scope.t_ms[:,], np.array(scope.dFF)[: , :n2plot] + sep*np.arange(n2plot), lw=.2,
↳c='w')
y_scale = 2
y0 = 0
x0 = -20
ax2.set(xlabel='time (ms)', title='F/F0 per ROI', xlim=(x0, tmax))
ax2.annotate('', xy=(x0, y0), xytext=(x0, y0 + y_scale), arrowprops=dict(arrowstyle='-'))
ax2.annotate(f'{y_scale * 100}% F/F', xy=(x0, y0), va='bottom', ha='right', rotation=90)
ax2.yaxis.set_visible(False)
ax2.spines['left'].set_visible(False)

```

../_readthedocs/jupyter_execute/all_optical_23_0.png

Results are as expected:

- We see one spike for every pulse, though for longer trains we might start to see missed spikes due to the opsin activation plateauing.
- If neurons were closer together, we'd expect to see some off-target effects; if there are any here, they are apparently subthreshold.

- Signal strength is proportional to the number of spikes, as expected, but varies somewhat due to heterogeneous expression.
- We also see that noise levels vary from ROI to ROI due to varying distances from the focal plane.

6.2.8 Video visualization

In this tutorial we'll see how to inject a video visualizer into a simulation.

Preamble:

```
from brian2 import *
import matplotlib.pyplot as plt

from cleo import *

utilities.style_plots_for_docs()

# numpy faster than cython for lightweight example
prefs.codegen.target = 'numpy'
# for reproducibility
np.random.seed(1866)

c_exc = 'xkcd:tomato'
c_inh = 'xkcd:cerulean blue'
```

Set up the simulation

Network

We'll use excitatory and inhibitory populations of *exponential integrate-and-fire neurons*.

```
n_e = 400
n_i = n_e // 4
def eif(n, name):
    ng = NeuronGroup(
        n,
        """
        dv/dt = -(v - E_L) + Delta_T*exp((v-theta)/Delta_T) + Rm*I) / tau_m : volt
        I : amp
        """,
        threshold="v>30*mV",
        reset="v=-55*mV",
        namespace={
            "tau_m": 20 * ms,
            "Rm": 500 * Mohm,
            "theta": -50 * mV,
            "Delta_T": 2 * mV,
            "E_L": -70*mV,
        },
        name=name,
    )
```

(continues on next page)

(continued from previous page)

```

    ng.v = -70 * mV
    return ng

exc = eif(n_e, "exc")
inh = eif(n_i, "inh")
W = 250
p_S = 0.3
S_ei = Synapses(exc, inh, on_pre="v_post+=W*mV/n_e")
S_ei.connect(p=p_S)
S_ie = Synapses(inh, exc, on_pre="v_post-=W*mV/n_i")
S_ie.connect(p=p_S)
S_ee = Synapses(exc, exc, on_pre="v_post+=W*mV/n_e")
S_ee.connect(condition='abs(i-j)<=20')

mon_e = SpikeMonitor(exc)
mon_i = SpikeMonitor(inh)

net = Network(exc, inh, S_ei, S_ie, S_ee, mon_e, mon_i)

```

Coordinates and optogenetics

Here we configure the coordinates and optogenetic stimulation. For more details, see the *“Optogenetic stimulation” tutorial*. Note that we save the arguments used in the plotting function for reuse later on when generating the video.

```

from cleo.coords import assign_coords_uniform_cylinder
from cleo.viz import plot

r = 1
assign_coords_uniform_cylinder(
    exc, xyz_start=(0, 0, 0.3), xyz_end=(0, 0, 0.4), radius=r
)
assign_coords_uniform_cylinder(
    inh, xyz_start=(0, 0, 0.3), xyz_end=(0, 0, 0.4), radius=r
)

from cleo.opto import chr2_4s, fiber473nm, Light

opsin = chr2_4s()
fibers = Light(
    coords=[[0, 0, 0], [700, 0, 0]] * umeter,
    light_model=fiber473nm(),
    max_Irr0_mW_per_mm2=30,
    save_history=True,
)

plotargs = {
    "colors": [c_exc, c_inh],
    "zlim": (0, 600),
    "scatterargs": {"s": 20}, # to adjust neuron marker size
    "axis_scale_unit": umeter,
}

```

(continues on next page)

(continued from previous page)

```
plot(
    exc,
    inh,
    **plotargs,
    devices=[fibers],
)
```

```
(<Figure size 640x480 with 1 Axes>,  
  <Axes3D: xlabel='x (um)', ylabel='y (um)', zlabel='z (um)'>)
```

../_readthedocs/jupyter_execute/video_visualization_5_1.png

Simulator, optogenetics injection

Here we create the simulator and inject the OptogeneticIntervention.

```
sim = CLSimulator(net)
sim.inject(opsin, exc, Iopto_var_name='I')
sim.inject(fibers, exc)
```

```

CLSimulator(io_processor=None, devices={Light(brian_objects=set(), sim=..., name='Light',
→ value=array([0., 0.]), save_history=True, light_model=OpticFiber(R0=100. * umetre,
→ NAfib=0.37, wavelength=0.473 * umetre, K=125. * metre ** -1, S=7370. * metre ** -1,
→ ntis=1.36), coords=array([[0. , 0. , 0. ],
→ [0.7, 0. , 0. ]]) * mmetre, direction=array([0., 0., 1.]), max_Irr0_mW_per_mm2=30,
→ max_Irr0_mW_per_mm2_viz=None, default_value=array([0., 0.])), FourStateOpsin(brian_
→ objects={NeuronGroup(clock=Clock(dt=100. * usecond, name='defaultclock'), when=start,
→ order=0, name='light_agg_ChR2_exc'), Synapses(clock=Clock(dt=100. * usecond, name=
→ 'defaultclock'), when=start, order=0, name='opto_syn_ChR2_exc')}, sim=..., name='ChR2',
→ action_spectrum=[(400, 0.34), (422, 0.65), (460, 0.96), (470, 1), (473, 1), (500, 0.
→ 57), (520, 0.22), (540, 0.06), (560, 0.01)], required_vars=[('Iopto', amp), ('v',
→ volt)], g0=114. * nsiemens, gamma=0.00742, phim=2.33e+23 * (second ** -1) / (meter **
→ 2), k1=4.15 * khertz, k2=0.868 * khertz, p=0.833, Gf0=37.3 * hertz, kf=58.1 * hertz,
→ Gb0=16.1 * hertz, kb=63. * hertz, q=1.94, Gd1=105. * hertz, Gd2=13.8 * hertz, Gr0=0.33
→ * hertz, E=0. * volt, v0=43. * mvolt, v1=17.1 * mvolt, model="\n          dC1/dt =
→ Gd1*O1 + Gr0*C2 - Ga1*C1 : 1 (clock-driven)\n          dO1/dt = Ga1*C1 + Gb*O2 -
→ (Gd1+Gf)*O1 : 1 (clock-driven)\n          dO2/dt = Ga2*C2 + Gf*O1 - (Gd2+Gb)*O2 : 1
→ (clock-driven)\n          C2 = 1 - C1 - O1 - O2 : 1\n\n          Theta = int(phi_pre >
→ 0*phi_pre) : 1\n          Hp = Theta * phi_pre**p/(phi_pre**p + phim**p) : 1\n
→ Ga1 = k1*Hp : hertz\n          Ga2 = k2*Hp : hertz\n          Hq = Theta * phi_pre**q/(phi_
→ pre**q + phim**q) : 1\n          Gf = kf*Hq + Gf0 : hertz\n          Gb = kb*Hq + Gb0 :
→ hertz\n\n          fphi = O1 + gamma*O2 : 1\n          # TODO: get this voltage dependence
→ right \n          # v1/v0 when v-E == 0 via l'Hopital's rule\n          # fv = (1 - exp(-
→ (V_VAR_NAME_post-E)/v0)) / -2 : 1\n          fv = f_unless_x0(\n          (1 - exp(-(V_
→ VAR_NAME_post-E)/v0)) / ((V_VAR_NAME_post-E)/v1),\n          V_VAR_NAME_post - E,\n
→ v1/v0\n          ) : 1\n\n          IOPTO_VAR_NAME_post = -g0*fphi*fV(V_VAR_
→ NAME_post-E)*rho_rel : ampere (summed)\n          rho_rel : 1", extra_namespace={'f_
→ unless_x0': <brian2.core.functions.Function object at 0x7fd1dc9ab4f0>}}))

```

(continued from previous page)

Processor

And we set up open-loop optogenetic stimulation:

```
from cleo.ioproc import LatencyIOProcessor
```

```
fibers.update([5, 0])
```

```
class OpenLoopOpto(LatencyIOProcessor):
```

```
    def process(self, state_dict, time_ms):
```

```
        # random walk stimulation
```

```
        fiber_intensities = fibers.value + np.random.randn(2)*.5
```

```
        fiber_intensities[fiber_intensities < 0] = 0
```

```
        if time_ms > 50:
```

```
            fiber_intensities = [0, 5]
```

```
        return ({"Light": fiber_intensities}, time_ms)
```

```
sim.set_io_processor(OpenLoopOpto(sample_period_ms=1))
```

```
CLSimulator(io_processor=<__main__.OpenLoopOpto object at 0x7fd1dc4be470>, devices=
→ {Light(brian_objects=set(), sim=..., name='Light', value=array([5, 0])), save_
→ history=True, light_model=OpticFiber(R0=100. * umetre, NAfib=0.37, wavelength=0.473 *
→ umetre, K=125. * metre ** -1, S=7370. * metre ** -1, ntis=1.36), coords=array([[0. , 0.
→ , 0. ],
→ [0.7, 0. , 0. ]]) * mmetre, direction=array([0., 0., 1.]), max_Irr0_mW_per_mm2=30,
→ max_Irr0_mW_per_mm2_viz=None, default_value=array([0., 0.])), FourStateOpsin(brian_
→ objects={NeuronGroup(clock=Clock(dt=100. * usecond, name='defaultclock'), when=start,
→ order=0, name='light_agg_ChR2_exc'), Synapses(clock=Clock(dt=100. * usecond, name=
→ 'defaultclock'), when=start, order=0, name='opto_syn_ChR2_exc')}, sim=..., name='ChR2',
→ action_spectrum=[(400, 0.34), (422, 0.65), (460, 0.96), (470, 1), (473, 1), (500, 0.
→ 57), (520, 0.22), (540, 0.06), (560, 0.01)], required_vars=[('Iopto', amp), ('v',
→ volt)], g0=114. * nsiemens, gamma=0.00742, phim=2.33e+23 * (second ** -1) / (meter **
→ 2), k1=4.15 * khertz, k2=0.868 * khertz, p=0.833, Gf0=37.3 * hertz, kf=58.1 * hertz,
→ Gb0=16.1 * hertz, kb=63. * hertz, q=1.94, Gd1=105. * hertz, Gd2=13.8 * hertz, Gr0=0.33
→ * hertz, E=0. * volt, v0=43. * mvolt, v1=17.1 * mvolt, model="\n          dC1/dt =
→ Gd1*O1 + Gr0*C2 - Ga1*C1 : 1 (clock-driven)\n          dO1/dt = Ga1*C1 + Gb*O2 -
→ (Gd1+Gf)*O1 : 1 (clock-driven)\n          dO2/dt = Ga2*C2 + Gf*O1 - (Gd2+Gb)*O2 : 1
→ (clock-driven)\n          C2 = 1 - C1 - O1 - O2 : 1\n\n          Theta = int(phi_pre >
→ 0*phi_pre) : 1\n          Hp = Theta * phi_pre**p/(phi_pre**p + phim**p) : 1\n
→ Ga1 = k1*Hp : hertz\n          Ga2 = k2*Hp : hertz\n          Hq = Theta * phi_pre**q/(phi_
→ pre**q + phim**q) : 1\n          Gf = kf*Hq + Gf0 : hertz\n          Gb = kb*Hq + Gb0 :
→ hertz\n\n          fphi = O1 + gamma*O2 : 1\n          # TODO: get this voltage dependence
→ right \n          # v1/v0 when v-E == 0 via l'Hopital's rule\n          # fv = (1 - exp(-
→ (V_VAR_NAME_post-E)/v0)) / -2 : 1\n          fv = f_unless_x0(\n          (1 - exp(-(V_
→ VAR_NAME_post-E)/v0)) / ((V_VAR_NAME_post-E)/v1),\n          V_VAR_NAME_post - E,\n
→ v1/v0\n          ) : 1\n\n          IOPTO_VAR_NAME_post = -g0*fphi*f_v*(V_VAR_
→ NAME_post-E)*rho_rel : ampere (summed)\n          rho_rel : 1", extra_namespace={'f_
→ unless_x0': <brian2.core.functions.Function object at 0x7fd1dc9ab4f0>}}))
```

Inject VideoVisualizer

A VideoVisualizer is an InterfaceDevice like recorders and stimulators and needs to be injected in order to properly interact with the Brian network. Keep in mind the following:

- It must be injected *after* all other devices for the devices='all' argument to work as expected.
- Similarly to recording and stimulation, you must specify the target neuron groups (to display, in this case) on injection
- The dt argument makes a huge difference on the amount of time it takes to generate the video. You may want to keep this high while experimenting and only lower it when you are ready to generate a high-quality video since the process is so slow.

```
from cleo.viz import VideoVisualizer
```

```
vv = VideoVisualizer(dt=1 * ms, devices="all")
sim.inject(vv, exc, inh)
```

```
CLSimulator(io_processor=<__main__.OpenLoopOpto object at 0x7fd1dc4be470>, devices=
→ {VideoVisualizer(brian_objects={<SpikeMonitor, recording from 'spikemonitor_2'>,
→ <SpikeMonitor, recording from 'spikemonitor_3'>}, sim=..., name='VideoVisualizer',
→ devices=[Light(brian_objects=set(), sim=..., name='Light', value=array([5, 0])), save_
→ history=True, light_model=OpticFiber(R0=100. * umetre, NAfib=0.37, wavelength=0.473 *
→ umetre, K=125. * metre ** -1, S=7370. * metre ** -1, ntis=1.36), coords=array([[0. , 0.
→ , 0. ]],
→ [0.7, 0. , 0. ]]) * mmetre, direction=array([0. , 0. , 1.]), max_Irr0_mW_per_mm2=30,
→ max_Irr0_mW_per_mm2_viz=None, default_value=array([0. , 0. ]]), dt=1. * msecond,
→ fig=None, ax=None, neuron_groups=[NeuronGroup(clock=Clock(dt=100. * usecond, name=
→ 'defaultclock'), when=start, order=0, name='exc'), NeuronGroup(clock=Clock(dt=100. *
→ usecond, name='defaultclock'), when=start, order=0, name='inh')], _spike_mons=[
→ <SpikeMonitor, recording from 'spikemonitor_2'>, <SpikeMonitor, recording from
→ 'spikemonitor_3'>], _num_old_spikes=[0, 0], _value_per_device_per_frame=[], _i_spikes_
→ per_ng_per_frame=[]), Light(brian_objects=set(), sim=..., name='Light', value=array([5,
→ 0])), save_history=True, light_model=OpticFiber(R0=100. * umetre, NAfib=0.37,
→ wavelength=0.473 * umetre, K=125. * metre ** -1, S=7370. * metre ** -1, ntis=1.36),
→ coords=array([[0. , 0. , 0. ]],
→ [0.7, 0. , 0. ]]) * mmetre, direction=array([0. , 0. , 1.]), max_Irr0_mW_per_mm2=30,
→ max_Irr0_mW_per_mm2_viz=None, default_value=array([0. , 0. ]]), FourStateOpsin(brian_
→ objects={NeuronGroup(clock=Clock(dt=100. * usecond, name='defaultclock'), when=start,
→ order=0, name='light_agg_ChR2_exc'), Synapses(clock=Clock(dt=100. * usecond, name=
→ 'defaultclock'), when=start, order=0, name='opto_syn_ChR2_exc')}, sim=..., name='ChR2',
→ action_spectrum=[[400, 0.34), (422, 0.65), (460, 0.96), (470, 1), (473, 1), (500, 0.
→ 57), (520, 0.22), (540, 0.06), (560, 0.01)], required_vars=[('Iopto', amp), ('v',
→ volt)], g0=114. * nsiemens, gamma=0.00742, phim=2.33e+23 * (second ** -1) / (meter **
→ 2), k1=4.15 * khertz, k2=0.868 * khertz, p=0.833, Gf0=37.3 * hertz, kf=58.1 * hertz,
→ Gb0=16.1 * hertz, kb=63. * hertz, q=1.94, Gd1=105. * hertz, Gd2=13.8 * hertz, Gr0=0.33
→ * hertz, E=0. * volt, v0=43. * mvolt, v1=17.1 * mvolt, model="\n
→ dC1/dt =
→ Gd1*O1 + Gr0*C2 - Ga1*C1 : 1 (clock-driven)\n
→ dO1/dt = Ga1*C1 + Gb*O2 -
→ (Gd1+Gf)*O1 : 1 (clock-driven)\n
→ dO2/dt = Ga2*C2 + Gf*O1 - (Gd2+Gb)*O2 : 1
→ (clock-driven)\n
→ C2 = 1 - C1 - O1 - O2 : 1\n\n
→ Theta = int(phi_pre >
→ 0*phi_pre) : 1\n
→ Hp = Theta * phi_pre**p/(phi_pre**p + phim**p) : 1\n
→ Ga1 = k1*Hp : hertz\n
→ Ga2 = k2*Hp : hertz\n
→ Hq = Theta * phi_pre**q/(phi_
→ pre**q + phim**q) : 1\n
→ Gf = kf*Hq + Gf0 : hertz\n
→ Gb = kb*Hq + Gb0 :
→ hertz\n\n
→ fphi = O1 + gamma*O2 : 1\n
→ # TODO: get this voltage dependence.
→ right \n
→ # v1/v0 when v-E == 0 via l'Hopital's rule\n
→ # fv (continues on next page)
→ (V_VAR_NAME_post-E)/v0)) / -2 : 1\n
→ fv = f_unless_x0(\n
→ (1 - exp(-(V_
→ VAR_NAME_post-E)/v0)) / ((V_VAR_NAME_post-E)/v1),\n
→ V_VAR_NAME_post - E\n
→ ) : 1\n\n
→ IOPTO_VAR_NAME_post = -g0*fphi*fv*(V_VAR_
→ NAME_post-E)*rho_rel : ampere (summed)\n
→ rho_rel : 1", extra_namespace={'f_
→ unless_x0': <brian2.core.functions.Function object at 0x7fd1dc9ab4f0>}}))
```

(continued from previous page)

Run simulation and visualize

Here we display a quick plot before generating the video:

```
T = 100
sim.run(T * ms)
```

```
WARNING      'T' is an internal variable of group 'light_prop_ChR2_exc', but also exists
↳ in the run namespace with the value 100. The internal variable will be used. [brian2.
↳ groups.group.Group.resolve.resolution_conflict]
INFO         No numerical integration method specified for group 'exc', using method 'euler
↳ ' (took 0.01s, trying other methods took 0.04s). [brian2.stateupdaters.base.method_
↳ choice]
INFO         No numerical integration method specified for group 'inh', using method 'euler
↳ ' (took 0.01s, trying other methods took 0.02s). [brian2.stateupdaters.base.method_
↳ choice]
```

```
INFO         No numerical integration method specified for group 'opto_syn_ChR2_exc',
↳ using method 'euler' (took 0.02s, trying other methods took 0.08s). [brian2.
↳ stateupdaters.base.method_choice]
```

```
fig, (ax1, ax2, ax3) = plt.subplots(3, 1, sharex=True)
stim_vals = np.array(fibers.values)
sptexc = mon_e.spike_trains()
ax1.eventplot([t/ms for t in sptexc.values()], lineoffsets=list(sptexc.keys()), color=c_
↳ exc)
ax1.set(ylabel="neuron index", title="exc spiking")
sptinh = mon_i.spike_trains()
ax2.eventplot([t/ms for t in sptinh.values()], lineoffsets=list(sptinh.keys()), color=c_
↳ inh)
ax2.set(ylabel="neuron index", title="inh spiking")
ax3.plot(fibers.t_ms, stim_vals[:, 0], c="#72a5f2", lw=2, label='fiber 1')
ax3.plot(fibers.t_ms, stim_vals[:, 1], c="#72c5f2", lw=2, label='fiber 2')
ax3.legend()
ax3.set(ylabel=r"$I_{irr,0}$ (nm/mW$^2$)", title="optogenetic stimulus", xlabel="time (ms)");
```

```
../_readthedocs/jupyter_execute/video_visualization_14_0.png
```


```
fibers.t_ms[0:5]
```

```
[0.0, 0.0, 1.0, 2.0, 3.0]
```

The VideoVisualizer stores the data it needs during the simulation, but hasn't yet produced any visual output. We first use the `generate_Animation()`, plugging in the arguments we used for the original plot.

Also, we set the `max_Irr0_mW_per_mm2_viz` attribute of the optogenetic intervention. This effectively scales how bright the light appears in the visualization. That is, a high maximum irradiance makes the stimulus values small in comparison and produces a faint light, while a low ceiling makes the values relatively large and produces a bright light in the resulting video.

```
fibers.max_Irr0_mW_per_mm2_viz = np.max(stim_vals)
ani = vv.generate_Animation(plotargs, slowdown_factor=10)
```

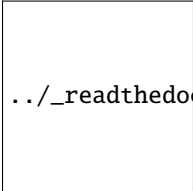

../_readthedocs/jupyter_execute/video_visualization_17_0.png

The `generate_Animation()` function returns a matplotlib `FuncAnimation` object, which you can then use however you want. You will probably want to [save a video](#).

Note that at this point the video still hasn't been rendered; that happens when you try and save or visualize the animation. This step takes a while if your temporal resolution is high, so we suggest you do this only after your experiment is finalized and after you've experimented with low framerate videos to finalize video parameters.

Here we embed the video using HTML so you can see the output:

```
from matplotlib import rc
rc('animation', html='jshtml')

ani
```

```
<matplotlib.animation.FuncAnimation at 0x7fd1d3f176d0>
```

6.2.9 Neo export

Cleo allows you to export data from the whole `CLSimulator` or from individual `InterfaceDevices` objects to `Neo` objects. This facilitates data analysis and visualization with [packages that take Neo as an input format](#), as well as export to various [open](#) and [proprietary](#) ephys data file formats, which could be useful for integration with existing pipelines developed for experimental data.

Setup

```
import brian2 as b2
from brian2 import np
import cleo
from cleo import opto, ephys, light
import neo

# numpy faster than cython for lightweight example
b2.prefs.codegen.target = 'numpy'
# for reproducibility
np.random.seed(33)

cleo.utilities.style_plots_for_docs()
```


We'll set up a basic simulation to see how export works for different devices and the whole simulation.

```

ng = b2.NeuronGroup(
    100,
    """dv/dt = ((-70*mV - v) + (500*Mohm)*Iopto) / (20*ms) : volt
    Iopto : amp""",
    threshold="v > -50*mV",
    reset="v = -70*mV",
)
ng.v = -70 * b2.mV
syn = b2.Synapses(ng, on_pre="v += 3*mV")
cleo.coords.assign_coords_rand_rect_prism(ng, (-.2, .2), (-.2, .2), (0, .4))

light = light.Light(
    coords=[[0.1, 0, 0], [0.1, 0, 0]] * b2.mm, light_model=light.fiber473nm()
)
probe = ephys.Probe(
    coords=[[0, 0, 50], [0, 0, 150], [0, 0, 250]] * b2.um,
    signals=[
        ephys.TKLFPSignal(),
        ephys.MultiUnitSpiking(
            r_perfect_detection=50 * b2.um, r_half_detection=100 * b2.um
        ),
        ephys.SortedSpiking(
            r_perfect_detection=50 * b2.um, r_half_detection=100 * b2.um
        ),
    ],
)
sim = cleo.CLSimulator(b2.Network(ng))

class IOProc(cleo.ioproc.LatencyIOProcessor):
    def process(self, state_dict, t_samp_ms):
        return {'Light': .5 * np.random.rand(light.n)}, t_samp_ms
sim.set_io_processor(IOProc(1))
sim.inject(light, ng).inject(probe, ng, tklf_type="exc")
sim.inject(opto.chr2_4s(), ng)

cleo.viz.plot(ng, colors=['orange'], sim=sim)

```

(<Figure size 640x480 with 1 Axes>,
 <Axes3D: xlabel='x (mm)', ylabel='y (mm)', zlabel='z (mm)'>)

../_readthedocs/jupyter_execute/neo_4_1.png

```
sim.run(50 * b2.ms)
```

INFO No numerical integration method specified for group 'neurongroup', using
 method 'exact' (took 0.23s). [brian2.stateupdaters.base.method_choice]

(continues on next page)

(continued from previous page)

```
INFO      No numerical integration method specified for group 'syn_ChR2_neurongroup',
↳ using method 'euler' (took 0.02s, trying other methods took 0.06s). [brian2.
↳ stateupdaters.base.method_choice]
```

Exporting individual devices

Depending on whether time values are regularly spaced, data is exported to either `AnalogSignal` or `IrregularlySampledSignal` objects, with appropriate annotations and per-channel annotations.

```
light.to_neo()
```

```
AnalogSignal with 2 channels of length 50; units mW/mm**2; datatype float64
name: 'Light'
description: 'Exported from Cleo Light device'
annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 181618)}
sampling rate: 1.0 1/ms
time: 0.0 ms to 50.0 ms
```

```
light.to_neo().array_annotations
```

```
{'x': array([-0.1,  0.1]) * mm,
'y': array([0.,  0.]) * mm,
'z': array([0.,  0.]) * mm,
'direction_x': array([0.,  0.]),
'direction_y': array([0.,  0.]),
'direction_z': array([1.,  1.]),
'i_channel': array([0, 1])}
```

Different signals from the same device are grouped together:

```
probe_neo = probe.to_neo()
probe_neo
```

```
Group with 2 groups, 1 analogsignals
name: 'Probe'
description: 'Exported from Cleo Probe device'
```

```
probe_neo.children
```

```
(AnalogSignal with 3 channels of length 50; units uV; datatype float64
name: 'Probe.TKLFPSignal'
description: 'Exported from Cleo TKLFPSignal object'
annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 239616)}
sampling rate: 1.0 1/ms
time: 0.0 ms to 50.0 ms,
Group with 3 spiketrains
name: 'Probe.MultiUnitSpiking'
description: 'Exported from Cleo MultiUnitSpiking object'
annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 241341)}),
```

(continues on next page)

(continued from previous page)

```

Group with 55 spiketrains
name: 'Probe.SortedSpiking'
description: 'Exported from Cleo SortedSpiking object'
annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 265177)}}

```

Whole-simulation export

Exporting the whole *CLSimulator* object generates a *Block* object with a single *Segment* object containing all the data from the simulation.

```
sim.to_neo()
```

```

Block with 1 segments, 1 groups
description: 'Exported from Cleo simulation'
annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 320150)}
# segments (N=1)
0: Segment with 2 analogsignals, 58 spiketrains
  # analogsignals (N=2)
  0: AnalogSignal with 2 channels of length 50; units mW/mm**2; datatype float64
    name: 'Light'
    description: 'Exported from Cleo Light device'
    annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 320674)}
    sampling rate: 1.0 1/ms
    time: 0.0 ms to 50.0 ms
  1: AnalogSignal with 3 channels of length 50; units uV; datatype float64
    name: 'Probe.TKLFPSignal'
    description: 'Exported from Cleo TKLFPSignal object'
    annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 321221)}
    sampling rate: 1.0 1/ms
    time: 0.0 ms to 50.0 ms

```

Trial structure

Using the convention that a *Segment* represents a single “trial,” there are a couple ways we can generate trial-structured Neo data:

1. Resetting the simulation and exporting the data from each trial individually, merging them into a single *Block* object.
2. Structuring multiple trials into a single call to *run()* and segmenting the data on export. This is not implemented, but could be in the future if there is sufficient demand.

Let’s try the first of these options:

```

block1 = sim.to_neo()
block1.segments[0].name = 'trial 1'
sim.reset()
sim.run(50 * b2.ms)
block2 = sim.to_neo()
block2.segments[0].name = 'trial 2'

```

(continues on next page)

(continued from previous page)

```

block = neo.Block()
block.segments.extend(block1.segments)
block.segments.extend(block2.segments)
block.groups.extend(block1.groups)
block.groups.extend(block2.groups)
block

```

```

Block with 2 segments, 2 groups
# segments (N=2)
0: Segment with 2 analogsignals, 58 spiketrains
  name: 'trial 1'
  # analogsignals (N=2)
  0: AnalogSignal with 2 channels of length 50; units mW/mm**2; datatype float64
    name: 'Light'
    description: 'Exported from Cleo Light device'
    annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 375717)}
    sampling rate: 1.0 1/ms
    time: 0.0 ms to 50.0 ms
  1: AnalogSignal with 3 channels of length 50; units uV; datatype float64
    name: 'Probe.TKLFPSignal'
    description: 'Exported from Cleo TKLFPSignal object'
    annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 33, 376227)}
    sampling rate: 1.0 1/ms
    time: 0.0 ms to 50.0 ms
1: Segment with 2 analogsignals, 63 spiketrains
  name: 'trial 2'
  # analogsignals (N=2)
  0: AnalogSignal with 2 channels of length 50; units mW/mm**2; datatype float64
    name: 'Light'
    description: 'Exported from Cleo Light device'
    annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 34, 138750)}
    sampling rate: 1.0 1/ms
    time: 0.0 ms to 50.0 ms
  1: AnalogSignal with 3 channels of length 50; units uV; datatype float64
    name: 'Probe.TKLFPSignal'
    description: 'Exported from Cleo TKLFPSignal object'
    annotations: {'export_datetime': datetime.datetime(2023, 9, 6, 9, 59, 34, 139149)}
    sampling rate: 1.0 1/ms
    time: 0.0 ms to 50.0 ms

```

The more attractive approach suggested by Neo's API of `block1.merge(block2)` does not work, but could potentially be fixed in the future.

Example analysis use case

Let's try `elephant`, following their statistics tutorial:

```
import matplotlib.pyplot as plt
from elephant.statistics import instantaneous_rate, time_histogram, mean_firing_rate
import quantities as pq

fig, ax = plt.subplots()
st1 = block.segments[0].spiketrains[0]
histogram_rate = time_histogram([st1], 5*pq.ms, output='rate')
inst_rate = instantaneous_rate(st1, sampling_period=1*pq.ms)

# plotting the original spiketrain
ax.plot(
    st1,
    [0] * len(st1),
    "r",
    marker=2,
    ms=25,
    markeredgewidth=2,
    lw=0,
    label="poisson spike times",
)

# mean firing rate
ax.hlines(
    mean_firing_rate(st1),
    xmin=st1.t_start,
    xmax=st1.t_stop,
    linestyle="--",
    label="mean firing rate",
)

# time histogram
ax.bar(
    histogram_rate.times,
    histogram_rate.magnitude.flatten(),
    width=histogram_rate.sampling_period,
    align="edge",
    alpha=0.3,
    label="time histogram (rate)",
)

# instantaneous rate
ax.plot(
    inst_rate.times.rescale(pq.ms),
    inst_rate.rescale(histogram_rate.dimensionality).magnitude.flatten(),
    label="instantaneous rate",
)

# axis labels and legend
ax.set(
```

(continues on next page)


(continued from previous page)

```

xlabel=f"time [{st1.times.dimensionality.latex}]",
ylabel=f"firing rate [{histogram_rate.dimensionality.latex}]",
xlim=(st1.t_start, st1.t_stop),
)
ax.legend()

```

```
<matplotlib.legend.Legend at 0x7f527c923e20>
```

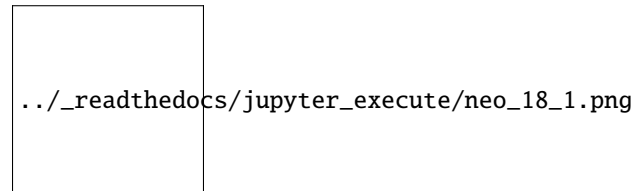

6.3 Reference

6.3.1 cleo module

class `cleo.CLSimulator`(*network*: *brian2.core.network.Network*)

Bases: `cleo.base.NeoExportable`

The centerpiece of cleo. Integrates simulation components and runs.

Method generated by attrs for class CLSimulator.

devices: `set[cleo.base.InterfaceDevice]`

get_state() → dict

Return current recorder measurements.

Returns A dictionary of *name: state* pairs for all recorders in the simulator.

Return type dict

inject(*device*: *cleo.base.InterfaceDevice*, **neuron_groups*: *brian2.groups.neurongroup.NeuronGroup*, ***kwargs*: *Any*) → *cleo.base.CLSimulator*

Inject `InterfaceDevice` into the network, connecting to specified neurons.

Calls `connect_to_neuron_group()` for each group with *kwargs* and adds the device's *brian_objects* to the simulator's *network*.

Parameters **device** (*InterfaceDevice*) – Device to inject

Returns self

Return type *CLSimulator*

io_processor: `cleo.base.IOProcessor`

network: `brian2.core.network.Network`

The Brian network forming the core model

recorders: `dict[str, cleo.base.Recorder]`

reset(**kwargs)

Reset the simulator to a neutral state

Restores the Brian Network to where it was when the CLSimulator was last modified (last injection, IO-Processor change). Calls reset() on devices and IOProcessor.

run(duration: *brian2.units.fundamentalunits.Quantity*, **kwparams) → None

Run simulation.

Parameters

- **duration** (*brian2 temporal Quantity*) – Length of simulation
- ****kwparams** (*additional arguments passed to brian2.run()*) – level has a default value of 1

set_io_processor(io_processor, communication_period=None) → *cleo.base.CLSimulator*

Set simulator IO processor

Will replace any previous IOProcessor so there is only one at a time. A Brian NetworkOperation is created to govern communication between the Network and the IOProcessor.

Parameters **io_processor** (*IOProcessor*) –

Returns self

Return type *CLSimulator*

stimulators: dict[str, *cleo.base.Stimulator*]

to_neo() → *neo.core.block.Block*

Return a neo.core.AnalogSignal object with the device's data

Returns Neo object representing exported data

Return type *neo.core.BaseNeo*

update_stimulators(ctrl_signals) → None

Update stimulators with output from the *IOProcessor*

Parameters **ctrl_signals** (*dict*) – {*stimulator_name: ctrl_signal*} dictionary with values to update each stimulator.

class *cleo.IOProcessor*

Bases: *abc.ABC*

Abstract class for implementing sampling, signal processing and control

This must be implemented by the user with their desired closed-loop use case, though most users will find the *LatencyIOProcessor()* class more useful, since delay handling is already defined.

abstract **get_ctrl_signal**(time) → dict

Get per-stimulator control signal from the *IOProcessor*.

Parameters **time** (*Brian 2 temporal Unit*) – Current timestep

Returns A {'stimulator_name': value} dictionary for updating stimulators.

Return type dict

abstract **is_sampling_now**(time) → bool

Determines whether the processor will take a sample at this timestep.

Parameters **time** (*Brian 2 temporal Unit*) – Current timestep.

Return type bool

abstract put_state(*state_dict: dict, time*) → None

Deliver network state to the [IOProcessor](#).

Parameters

- **state_dict** (*dict*) – A dictionary of recorder measurements, as returned by [get_state\(\)](#)
- **time** (*brian2 temporal Unit*) – The current simulation timestep. Essential for simulating control latency and for time-varying control.

reset(***kwargs*) → None

sample_period_ms: float

Determines how frequently the processor takes samples

class cleo.**InterfaceDevice**(**, name: str = NOTHING, save_history: bool = True*)

Bases: abc.ABC

Base class for devices to be injected into the network

Method generated by attrs for class InterfaceDevice.

add_self_to_plot(*ax: mpl_toolkits.mplot3d.axes3d.Axes3D, axis_scale_unit: brian2.units.fundamentalunits.Unit, **kwargs*) → list[matplotlib.artist.Artist]

Add device to an existing plot

Should only be called by [plot\(\)](#).

Parameters

- **ax** (*Axes3D*) – The existing matplotlib Axes object
- **axis_scale_unit** (*Unit*) – The unit used to label axes and define chart limits
- ****kwargs** (*optional*) –

Returns A list of artists used to render the device. Needed for use in conjunction with [VideoVisualizer](#).

Return type list[Artist]

brian_objects: set

All the Brian objects added to the network by this device. Must be kept up-to-date in [connect_to_neuron_group\(\)](#) and other functions so that those objects can be automatically added to the network when the device is injected.

abstract connect_to_neuron_group(*neuron_group: brian2.groups.neurongroup.NeuronGroup, **kwargs*) → None

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a NeuronGroup, Synapses, or Monitor, make sure to add these to *self.brian_objects*.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwargs** (*optional, passed from inject or*) – *inject*

init_for_simulator(*simulator*: [cleo.base.CLSimulator](#)) → None

Initialize device for simulator on initial injection

This function is called only the first time a device is injected into a simulator and performs any operations that are independent of the individual neuron groups it is connected to.

Parameters **simulator** ([CLSimulator](#)) – simulator being injected into

name: **str**

Unique identifier for device, used in sampling, plotting, etc. Name of the class by default.

reset(***kwargs*) → None

Reset the device to a neutral state

save_history: **bool**

Determines whether times and inputs/outputs are recorded.

For stimulators, this is when [update\(\)](#) is called. For recorders, it is when [get_state\(\)](#) is called.

sim: [cleo.base.CLSimulator](#)

The simulator the device is injected into

update_artists(*artists*: [list\[matplotlib.artist.Artist\]](#), **args*, ***kwargs*) → [list\[matplotlib.artist.Artist\]](#)

Update the artists used to render the device

Used to set the artists' state at every frame of a video visualization. The current state would be passed in **args* or ***kwargs*

Parameters **artists** ([list\[Artist\]](#)) – the artists used to render the device originally, i.e., which were returned from the first [add_self_to_plot\(\)](#) call.

Returns The artists that were actually updated. Needed for efficient blit rendering, where only updated artists are re-rendered.

Return type [list\[Artist\]](#)

class [cleo.Recorder](#)(**name*: [str](#) = [NOTHING](#), *save_history*: [bool](#) = [True](#))

Bases: [cleo.base.InterfaceDevice](#)

Device for taking measurements of the network.

Method generated by attrs for class Recorder.

abstract **get_state**() → Any

Return current measurement.

class [cleo.Stimulator](#)(*default_value*: Any = 0, **name*: [str](#) = [NOTHING](#), *save_history*: [bool](#) = [True](#))

Bases: [cleo.base.InterfaceDevice](#), [cleo.base.NeoExportable](#)

Device for manipulating the network

Method generated by attrs for class Stimulator.

default_value: Any

The default value of the device—used on initialization and on [reset\(\)](#)

reset(***kwargs*) → None

Reset the stimulator device to a neutral state

t_ms: [list\[float\]](#)

Times stimulator was updated, stored if [save_history](#)

to_neo()

Return a `neo.core.AnalogSignal` object with the device's data

Returns Neo object representing exported data

Return type `neo.core.BaseNeo`

update(ctrl_signal) → None

Set the stimulator value.

By default this simply sets *value* to *ctrl_signal*. You will want to implement this method if your stimulator requires additional logic. Use `super.update(self, value)` to preserve the `self.value` attribute logic

Parameters **ctrl_signal** (any) – The value the stimulator is to take.

value: Any

The current value of the stimulator device

values: list[Any]

Values taken by the stimulator at each `update()` call, stored if `save_history`

class `cleo.SynapseDevice`(*extra_namespace: dict = NOTHING, *, name: str = NOTHING, save_history: bool = True*)

Bases: `cleo.base.InterfaceDevice`

Base class for devices that record from/stimulate neurons via a Synapses object with device-specific model. Used for opsin and indicator classes

Method generated by attrs for class SynapseDevice.

connect_to_neuron_group(*neuron_group: brian2.groups.neurongroup.NeuronGroup, **kwargs*) → None

Transfect neuron group with device.

Parameters **neuron_group** (*NeuronGroup*) – The neuron group to transform

Keyword Arguments

- **p_expression** (*float*) – Probability ($0 \leq p \leq 1$) that a given neuron in the group will express the protein. 1 by default.
- **i_targets** (*array-like*) – Indices of neurons in the group to transfect. recommended for efficiency when stimulating or imaging a small subset of the group. Incompatible with `p_expression`.
- **rho_rel** (*float*) – The expression level, relative to the standard model fit, of the protein. 1 by default. For heterogeneous expression, this would have to be modified in the light-dependent synapse post-injection, e.g., `opsin.syns["neuron_group_name"].rho_rel = ...`
- **[default_name]_var_name** (*str*) – See [required_vars](#). Allows for custom variable names.

extra_namespace: dict

Additional items (beyond parameters) to be added to the opto synapse namespace

init_syn_vars(*syn: brian2.synapses.synapses.Synapses*) → None

Initializes appropriate variables in Synapses implementing the model

Can also be used to reset the variables.

Parameters **syn** (*Synapses*) – The synapses object implementing this model

model: str

Basic Brian model equations string.

Should contain a *rho_rel* term reflecting relative expression levels. Will likely also contain special NeuronGroup-dependent symbols such as V_VAR_NAME to be replaced on injection in [modify_model_and_params_for_ng\(\)](#).

modify_model_and_params_for_ng(*neuron_group*: [brian2.groups.neurongroup.NeuronGroup](#),
injet_params: dict) → Tuple[[brian2.equations.equations.Equations](#), dict]

Adapt model for given neuron group on injection

This enables the specification of variable names differently for each neuron group, allowing for custom names and avoiding conflicts.

Parameters

- **neuron_group** (*NeuronGroup*) – NeuronGroup this opsin model is being connected to
- **injet_params** (*dict*) – kwargs passed in on injection, could contain variable names to plug into the model

Keyword Arguments **model** (*str*, *optional*) – Model to start with, by default that defined for the class. This allows for prior string manipulations before it can be parsed as an *Equations* object.

Returns A tuple containing an *Equations* object and a parameter dictionary, constructed from [model](#) and [params](#), respectively, with modified names for use in *synapses*

Return type *Equations*, dict

on_pre: str

Model string for Synapses reacting to spikes.

property params: dict

Returns a dictionary of parameters for the model

per_ng_unit_replacements: list[Tuple[str, str]]

List of (UNIT_NAME, neuron_group_specific_unit_name) tuples to be substituted in the model string on injection and before checking required variables.

required_vars: list[Tuple[str, [brian2.units.fundamentalunits.Unit](#)]]

Default names of state variables required in the neuron group, along with units, e.g., [('Iopto', amp)].

It is assumed that non-default values can be passed in on injection as a keyword argument [default_name]_var_name=[non_default_name] and that these are found in the model string as [DEFAULT_NAME]_VAR_NAME before replacement.

reset(**kwargs)

Reset the device to a neutral state

source_ngs: dict[str, [brian2.groups.neurongroup.NeuronGroup](#)]

{target_ng.name: source_ng} dict of source neuron groups.

The source is the target itself by default or light aggregator neurons for [LightDependent](#).

synapses: dict[str, [brian2.synapses.synapses.Synapses](#)]

Stores the synapse objects implementing the model, connecting from source (light aggregator neurons or the target group itself) to target neuron groups, with NeuronGroup name keys and Synapses values.

6.3.2 cleo.coords module

Contains functions for assigning neuron coordinates and visualizing

`cleo.coords.assign_coords(neuron_group: brian2.groups.neurongroup.NeuronGroup, coords: brian2.units.fundamentalunits.Quantity)`

`cleo.coords.assign_coords_grid_rect_prism(neuron_group: brian2.groups.neurongroup.NeuronGroup, xlim: Tuple[float, float], ylim: Tuple[float, float], zlim: Tuple[float, float], shape: Tuple[int, int, int], unit: brian2.units.fundamentalunits.Unit = mmetre) → None`

Assign grid coordinates to neurons in a rectangular grid

Parameters

- **neuron_group** (*NeuronGroup*) – The neuron group to assign coordinates to
- **xlim** (*Tuple[float, float]*) – xmin, xmax, with no unit
- **ylim** (*Tuple[float, float]*) – ymin, ymax, with no unit
- **zlim** (*Tuple[float, float]*) – zmin, zmax with no unit
- **shape** (*Tuple[int, int, int]*) – n_x, n_y, n_z tuple representing the shape of the resulting grid
- **unit** (*Unit, optional*) – Brian unit determining what scale to use for coordinates, by default mm

Raises **ValueError** – When the shape is incompatible with the number of neurons in the group

`cleo.coords.assign_coords_rand_cylinder(neuron_group: brian2.groups.neurongroup.NeuronGroup, xyz_start: Tuple[float, float, float], xyz_end: Tuple[float, float, float], radius: float, unit: brian2.units.fundamentalunits.Unit = mmetre) → None`

Assign random coordinates within a cylinder.

Parameters

- **neuron_group** (*NeuronGroup*) – neurons to assign coordinates to
- **xyz_start** (*Tuple[float, float, float]*) – starting position of cylinder without unit
- **xyz_end** (*Tuple[float, float, float]*) – ending position of cylinder without unit
- **radius** (*float*) – radius of cylinder without unit
- **unit** (*Unit, optional*) – Brian unit to scale other params, by default mm

`cleo.coords.assign_coords_rand_rect_prism(neuron_group: brian2.groups.neurongroup.NeuronGroup, xlim: Tuple[float, float], ylim: Tuple[float, float], zlim: Tuple[float, float], unit: brian2.units.fundamentalunits.Unit = mmetre) → None`

Assign random coordinates to neurons within a rectangular prism

Parameters

- **neuron_group** (*NeuronGroup*) – neurons to assign coordinates to
- **xlim** (*Tuple[float, float]*) – xmin, xmax without unit
- **ylim** (*Tuple[float, float]*) – ymin, ymax without unit
- **zlim** (*Tuple[float, float]*) – zmin, zmax without unit

- **unit** (*Unit*, *optional*) – Brian unit to specify scale implied in limits, by default mm

```
cleo.coords.assign_coords_uniform_cylinder(neuron_group: brian2.groups.neurongroup.NeuronGroup,
                                           xyz_start: Tuple[float, float, float], xyz_end: Tuple[float,
float, float], radius: float, unit:
                                           brian2.units.fundamentalunits.Unit = mmetre) → None
```

Assign uniformly spaced coordinates within a cylinder.

Parameters

- **neuron_group** (*NeuronGroup*) – neurons to assign coordinates to
- **xyz_start** (*Tuple[float, float, float]*) – starting position of cylinder without unit
- **xyz_end** (*Tuple[float, float, float]*) – ending position of cylinder without unit
- **radius** (*float*) – radius of cylinder without unit
- **unit** (*Unit*, *optional*) – Brian unit to scale other params, by default mm

```
cleo.coords.assign_xyz(neuron_group: brian2.groups.neurongroup.NeuronGroup, x: numpy.ndarray, y:
numpy.ndarray, z: numpy.ndarray, unit: brian2.units.fundamentalunits.Unit =
mmetre)
```

Assign arbitrary coordinates to neuron group.

Parameters

- **neuron_group** (*NeuronGroup*) – neurons to be assigned coordinates
- **x** (*np.ndarray*) – x positions to assign (preferably 1D with no unit)
- **y** (*np.ndarray*) – y positions to assign (preferably 1D with no unit)
- **z** (*np.ndarray*) – z positions to assign (preferably 1D with no unit)
- **unit** (*Unit*, *optional*) – Brian unit determining what scale to use for coordinates, by default mm

```
cleo.coords.concat_coords(*coords: brian2.units.fundamentalunits.Quantity) →
brian2.units.fundamentalunits.Quantity
```

Combine multiple coordinate Quantity arrays into one

Parameters **coords* (*Quantity*) – Multiple coordinate n x 3 Quantity arrays to combine

Returns A single n x 3 combined Quantity array

Return type Quantity

```
cleo.coords.coords_from_ng(ng: brian2.groups.neurongroup.NeuronGroup) →
brian2.units.fundamentalunits.Quantity
```

Get nx3 coordinate array from NeuronGroup.

```
cleo.coords.coords_from_xyz(x: brian2.units.fundamentalunits.Quantity, y:
brian2.units.fundamentalunits.Quantity, z:
brian2.units.fundamentalunits.Quantity) →
brian2.units.fundamentalunits.Quantity
```

Get ...x3 coordinate array from x, y, z arrays (with units).

6.3.3 cleo.ephys module

Contains probes, coordinate convenience functions, signals, spiking, and LFP

```
class cleo.ephys.MultiUnitSpiking(r_perfect_detection: brian2.units.fundamentalunits.Quantity,  
                                r_half_detection: brian2.units.fundamentalunits.Quantity,  
                                cutoff_probability: float = 0.01, *, name: str = NOTHING)
```

Bases: [cleo.ephys.spiking.Spiking](#)

Detects spikes per channel, that is, unsorted.

Method generated by attrs for class MultiUnitSpiking.

```
connect_to_neuron_group(neuron_group: brian2.groups.neurongroup.NeuronGroup, **kwargs) →  
None
```

Configure signal to record from specified neuron group

Parameters *neuron_group* (*NeuronGroup*) – group to record from

```
get_state() → tuple[NDArray[Any, ..., UInt[64]], NDArray[Any, ..., Float[64]], NDArray[Any, ...,  
                    UInt[64]]]
```

Return spikes since method was last called (*i*, *t_ms*, *y*)

Returns (*i*, *t_ms*, *y*) where *i* is channel (for multi-unit) or neuron (for sorted) spike indices, *t_ms* is spike times, and *y* is a spike count vector suitable for control- theoretic uses—i.e., a 0 for every channel/neuron that hasn’t spiked and a 1 for a single spike.

Return type tuple[NDArray[np.uint], NDArray[float], NDArray[np.uint]]

```
to_neo() → neo.core.group.Group
```

Return a neo.core.AnalogSignal object with the device’s data

Returns Neo object representing exported data

Return type neo.core.BaseNeo

```
class cleo.ephys.Probe(coords: brian2.units.fundamentalunits.Quantity, signals:  
                     list[cleo.ephys.probes.Signal] = NOTHING, *, name: str = NOTHING, save_history:  
                     bool = True)
```

Bases: [cleo.base.Recorder](#), [cleo.base.NeoExportable](#)

Picks up specified signals across an array of electrodes.

Visualization kwargs

- **marker** (*str, optional*) – The marker used to represent each contact. “x” by default.
- **size** (*float, optional*) – The size of each contact marker. 40 by default.
- **color** (*Any, optional*) – The color of contact markers. “xkcd:dark gray” by default.

Method generated by attrs for class Probe.

```
add_self_to_plot(ax: mpl_toolkits.mplot3d.axes3d.Axes3D, axis_scale_unit:  
                 brian2.units.fundamentalunits.Unit, **kwargs) → list[matplotlib.artist.Artist]
```

Add device to an existing plot

Should only be called by [plot\(\)](#).

Parameters

- **ax** (*Axes3D*) – The existing matplotlib Axes object
- **axis_scale_unit** (*Unit*) – The unit used to label axes and define chart limits

- ****kwargs** (*optional*) –

Returns A list of artists used to render the device. Needed for use in conjunction with [VideoVisualizer](#).

Return type list[Artist]

add_signals(*signals: [cleo.ephys.probes.Signal](#)) → None

Add signals to the probe for recording

Parameters *signals ([Signal](#)) – signals to add

connect_to_neuron_group(neuron_group: [brian2.groups.neurongroup.NeuronGroup](#), **kwargs: Any) → None

Configure probe to record from given neuron group

Will call [Signal.connect_to_neuron_group\(\)](#) for each signal

Parameters

- **neuron_group** ([NeuronGroup](#)) – neuron group to connect to, i.e., record from
- ****kwargs** (Any) – Passed in to signals' connect functions, needed for some signals

coords: [Quantity](#)

Coordinates of n electrodes. Must be an n x 3 array (with unit) where columns represent x, y, and z

get_state() → dict

Get current state from probe, i.e., all signals

Returns {'signal_name': value} dict with signal states

Return type dict

property n

Number of electrode contacts in the probe

probe: [Probe](#)

reset(**kwargs)

Reset the probe to a neutral state

Calls reset() on each signal

signals: list[[Signal](#)]

Signals recorded by the probe. Can be added to post-init with [add_signals\(\)](#).

to_neo() → [neo.core.group.Group](#)

Return a neo.core.AnalogSignal object with the device's data

Returns Neo object representing exported data

Return type neo.core.BaseNeo

property xs: [brian2.units.fundamentalunits.Quantity](#)

x coordinates of recording contacts

Returns x coordinates represented as a Brian quantity, that is, including units. Should be like a 1D array.

Return type Quantity

property ys: `brian2.units.fundamentalunits.Quantity`

y coordinates of recording contacts

Returns y coordinates represented as a Brian quantity, that is, including units. Should be like a 1D array.

Return type Quantity

property zs: `brian2.units.fundamentalunits.Quantity`

z coordinates of recording contacts

Returns z coordinates represented as a Brian quantity, that is, including units. Should be like a 1D array.

Return type Quantity

class `cleo.ephys.Signal(*, name: str = NOTHING)`

Bases: `abc.ABC`

Base class representing something an electrode can record

Method generated by attrs for class Signal.

brian_objects: `set`

All Brian objects created by the signal. Must be kept up-to-date for automatic injection into the network

abstract connect_to_neuron_group(*neuron_group*: `brian2.groups.neurongroup.NeuronGroup`,
***kwparams*)

Configure signal to record from specified neuron group

Parameters *neuron_group* (`NeuronGroup`) – group to record from

abstract get_state() → Any

Get the signal's current value

init_for_probe(*probe*: `cleo.ephys.probes.Probe`) → None

Called when attached to a probe.

Ensures signal can access probe and is only attached to one

Parameters *probe* (`Probe`) – Probe to attach to

Raises `ValueError` – When signal already attached to another probe

name: `str`

Unique identifier used to organize probe output. Name of the class by default.

probe: `cleo.ephys.probes.Probe`

The probe the signal is configured to record for.

reset(***kwargs*) → None

Reset signal to a neutral state

class `cleo.ephys.SortedSpiking`(*r_perfect_detection*: `brian2.units.fundamentalunits.Quantity`,
r_half_detection: `brian2.units.fundamentalunits.Quantity`,
cutoff_probability: `float = 0.01`, *, *name*: `str = NOTHING`)

Bases: `cleo.ephys.spiking.Spiking`

Detect spikes identified by neuron indices.

The indices used by the probe do not correspond to those coming from neuron groups, since the probe must consider multiple potential groups and within a group ignores those neurons that are too far away to be easily detected.

Method generated by attrs for class SortedSpiking.

connect_to_neuron_group(*neuron_group*: *brian2.groups.neurongroup.NeuronGroup*, ***kwargs*) → None

Configure sorted spiking signal to record from given neuron group

Parameters *neuron_group* (*NeuronGroup*) – group to record from

get_state() → tuple[NDArray[Any, ..., UInt[64]], NDArray[Any, ..., Float[64]], NDArray[Any, ..., UInt[64]]]

Return spikes since method was last called (i, t_ms, y)

Returns (i, t_ms, y) where i is channel (for multi-unit) or neuron (for sorted) spike indices, t_ms is spike times, and y is a spike count vector suitable for control- theoretic uses—i.e., a 0 for every channel/neuron that hasn't spiked and a 1 for a single spike.

Return type tuple[NDArray[np.uint], NDArray[float], NDArray[np.uint]]

class cleo.ephys.**Spiking**(*r_perfect_detection*: *brian2.units.fundamentalunits.Quantity*, *r_half_detection*: *brian2.units.fundamentalunits.Quantity*, *cutoff_probability*: float = 0.01, ***, *name*: str = NOTHING)

Bases: *cleo.ephys.probes.Signal*, *cleo.base.NeoExportable*

Base class for probabilistically detecting spikes

Method generated by attrs for class Spiking.

connect_to_neuron_group(*neuron_group*: *brian2.groups.neurongroup.NeuronGroup*, ***kwargs*) → numpy.ndarray

Configure signal to record from specified neuron group

Parameters *neuron_group* (*NeuronGroup*) – Neuron group to record from

Returns num_neurons_to_consider x num_channels array of spike detection probabilities, for use in subclasses

Return type np.ndarray

cutoff_probability: float

Spike detection probability below which neurons will not be considered. For computational efficiency.

abstract get_state() → tuple[NDArray[Any, ..., UInt[64]], NDArray[Any, ..., Float[64]], NDArray[Any, ..., UInt[64]]]

Return spikes since method was last called (i, t_ms, y)

Returns (i, t_ms, y) where i is channel (for multi-unit) or neuron (for sorted) spike indices, t_ms is spike times, and y is a spike count vector suitable for control- theoretic uses—i.e., a 0 for every channel/neuron that hasn't spiked and a 1 for a single spike.

Return type tuple[NDArray[np.uint], NDArray[float], NDArray[np.uint]]

i: NDArray[Any, np.uint]

Channel (for multi-unit) or neuron (for sorted) indices of spikes, stored if *save_history* on *probe*

i_probe_by_i_ng: bidict

(*neuron_group*, *i_ng*) keys, *i_probe* values. bidict for converting between neuron group indices and the indices the probe uses

r_half_detection: Quantity

Radius (with Brian unit) within which half of all spikes are detected

r_perfect_detection: Quantity
Radius (with Brian unit) within which all spikes are detected

reset(kwargs)** → None
Reset signal to a neutral state

t_ms: `NDArray[Any, float]`
Spike times in ms, stored if `save_history` on `probe`

t_samp_ms: `NDArray[Any, float]`
Sample times in ms when each spike was recorded, stored if `save_history` on `probe`

to_neo() → `neo.core.group.Group`
Return a `neo.core.AnalogSignal` object with the device's data
Returns Neo object representing exported data
Return type `neo.core.BaseNeo`

class `cleo.ephys.TKLFPSignal(uLFP_threshold_uV: float = 0.001, *, name: str = NOTHING)`
Bases: `cleo.ephys.probes.Signal`, `cleo.base.NeoExportable`
Records the Teleńczuk kernel LFP approximation.
Requires `tklfp_type='exc' | 'inh'` to specify cell type on injection.
An orientation keyword argument can also be specified on injection, which should be an array of shape `(n_neurons, 3)` representing which way is “up,” that is, towards the surface of the cortex, for each neuron. If a single vector is given, it is taken to be the orientation for all neurons in the group. `[0, 0, -1]` is the default, meaning the negative z axis is “up.” As stated elsewhere, Cleo’s convention is that `z=0` corresponds to the cortical surface and increasing z values represent increasing depth.
TKLFP is computed from spikes using the `tklfp` package.
Method generated by `attrs` for class `TKLFPSignal`.

connect_to_neuron_group(neuron_group: `brian2.groups.neurongroup.NeuronGroup`, **kwparams)
Configure signal to record from specified neuron group
Parameters `neuron_group` (`NeuronGroup`) – group to record from

get_state() → `numpy.ndarray`
Get the signal’s current value

init_for_probe(probe: `cleo.ephys.probes.Probe`)
Called when attached to a probe.
Ensures signal can access probe and is only attached to one
Parameters `probe` (`Probe`) – Probe to attach to
Raises `ValueError` – When signal already attached to another probe

lfp_uV: `numpytypes._ndarray.NDArray[Any, Any, numpytypes._number.Float]`
Approximated LFP from every call to `get_state()`. Shape is `(n_samples, n_channels)`. Stored if `save_history` on `probe`

reset(kwargs)** → None
Reset signal to a neutral state

t_ms: `numpytypes._ndarray.NDArray[Any, numpytypes._number.Float]`
Times at which LFP is recorded, in ms, stored if `save_history` on `probe`

to_neo() → neo.AnalogSignal

Return a neo.core.AnalogSignal object with the device's data

Returns Neo object representing exported data

Return type neo.core.BaseNeo

uLFP_threshold_uV: float

Threshold, in microvolts, above which the uLFP for a single spike is guaranteed to be considered, by default 1e-3. This determines the buffer length of past spikes, since the uLFP from a long-past spike becomes negligible and is ignored.

cleo.ephys.linear_shank_coords(*array_length: brian2.units.fundamentalunits.Quantity*, *channel_count: int*, *start_location: brian2.units.fundamentalunits.Quantity = array([0., 0., 0.]) * metre*, *direction: Tuple[float, float, float] = (0, 0, 1)*) → *brian2.units.fundamentalunits.Quantity*

Generate coordinates in a linear pattern

Parameters

- **array_length** (*Quantity*) – Distance from the first to the last contact (with a Brian unit)
- **channel_count** (*int*) – Number of coordinates to generate, i.e. electrode contacts
- **start_location** (*Quantity*, *optional*) – x, y, z coordinate (with unit) for the start of the electrode array, by default (0, 0, 0)*mm
- **direction** (*Tuple[float, float, float]*, *optional*) – x, y, z vector indicating the direction in which the array extends, by default (0, 0, 1), meaning pointing straight down

Returns channel_count x 3 array of coordinates, where the 3 columns represent x, y, and z

Return type Quantity

cleo.ephys.poly2_shank_coords(*array_length: brian2.units.fundamentalunits.Quantity*, *channel_count: int*, *intercol_space: brian2.units.fundamentalunits.Quantity*, *start_location: brian2.units.fundamentalunits.Quantity = array([0., 0., 0.]) * metre*, *direction: Tuple[float, float, float] = (0, 0, 1)*) → *brian2.units.fundamentalunits.Quantity*

Generate NeuroNexus-style Poly2 array coordinates

Poly2 refers to 2 parallel columns with staggered contacts. See <https://www.neuronexus.com/products/electrode-arrays/up-to-15-mm-depth> for more detail.

Parameters

- **array_length** (*Quantity*) – Length from the beginning to the end of the two-column array, as measured in the center
- **channel_count** (*int*) – Total (not per-column) number of coordinates (recording contacts) desired
- **intercol_space** (*Quantity*) – Distance between columns (with Brian unit)
- **start_location** (*Quantity*, *optional*) – Where to place the beginning of the array, by default (0, 0, 0)*mm
- **direction** (*Tuple[float, float, float]*, *optional*) – x, y, z vector indicating the direction in which the two columns extend; by default (0, 0, 1), meaning straight down.

Returns channel_count x 3 array of coordinates, where the 3 columns represent x, y, and z

Return type Quantity

```
cleo.ephys.poly3_shank_coords(array_length: brian2.units.fundamentalunits.Quantity, channel_count: int,
                             intercol_space: brian2.units.fundamentalunits.Quantity, start_location:
                             brian2.units.fundamentalunits.Quantity = array([0., 0., 0.]) * metre,
                             direction: Tuple[float, float, float] = (0, 0, 1)) →
                             brian2.units.fundamentalunits.Quantity
```

Generate NeuroNexus Poly3-style array coordinates

Poly3 refers to three parallel columns of electrodes. The middle column will be longest if the channel count isn't divisible by three and the side columns will be centered vertically with respect to the middle.

Parameters

- **array_length** (*Quantity*) – Length from beginning to end of the array as measured along the center column
- **channel_count** (*int*) – Total (not per-column) number of coordinates to generate (i.e., electrode contacts)
- **intercol_space** (*Quantity*) – Spacing between columns, with Brian unit
- **start_location** (*Quantity*, *optional*) – Location of beginning of the array, that is, the first contact in the center column, by default (0, 0, 0)*mm
- **direction** (*Tuple[float, float, float]*, *optional*) – x, y, z vector indicating the direction along which the array extends, by default (0, 0, 1), meaning straight down

Returns channel_count x 3 array of coordinates, where the 3 columns represent x, y, and z

Return type Quantity

```
cleo.ephys.tetrode_shank_coords(array_length: brian2.units.fundamentalunits.Quantity, tetrode_count: int,
                                start_location: brian2.units.fundamentalunits.Quantity = array([0., 0., 0.])
                                * metre, direction: Tuple[float, float, float] = (0, 0, 1), tetrode_width:
                                brian2.units.fundamentalunits.Quantity = 25. * umetre) →
                                brian2.units.fundamentalunits.Quantity
```

Generate coordinates for a linear array of tetrodes

See <https://www.neuronexus.com/products/electrode-arrays/up-to-15-mm-depth> to visualize NeuroNexus-style arrays.

Parameters

- **array_length** (*Quantity*) – Distance from the center of the first tetrode to the last (with a Brian unit)
- **tetrode_count** (*int*) – Number of tetrodes desired
- **start_location** (*Quantity*, *optional*) – Center location of the first tetrode in the array, by default (0, 0, 0)*mm
- **direction** (*Tuple[float, float, float]*, *optional*) – x, y, z vector determining the direction in which the linear array extends, by default (0, 0, 1), meaning straight down.
- **tetrode_width** (*Quantity*, *optional*) – Distance between contacts in a single tetrode. Not the diagonal distance, but the length of one side of the square. By default 25*umeter, as in NeuroNexus probes.

Returns (tetrode_count*4) x 3 array of coordinates, where 3 columns represent x, y, and z

Return type Quantity

`cleo.ephys.tile_coords`(*coords: brian2.units.fundamentalunits.Quantity, num_tiles: int, tile_vector: brian2.units.fundamentalunits.Quantity*) → *brian2.units.fundamentalunits.Quantity*

Tile (repeat) coordinates to produce multi-shank/matrix arrays

Parameters

- **coords** (*Quantity*) – The $n \times 3$ coordinates array to tile
- **num_tiles** (*int*) – Number of times to tile (repeat) the coordinates. For example, if you are tiling linear shank coordinates to produce multi-shank coordinates, this would be the desired number of shanks
- **tile_vector** (*Quantity*) – x, y, z array with Brian unit determining both the length and direction of the tiling

Returns $(n * \text{num_tiles}) \times 3$ array of coordinates, where the 3 columns represent x, y, and z

Return type *Quantity*

6.3.4 cleo.light module

class `cleo.light.GaussianEllipsoid`(*sigma_axial: brian2.units.fundamentalunits.Quantity = 18. * umetre, sigma_lateral: brian2.units.fundamentalunits.Quantity = 8. * umetre*)

Bases: `cleo.light.light.LightModel`

Method generated by attrs for class GaussianEllipsoid.

sigma_axial: *brian2.units.fundamentalunits.Quantity*

Standard deviation distance along the focal axis.

Standard deviations estimated by taking point where response is ~60% of peak:

Publication	axial	lateral	measure
Prakash et al., 2012	13 m	7 m	photocurrent
Rickgauer et al., 2014	8 m	4 m	Ca2+ dF/F response
Packer et al., 2015	18 m	8 m	AP probability
Chen et al., 2019	18/11 m	8/? m	AP probability/photocurrent

sigma_lateral: *brian2.units.fundamentalunits.Quantity*

Standard deviation distance along the focal plane.

transmittance(*source_coords: brian2.units.fundamentalunits.Quantity, source_dir_uvec: nptyping.types._ndarray.NDArray[Any, 3, Any], target_coords: brian2.units.fundamentalunits.Quantity*) → *nptyping.types._ndarray.NDArray[Any, Any, nptyping.types._number.Float]*

Output must be between 0 and shape (n_sources, n_targets).

viz_params(*coords: brian2.units.fundamentalunits.Quantity, direction: nptyping.types._ndarray.NDArray[Any, 3, Any], T_threshold: float, n_points_per_source: int = 4000, **kwargs*) → *brian2.units.fundamentalunits.Quantity*

Outputs $(m \times n_points_per_source \times 3)$ viz_points array, markersize_um, intensity_scale)

class `cleo.light.Light`(*coords=array([0., 0., 0.]) * metre, direction: brian2.units.fundamentalunits.Quantity = (0, 0, 1), *, name: str = NOTHING, save_history: bool = True, light_model: cleo.light.light.LightModel, wavelength: brian2.units.fundamentalunits.Quantity = 0.473 * umetre, max_Irr0_mW_per_mm2: Optional[float] = None, max_Irr0_mW_per_mm2_viz: Optional[float] = None, default_value: nptyping.types._ndarray.NDArray[Any, nptyping.types._number.Float] = NOTHING*)

Bases: `cleo.base.Stimulator`

Delivers photostimulation of the network.

Essentially “transfects” neurons and provides a light source. Under the hood, it delivers current via a Brian `Synapses` object.

Requires neurons to have 3D spatial coordinates already assigned. Also requires that the neuron model has a current term (by default `Iopto`) which is assumed to be positive (unlike the convention in many opsin modeling papers, where the current is described as negative).

See `connect_to_neuron_group()` for optional keyword parameters that can be specified when calling `cleo.CLSimulator.inject()`.

Visualization kwargs

- **n_points_per_source** (*int, optional*) – The number of points per light source used to represent light intensity in space. Default varies by `light_model`. Alias `n_points`.
- **T_threshold** (*float, optional*) – The transmittance below which no points are plotted. By default `1e-3`.
- **intensity** (*float, optional*) – How bright the light appears, should be between 0 and 1. By default `0.5`.
- **rasterized** (*bool, optional*) – Whether to render as rasterized in vector output, True by default. Useful since so many points makes later rendering and editing slow.

Method generated by attrs for class `Light`.

add_self_to_plot (*ax, axis_scale_unit, **kwargs*) → `list[matplotlib.collections.PathCollection]`

Add device to an existing plot

Should only be called by `plot()`.

Parameters

- **ax** (*Axes3D*) – The existing matplotlib Axes object
- **axis_scale_unit** (*Unit*) – The unit used to label axes and define chart limits
- ****kwargs** (*optional*) –

Returns A list of artists used to render the device. Needed for use in conjunction with `VideoVisualizer`.

Return type `list[Artist]`

connect_to_neuron_group (*neuron_group: brian2.groups.neurongroup.NeuronGroup, **kwparams: Any*) → `None`

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a `NeuronGroup`, `Synapses`, or `Monitor`, make sure to add these to `self.brian_objects`.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwparams** (*optional, passed from inject or*) – *inject*

coords: **Quantity**

(x, y, z) coords with Brian unit specifying where to place the base of the light source, by default (0, 0, 0)*mm. Can also be an nx3 array for multiple sources.

default_value: `NDArray[Any, float]`

The default value of the device—used on initialization and on `reset()`

direction: `NDArray[Any, 3, Any]`

(x, y, z) vector specifying direction in which light source is pointing, by default (0, 0, 1).

Will be converted to unit magnitude.

init_for_simulator(*sim*: `CLSimulator`) → `None`

Initialize device for simulator on initial injection

This function is called only the first time a device is injected into a simulator and performs any operations that are independent of the individual neuron groups it is connected to.

Parameters `simulator` (`CLSimulator`) – simulator being injected into

light_model: `LightModel`

LightModel object defining how light is emitted. See [OpticFiber](#) for an example.

max_Irr0_mW_per_mm2: `float`

The maximum irradiance the light source can emit.

Usually determined by hardware in a real experiment.

max_Irr0_mW_per_mm2_viz: `float`

Maximum irradiance for visualization purposes.

i.e., the level at or above which the light appears maximally bright. Only relevant in video visualization.

property `n`

Number of light sources

property `source:` `brian2.groups.subgroup.Subgroup`

to_neo()

Return a `neo.core.AnalogSignal` object with the device's data

Returns Neo object representing exported data

Return type `neo.core.BaseNeo`

transmittance(*target_coords*) → `numpy.ndarray`

update(*value*: `Union[float, numpy.ndarray]`) → `None`

Set the light intensity, in mW/mm2 (without unit) for 1P excitation or laser power (mW) for 2P excitation (GaussianEllipsoid light model).

Parameters `Irr0_mW_per_mm2` (`float`) – Desired light intensity for light source

update_artists(*artists*: `list[matplotlib.artist.Artist]`, *value*, **args*, ***kwargs*) → `list[matplotlib.artist.Artist]`

Update the artists used to render the device

Used to set the artists' state at every frame of a video visualization. The current state would be passed in **args* or ***kwargs*

Parameters `artists` (`list[Artist]`) – the artists used to render the device originally, i.e., which were returned from the first [add_self_to_plot\(\)](#) call.

Returns The artists that were actually updated. Needed for efficient blit rendering, where only updated artists are re-rendered.

Return type `list[Artist]`

wavelength: Quantity

light wavelength with unit (usually nmeter)

```
class cleo.light.LightDependent(spectrum: list[tuple[float, float]] = NOTHING, spectrum_interpolator:
                                typing.Callable = <function cubic_interpolator>)
```

Bases: object

Mix-in class for opsin and indicator. TODO

We approximate dynamics under multiple wavelengths using a weighted sum of photon fluxes, where the factor indicates the activation relative to the peak-sensitivity wavelength for a equivalent power, which most papers report. When they report the action spectrum for equivalent photon flux instead (see Mager et al, 2018), use `equal_photon_flux_spectrum()`. This weighted sum is an approximation of a nonlinear peak-non-peak wavelength relation; see `notebooks/multi_wavelength_model.ipynb` for details.

Method generated by attrs for class LightDependent.

epsilon(*lambda_new*) → float

Returns the epsilon value for a given lambda (in nm) representing the relative sensitivity of the opsin to that wavelength.

property light_agg_ngs

spectrum: list[tuple[float, float]]

List of (wavelength, epsilon) tuples representing the action (opsin) or excitation (indicator) spectrum.

spectrum_interpolator: Callable

Function of signature (*lambdas_nm, epsilons, lambda_new_nm*) that interpolates the action spectrum data and returns $\varepsilon \in [0, 1]$ for the new wavelength.

```
class cleo.light.LightModel
```

Bases: abc.ABC

Method generated by attrs for class LightModel.

```
abstract transmittance(source_coords: brian2.units.fundamentalunits.Quantity, source_direction:
                        nptyping.types._ndarray.NDArray[Any, 3, Any], target_coords:
                        brian2.units.fundamentalunits.Quantity) →
                        nptyping.types._ndarray.NDArray[Any, Any, nptyping.types._number.Float]
```

Output must be between 0 and shape (n_sources, n_targets).

```
abstract viz_params(coords: brian2.units.fundamentalunits.Quantity, direction:
                      nptyping.types._ndarray.NDArray[Any, 3, Any], T_threshold: float,
                      n_points_per_source: Optional[int] = None, **kwargs) →
                      brian2.units.fundamentalunits.Quantity
```

Outputs (m x n_points_per_source x 3 viz_points array, markersize_um, intensity_scale)

```
class cleo.light.OpticFiber(R0: brian2.units.fundamentalunits.Quantity = 100. * umetre, NAfib:
                             brian2.units.fundamentalunits.Quantity = 0.37, K:
                             brian2.units.fundamentalunits.Quantity = 125. * metre ** - 1, S:
                             brian2.units.fundamentalunits.Quantity = 7370. * metre ** - 1, ntis:
                             brian2.units.fundamentalunits.Quantity = 1.36)
```

Bases: `cleo.light.light.LightModel`

Optic fiber light model from Foutz et al., 2012.

Defaults are from paper for 473 nm wavelength.

Method generated by attrs for class OpticFiber.

K: `brian2.units.fundamentalunits.Quantity`

absorbance coefficient (wavelength/tissue dependent)

NAfib: `brian2.units.fundamentalunits.Quantity`

optical fiber numerical aperture

R0: `brian2.units.fundamentalunits.Quantity`

optical fiber radius

S: `brian2.units.fundamentalunits.Quantity`

scattering coefficient (wavelength/tissue dependent)

ntis: `brian2.units.fundamentalunits.Quantity`

tissue index of refraction (wavelength/tissue dependent)

transmittance(*source_coords: brian2.units.fundamentalunits.Quantity, source_dir_uvec: nptyping.types._ndarray.NDArray[Any, 3, Any], target_coords: brian2.units.fundamentalunits.Quantity*) → `nptyping.types._ndarray.NDArray[Any, Any, nptyping.types._number.Float]`

Output must be between 0 and shape (n_sources, n_targets).

viz_params(*coords: brian2.units.fundamentalunits.Quantity, direction: nptyping.types._ndarray.NDArray[Any, 3, Any], T_threshold: float, n_points_per_source: int = 4000, **kwargs*) → `brian2.units.fundamentalunits.Quantity`

Outputs (m x n_points_per_source x 3 viz_points array, markersize_um, intensity_scale)

`cleo.light.cubic_interpolator(lambdas_nm, epsilons, lambda_new_nm)`

`cleo.light.equal_photon_flux_spectrum(spectrum: list[tuple[float, float]]) → list[tuple[float, float]]`

Converts an equival photon flux spectrum to an equal power density spectrum.

`cleo.light.fiber473nm(R0=100. * umetre, NAfib=0.37, K=125. * metre ** - 1, S=7370. * metre ** - 1, ntis=1.36) → cleo.light.light.OpticFiber`

Light parameters for 473 nm wavelength delivered via an optic fiber.

From Foutz et al., 2012. See [OpticFiber](#) for parameter descriptions.

`cleo.light.linear_interpolator(lambdas_nm, epsilons, lambda_new_nm)`

`cleo.light.plot_spectra(*lids: cleo.light.light_dependence.LightDependent)`

`cleo.light.tp_light_from_scope(scope, wavelength=1.06 * umetre, **kwargs)`

6.3.5 cleo.opto module

Contains opsin models, light sources, and some parameters

```

class cleo.opto.BansalFourStateOpsin(extra_namespace: dict = NOTHING, spectrum: list[tuple[float,
float]] = NOTHING, spectrum_interpolator: typing.Callable =
<function cubic_interpolator>, Gd1:
brian2.units.fundamentalunits.Quantity = 66. * hertz, Gd2:
brian2.units.fundamentalunits.Quantity = 10. * hertz, Gr0:
brian2.units.fundamentalunits.Quantity = 0.333 * hertz, g0:
brian2.units.fundamentalunits.Quantity = 3.2 * nsiemens, phim:
brian2.units.fundamentalunits.Quantity = 1.e+22 * metre ** -2 *
second ** -1, k1: brian2.units.fundamentalunits.Quantity = 0.4 *
khertz, k2: brian2.units.fundamentalunits.Quantity = 120. * hertz,
Gf0: brian2.units.fundamentalunits.Quantity = 18. * hertz, Gb0:
brian2.units.fundamentalunits.Quantity = 8. * hertz, kf:
brian2.units.fundamentalunits.Quantity = 10. * hertz, kb:
brian2.units.fundamentalunits.Quantity = 8. * hertz, gamma:
brian2.units.fundamentalunits.Quantity = 0.05, p:
brian2.units.fundamentalunits.Quantity = 1, q:
brian2.units.fundamentalunits.Quantity = 1, E:
brian2.units.fundamentalunits.Quantity = 0. * volt, *, name: str =
NOTHING, save_history: bool = True)

```

Bases: `cleo.opto.opsins.MarkovOpsin`

4-state model from Bansal et al. 2020.

The difference from the PyRhO model is that there is no voltage dependence.

`rho_rel` is channel density relative to standard model fit; modifying it post-injection allows for heterogeneous opsin expression.

`IOPTO_VAR_NAME` and `V_VAR_NAME` are substituted on injection.

Method generated by attrs for class `BansalFourStateOpsin`.

E: Quantity

Gb0: Quantity

Gd1: Quantity

Gd2: Quantity

Gf0: Quantity

Gr0: Quantity

g0: Quantity

gamma: Quantity

init_syn_vars(*opto_syn*: `brian2.synapses.synapses.Synapses`) → None

Initializes appropriate variables in `Synapses` implementing the model

Can also be used to reset the variables.

Parameters `syn` (`Synapses`) – The `synapses` object implementing this model

k1: Quantity

k2: Quantity

kb: Quantity

kf: Quantity

model: str

Basic Brian model equations string.

Should contain a *rho_rel* term reflecting relative expression levels. Will likely also contain special NeuronGroup-dependent symbols such as V_VAR_NAME to be replaced on injection in `modify_model_and_params_for_ng()`.

p: Quantity

phim: Quantity

q: Quantity

```
class cleo.opto.FourStateOpsin(spectrum: list[tuple[float, float]] = NOTHING, spectrum_interpolator:
    typing.Callable = <function cubic_interpolator>, g0:
    brian2.units.fundamentalunits.Quantity = 114. * nsiemens, gamma:
    brian2.units.fundamentalunits.Quantity = 0.00742, phim:
    brian2.units.fundamentalunits.Quantity = 2.33e+23 * metre ** -2 * second
    ** -1, k1: brian2.units.fundamentalunits.Quantity = 4.15 * khertz, k2:
    brian2.units.fundamentalunits.Quantity = 0.868 * khertz, p:
    brian2.units.fundamentalunits.Quantity = 0.833, Gf0:
    brian2.units.fundamentalunits.Quantity = 37.3 * hertz, kf:
    brian2.units.fundamentalunits.Quantity = 58.1 * hertz, Gb0:
    brian2.units.fundamentalunits.Quantity = 16.1 * hertz, kb:
    brian2.units.fundamentalunits.Quantity = 63. * hertz, q:
    brian2.units.fundamentalunits.Quantity = 1.94, Gd1:
    brian2.units.fundamentalunits.Quantity = 105. * hertz, Gd2:
    brian2.units.fundamentalunits.Quantity = 13.8 * hertz, Gr0:
    brian2.units.fundamentalunits.Quantity = 0.33 * hertz, E:
    brian2.units.fundamentalunits.Quantity = 0. * volt, v0:
    brian2.units.fundamentalunits.Quantity = 43. * mvolt, v1:
    brian2.units.fundamentalunits.Quantity = 17.1 * mvolt, *, name: str =
    NOTHING, save_history: bool = True)
```

Bases: `cleo.opto.opsins.MarkovOpsin`

4-state model from PyRhO (Evans et al. 2016).

rho_rel is channel density relative to standard model fit; modifying it post-injection allows for heterogeneous opsin expression.

IOPTO_VAR_NAME and V_VAR_NAME are substituted on injection.

Defaults are for Chr2.

Method generated by attrs for class `FourStateOpsin`.

E: Quantity

Gb0: Quantity

Gd1: Quantity

Gd2: Quantity

Gf0: Quantity

Gr0: Quantity

extra_namespace: `dict[str, Any]`

Additional items (beyond parameters) to be added to the opto synapse namespace

g0: `Quantity`

gamma: `Quantity`

init_syn_vars(*opto_syn*: `brian2.synapses.synapses.Synapses`) \rightarrow `None`

Initializes appropriate variables in Synapses implementing the model

Can also be used to reset the variables.

Parameters `syn` (`Synapses`) – The synapses object implementing this model

k1: `Quantity`

k2: `Quantity`

kb: `Quantity`

kf: `Quantity`

model: `str`

Basic Brian model equations string.

Should contain a *rho_rel* term reflecting relative expression levels. Will likely also contain special NeuronGroup-dependent symbols such as `V_VAR_NAME` to be replaced on injection in `modify_model_and_params_for_ng()`.

p: `Quantity`

phim: `Quantity`

q: `Quantity`

v0: `Quantity`

v1: `Quantity`

class `cleo.opto.Opsin`(*extra_namespace*: `dict = NOTHING`, *spectrum*: `list[tuple[float, float]] = NOTHING`, *spectrum_interpolator*: `typing.Callable = <function cubic_interpolator>`, ***, *name*: `str = NOTHING`, *save_history*: `bool = True`)

Bases: `cleo.light.light_dependence.LightDependent`, `cleo.base.SynapseDevice`

Base class for opsin model.

We approximate dynamics under multiple wavelengths using a weighted sum of photon fluxes, where the factor indicates the activation relative to the peak-sensitivity wavelength for an equivalent number of photons (see Mager et al, 2018). This weighted sum is an approximation of a nonlinear peak-non-peak wavelength relation; see `notebooks/multi_wavelength_model.ipynb` for details.

Method generated by attrs for class Opsin.

property `action_spectrum`

Alias for `light_receptor.spectrum`

spectrum: `list[tuple[float, float]]`

List of (wavelength, epsilon) tuples representing the action (opsin) or excitation (indicator) spectrum.

spectrum_interpolator: Callable

Function of signature (lambdas_nm, epsilons, lambda_new_nm) that interpolates the action spectrum data and returns $\varepsilon \in [0, 1]$ for the new wavelength.

```
class cleo.opto.ProportionalCurrentOpsin(extra_namespace: dict = NOTHING, spectrum:  
                                         list[tuple[float, float]] = NOTHING, spectrum_interpolator:  
                                         typing.Callable = <function cubic_interpolator>, *, name: str  
                                         = NOTHING, save_history: bool = True, I_per_Irr:  
                                         brian2.units.fundamentalunits.Quantity)
```

Bases: `cleo.opto.opsins.Opsin`

A simple model delivering current proportional to light intensity

Method generated by attrs for class ProportionalCurrentOpsin.

I_per_Irr: Quantity

How much current (in amps or unitless, depending on neuron model) to deliver per mW/mm2.

model: str

Basic Brian model equations string.

Should contain a *rho_rel* term reflecting relative expression levels. Will likely also contain special NeuronGroup-dependent symbols such as V_VAR_NAME to be replaced on injection in `modify_model_and_params_for_ng()`.

required_vars: list[Tuple[str, Unit]]

Default names of state variables required in the neuron group, along with units, e.g., [('Iopto', amp)].

It is assumed that non-default values can be passed in on injection as a keyword argument [`default_name`]_{var_name}=[`non_default_name`] and that these are found in the model string as [`DEFAULT_NAME`]_{VAR_NAME} before replacement.

`cleo.opto.chr2_4s()` → `cleo.opto.opsins.FourStateOpsin`

Returns a 4-state Chr2 model.

Params taken from [try.projectpyrho.org](https://tryprojectpyrho.org)'s default 4-state configuration. Action spectrum from Nagel et al., 2003, Fig. 4a, extracted using [Plot Digitizer](#).

Parameters can be changed after initialization but *before injection*.

`cleo.opto.chr2_b4s()` → `cleo.opto.opsins.BansalFourStateOpsin`

Returns a 4-state Vf-Chrimson model.

Params given in Bansal et al., 2020. Action spectrum from Nagel et al., 2003, Fig. 4a, extracted using [Plot Digitizer](#).

Parameters can be changed after initialization but *before injection*.

`cleo.opto.chrimson_4s()` → `cleo.opto.opsins.BansalFourStateOpsin`

Returns a 4-state Chrimson model.

Params given in Bansal et al., 2020. Action spectrum from Mager et al., 2018, Supp. Fig. 1a, extracted using [Plot Digitizer](#).

Parameters can be changed after initialization but *before injection*.

`cleo.opto.gtacr2_4s()` → `cleo.opto.opsins.BansalFourStateOpsin`

Returns a 4-state model of GtACR2, an anion channel.

Params given in Bansal et al., 2020. Action spectra from Govorunova et al., 2015, Fig. 1f, extracted using [Plot Digitizer](#).

Parameters can be changed after initialization but *before injection*.

`cleo.opto.vfchromson_4s()` → `cleo.opto.opsins.BansalFourStateOpsin`

Returns a 4-state Vf-Chrimson model.

Params given in Bansal et al., 2020. Action spectrum from Mager et al., 2018, Supp. Fig. 1a, extracted using [Plot Digitizer](#).

Parameters can be changed after initialization but *before injection*.

6.3.6 cleo.imaging module

class `cleo.imaging.CalBindingActivationModel(model: str)`

Bases: `object`

Base class for modeling calcium binding/activation

Method generated by attrs for class `CalBindingActivationModel`.

model: `str`

class `cleo.imaging.CalciumModel(model: str)`

Bases: `object`

Base class for how GECI computes calcium concentration

Must provide variable Ca (molar) in model.

Method generated by attrs for class `CalciumModel`.

model: `str`

on_pre: `str`

class `cleo.imaging.DoubExpCalBindingActivation(*, A: float, tau_on: brian2.units.fundamentalunits.Quantity, tau_off: brian2.units.fundamentalunits.Quantity, Ca_rest: brian2.units.fundamentalunits.Quantity)`

Bases: `cleo.imaging.sensors.CalBindingActivationModel`

Double exponential kernel convolution representing CaB binding/activation.

Convolution is implemented via ODEs; see `notebooks/double_exp_conv_as_ode.ipynb` for derivation.

A, tau_on, and tau_off are the versions with proper scale and units of NAOMi's ca_amp, t_on, and t_off.

Some parameters found [here](#). Parameters fit [here](#).

Method generated by attrs for class `DoubExpCalBindingActivation`.

A: `float`

amplitude of double exponential kernel

Ca_rest: `brian2.units.fundamentalunits.Quantity`

Resting Ca²⁺ concentration (molar).

init_syn_vars(syn: `brian2.synapses.synapses.Synapses`) → `None`

model: `str`

tau_off: `brian2.units.fundamentalunits.Quantity`
CaB unbinding/deactivation time constant (sec)

tau_on: `brian2.units.fundamentalunits.Quantity`
CaB binding/activation time constant (sec)

class `cleo.imaging.DynamicCalcium(*, Ca_rest: brian2.units.fundamentalunits.Quantity, gamma: brian2.units.fundamentalunits.Quantity, B_T: brian2.units.fundamentalunits.Quantity, kappa_S: float, dCa_T: brian2.units.fundamentalunits.Quantity)`

Bases: `cleo.imaging.sensors.CalciumModel`

Pieced together from Lütke et al., 2013; Helmchen and Tank, 2015; and Song et al., 2021. ([code](#))

Method generated by attrs for class DynamicCalcium.

B_T: `brian2.units.fundamentalunits.Quantity`
total indicator (buffer) concentration (molar)

Ca_rest: `brian2.units.fundamentalunits.Quantity`
resting Ca2+ concentration (molar)

dCa_T: `brian2.units.fundamentalunits.Quantity`
total Ca2+ concentration increase per spike (molar)

gamma: `brian2.units.fundamentalunits.Quantity`
clearance/extrusion rate (1/sec)

init_syn_vars(*syn*: `brian2.synapses.synapses.Synapses`) → None

kappa_S: `float`
Ca2+ binding ratio of the endogenous buffer

model: `str`
from eq 8 in Lütke et al., 2013

on_pre: `str`
from eq 9 in Lütke et al., 2013

class `cleo.imaging.ExcitationModel(model: str)`

Bases: `object`

Defines `exc_factor`

Method generated by attrs for class ExcitationModel.

model: `str`

class `cleo.imaging.GECI(extra_namespace: dict = NOTHING, *, name: str = NOTHING, save_history: bool = True, sigma_noise: float, dFF_IAP: float, cal_model: cleo.imaging.sensors.CalciumModel, bind_act_model: cleo.imaging.sensors.CalBindingActivationModel, exc_model: cleo.imaging.sensors.ExcitationModel, K_d: brian2.units.fundamentalunits.Quantity, n_H: float, dFF_max: float)`

Bases: `cleo.imaging.sensors.Sensor`

GECI model based on Song et al., 2021, with interchangeable components.

See `geci()` for a convenience function for creating GECI models.

As a potentially simpler alternative for future work, see the phenomenological S2F model from [Zhang et al., 2023](#). While parameter count looks similar, at least they have parameters fit already, and directly to data, rather than to biophysical processes before the data.

Method generated by attrs for class GECL.

K_d: Quantity

indicator dissociation constant (binding affinity) (molar)

bind_act_model: [*CalBindingActivationModel*](#)

cal_model: [*CalciumModel*](#)

dFF_max: float

amplitude of Hill equation for conversion from Ca2+ to F/F, Fmax/F0. May be approximated from ‘dynamic range’ in literature Fmax/Fmin

exc_model: [*ExcitationModel*](#)

fluor_model: str

Uses a Hill equation to convert from Ca2+ to F/F, as in Song et al., 2021

get_state() → dict[brian2.groups.neurongroup.NeuronGroup, numpy.ndarray]

Returns a list of arrays in the order neuron groups/targets were received.

Signals should be normalized to baseline of 0 and 1 corresponding to an action potential peak.

init_syn_vars(syn: [*brian2.synapses.synapses.Synapses*](#)) → None

Initializes appropriate variables in Synapses implementing the model

Can also be used to reset the variables.

Parameters **syn** (*Synapses*) – The synapses object implementing this model

location: str

cytoplasm or membrane

n_H: float

Hill coefficient for conversion from Ca2+ to F/F

property params: dict

Returns a dictionary of parameters for the model

class cleo.imaging.**LightExcitation**

Bases: [*cleo.imaging.sensors.ExcitationModel*](#)

Method generated by attrs for class LightExcitation.

model: str

class cleo.imaging.**NullBindingActivation**

Bases: [*cleo.imaging.sensors.CalBindingActivationModel*](#)

Doesn’t model binding/activation; i.e., goes straight from [Ca2+] to F/F

Method generated by attrs for class NullBindingActivation.

model: str

class cleo.imaging.NullExcitation

Bases: [cleo.imaging.sensors.ExcitationModel](#)

Method generated by attrs for class NullExcitation.

model: str

class cleo.imaging.PreexistingCalcium

Bases: [cleo.imaging.sensors.CalciumModel](#)

Calcium concentration is pre-existing in neuron model

Method generated by attrs for class PreexistingCalcium.

model: str

class cleo.imaging.Scope(*sensor: cleo.imaging.sensors.Sensor, img_width: brian2.units.fundamentalunits.Quantity, focus_depth: Optional[brian2.units.fundamentalunits.Quantity] = None, location: brian2.units.fundamentalunits.Quantity = array([0., 0., 0.]) * metre, direction: brian2.units.fundamentalunits.Quantity = (0, 0, 1), soma_radius: brian2.units.fundamentalunits.Quantity = 10. * umetre, snr_cutoff: float = 1, rand_seed: Optional[int] = None, *, name: str = NOTHING, save_history: bool = True*)

Bases: [cleo.base.Recorder](#)

rho_rel_generator : callable, takes n, outputs float array

Method generated by attrs for class Scope.

add_self_to_plot(*ax: mpl_toolkits.mplot3d.axes3d.Axes3D, axis_scale_unit: brian2.units.fundamentalunits.Unit, **kwargs*) → list[matplotlib.artist.Artist]

Add device to an existing plot

Should only be called by [plot\(\)](#).

Parameters

- **ax** (*Axes3D*) – The existing matplotlib Axes object
- **axis_scale_unit** (*Unit*) – The unit used to label axes and define chart limits
- ****kwargs** (*optional*) –

Returns A list of artists used to render the device. Needed for use in conjunction with [VideoVisualizer](#).

Return type list[Artist]

connect_to_neuron_group(*neuron_group: brian2.groups.neurongroup.NeuronGroup, **kwparams*) → None

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a NeuronGroup, Synapses, or Monitor, make sure to add these to *self.brian_objects*.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwparams** (*optional, passed from inject or*) – *inject*

dFF: `list[NDArray[Any, float]]`
 F/F from every call to `get_state()`. Shape is (n_samples, n_ROIs). Stored if `save_history`

property dFF_1AP
 gets dFF_1AP for all targets, in order injected. Varies with expression levels.

direction: `NDArray[3, float]`

focus_coords_per_injct: `list[NDArray[Any, float]]`

focus_depth: `Quantity`

get_state() \rightarrow `nptyping.types._ndarray.NDArray[Any, nptyping.types._number.Float]`
 Return current measurement.

i_targets_for_neuron_group(*neuron_group*)
 can handle multiple injections into same ng

i_targets_per_injct: `list[NDArray[Any, int]]`

img_width: `Quantity`

inject_sensor_for_targets(***kwparams*) \rightarrow `None`

location: `Quantity`

property n

neuron_groups: `list[NeuronGroup]`

rand_seed: `int`

reset(***kwargs*) \rightarrow `None`
 Reset the device to a neutral state

rho_rel_per_injct: `list[NDArray[Any, float]]`

sensor: `Sensor`

property sigma_noise
 gets noise for all targets, in order injected.

sigma_per_injct: `list[NDArray[Any, float]]`

snr_cutoff: `float`
 applied only when not focus_depth is not None

soma_radius: `Quantity`

t_ms: `list[float]`
 Times at which sensor traces are recorded, in ms, stored if `save_history`

target_neurons_in_plane(*ng*, *focus_depth: Optional[brian2.units.fundamentalunits.Quantity] = None*,
soma_radius: Optional[brian2.units.fundamentalunits.Quantity] = None) \rightarrow
`tuple[NDArray[Any, Int[64]], NDArray[Any, Float[64]]]`

```
class cleo.imaging.Sensor(extra_namespace: dict = NOTHING, *, name: str = NOTHING, save_history:  
                        bool = True, sigma_noise: float, dFF_1AP: float, location: str)
```

Bases: [cleo.base.SynapseDevice](#)

Base class for sensors

Method generated by attrs for class Sensor.

dFF_1AP: float

F/F for 1 AP, only used for scope SNR cutoff

property exc_spectrum: list[tuple[float, float]]

Excitation spectrum, alias for spectrum

get_state() → dict[brian2.groups.neurongroup.NeuronGroup, numpy.ndarray]

Returns a list of arrays in the order neuron groups/targets were received.

Signals should be normalized to baseline of 0 and 1 corresponding to an action potential peak.

location: str

cytoplasm or membrane

sigma_noise: float

standard deviation of Gaussian noise in F/F measurement

property snr: float

Signal-to-noise ratio for 1 AP

```
cleo.imaging.gcamp3(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=287. * nmolar,  
                  n_H=2.52, dFF_max=12, sigma_noise=0.018571428571428572,  
                  dFF_1AP=0.039, ca_amp=0.05, t_on=1, t_off=1, Ca_rest=50. * nmolar, kappa_S=110,  
                  gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar, name='gcamp3',  
                  **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is None, then we don't have it fit yet.

```
cleo.imaging.gcamp6f(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=290. * nmolar,  
                   n_H=2.7, dFF_max=25.2, sigma_noise=0.03748181818181818,  
                   dFF_1AP=0.09775500000000001, ca_amp=76.1251, t_on=0.8535, t_off=98.6173,  
                   Ca_rest=50. * nmolar, kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar,  
                   dCa_T=7.6 * umolar, name='gcamp6f', **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is None, then we don't have it fit yet.

```
cleo.imaging.gcamp6rs06(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=0.32 * umolar,  
                      n_H=3, dFF_max=15, sigma_noise=0.030227272727272728,  
                      dFF_1AP=None, ca_amp=None, t_on=None, t_off=None, Ca_rest=50. * nmolar,  
                      kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar,  
                      name='gcamp6rs06', **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

Don't know sigma_noise

```
cleo.imaging.gcamp6rs09(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=0.52 *
    umolar, n_H=3.2, dFF_max=25, sigma_noise=0.030227272727272728,
    dFF_IAP=None, ca_amp=None, t_on=None, t_off=None, Ca_rest=50. * nmolar,
    kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar,
    name='gcamp6rs09', **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

Don't know sigma_noise

```
cleo.imaging.gcamp6s(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=147. *
    nmolar, n_H=2.45, dFF_max=27.2, sigma_noise=0.030227272727272728,
    dFF_IAP=0.133, ca_amp=54.6943, t_on=0.4526, t_off=68.5461, Ca_rest=50. *
    nmolar, kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 *
    umolar, name='gcamp6s', **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

```
cleo.imaging.geci(light_dependent, doub_exp_conv, pre_existing_cal, **kwargs)
```

```
cleo.imaging.jgcamp7b(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=82. *
    nmolar, n_H=3.06, dFF_max=22.1, sigma_noise=0.007556818181818182,
    dFF_IAP=None, ca_amp=None, t_on=None, t_off=None, Ca_rest=50. * nmolar,
    kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar,
    name='jgcamp7b', **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

```
cleo.imaging.jgcamp7c(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=298. *
    nmolar, n_H=2.44, dFF_max=145.6, sigma_noise=0.011788636363636364,
    dFF_IAP=None, ca_amp=None, t_on=None, t_off=None, Ca_rest=50. * nmolar,
    kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar,
    name='jgcamp7c', **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

```
cleo.imaging.jgcamp7f(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=174. *
    nmolar, n_H=2.3, dFF_max=30.2, sigma_noise=0.021763636363636363,
    dFF_IAP=None, ca_amp=None, t_on=None, t_off=None, Ca_rest=50. * nmolar,
    kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar,
    name='jgcamp7f', **kwargs)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

```
cleo.imaging.jgcamp7s(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=68. *
                      nmolar, n_H=2.49, dFF_max=40.4, sigma_noise=0.009975000000000001,
                      dFF_IAP=None, ca_amp=None, t_on=None, t_off=None, Ca_rest=50. * nmolar,
                      kappa_S=110, gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar,
                      name='jgcamp7s', **kwparams)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

```
cleo.imaging.ogb1(light_dependent=False, doub_exp_conv=True, pre_existing_cal=False, K_d=250. * nmolar,
                  n_H=1, dFF_max=14, sigma_noise=0.030227272727272728, dFF_IAP=None,
                  ca_amp=None, t_on=None, t_off=None, Ca_rest=50. * nmolar, kappa_S=110,
                  gamma=292.3 * hertz, B_T=200. * umolar, dCa_T=7.6 * umolar, name='ogb1',
                  **kwparams)
```

Returns a GECI model with parameters given in [NAOMi's code](#) (Song et al., 2021) as well as Dana et al., 2019 and Zhang et al., 2023.

Only those parameters used in chosen model components apply. If the default is `None`, then we don't have it fit yet.

Don't know sigma_noise

```
cleo.imaging.target_neurons_in_plane(ng, scope_focus_depth, scope_img_width,
                                     scope_location=array([0., 0., 0.]) * metre, scope_direction=(0, 0,
                                     1), soma_radius=10. * umetre, sensor_location='cytoplasm')
```

Returns a tuple of (i_targets, noise_focus_factor, coords_on_plane)

6.3.7 cleo.ioproc module

Classes and functions for constructing and configuring an *IOProcessor*.

class cleo.ioproc.ConstantDelay(delay_ms: float)

Bases: [cleo.ioproc.delays.Delay](#)

Simply adds a constant delay to the computation

Parameters delay_ms (float) – Desired delay in milliseconds

compute()

Compute delay.

class cleo.ioproc.Delay

Bases: abc.ABC

Abstract base class for computing delays.

abstract compute() → float

Compute delay.

class cleo.ioproc.FiringRateEstimator(*tau_ms: float, sample_period_ms: float, **kwargs*)

Bases: [cleo.ioproc.base.ProcessingBlock](#)

Exponential filter to estimate firing rate.

Requires *sample_time_ms* kwarg when process is called.

Parameters

- **tau_ms** (*float*) – Time constant of filter
- **sample_period_ms** (*float*) – Sampling period in milliseconds

compute_output (*input: nptyping.types._ndarray.NDArray[Any, nptyping.types._number.UInt], **kwargs*)
→ *nptyping.types._ndarray.NDArray[Any, nptyping.types._number.Float]*

Estimate firing rate given past and current spikes.

Parameters *input* (*NDArray[(n,), np.uint]*) – n-length vector of spike counts

Keyword Arguments *sample_time_ms* (*float*) – Time measurement was taken in milliseconds

Returns n-length vector of firing rates

Return type *NDArray[(n,), float]*

delay: [Delay](#)

The delay object determining compute latency for the block

save_history: **bool**

Whether to record *t_in_ms*, *t_out_ms*, and *values* with every timestep

t_in_ms: *list[float]*

The walltime the block received each input. Only recorded if *save_history*

t_out_ms: *list[float]*

The walltime of each of the block's outputs. Only recorded if *save_history*

values: *list[Any]*

Each of the block's outputs. Only recorded if *save_history*

class cleo.ioproc.GaussianDelay(*loc: float, scale: float*)

Bases: [cleo.ioproc.delays.Delay](#)

Generates normal-distributed delay.

Will return 0 when a negative value is sampled.

Parameters

- **loc** (*float*) – Center of distribution
- **scale** (*float*) – Standard deviation of delay distribution

compute() → *float*

Compute delay.

class cleo.ioproc.LatencyIOProcessor(*sample_period_ms: float, **kwargs*)

Bases: [cleo.base.IOProcessor](#)

IOProcessor capable of delivering stimulation some time after measurement.

Parameters *sample_period_ms* (*float*) – Determines how frequently samples are taken from the network.

Keyword Arguments

- **sampling** (*str*) – “fixed” or “when idle”; “fixed” by default
 “fixed” sampling means samples are taken on a fixed schedule, with no exceptions.
 “when idle” sampling means no samples are taken before the previous sample’s output has been delivered. A sample is taken ASAP after an over-period computation: otherwise remains on schedule.
- **processing** (*str*) – “parallel” or “serial”; “parallel” by default
 “parallel” computes the output time by adding the delay for a sample onto the sample time, so if the delay is 2 ms, for example, while the sample period is only 1 ms, some of the processing is happening in parallel. Output order matches input order even if the computed output time for a sample is sooner than that for a previous sample.
 “serial” computes the output time by adding the delay for a sample onto the output time of the previous sample, rather than the sampling time. Note this may be of limited utility because it essentially means the *entire* round trip cannot be in parallel at all. More realistic is that simply each block or phase of computation must be serial. If anyone cares enough about this, it will have to be implemented in the future.

Note: It doesn’t make much sense to combine parallel computation with “when idle” sampling, because “when idle” sampling only produces one sample at a time to process.

Raises **ValueError** – For invalid *sampling* or *processing* kwargs

get_ctrl_signal(*query_time_ms*)

Get per-stimulator control signal from the *IOProcessor*.

Parameters **time** (*Brian 2 temporal Unit*) – Current timestep

Returns A {‘stimulator_name’: value} dictionary for updating stimulators.

Return type dict

is_sampling_now(*query_time_ms*)

Determines whether the processor will take a sample at this timestep.

Parameters **time** (*Brian 2 temporal Unit*) – Current timestep.

Return type bool

abstract process(*state_dict: dict, sample_time_ms: float*) → Tuple[dict, float]

Process network state to generate output to update stimulators.

This is the function the user must implement to define the signal processing pipeline.

Parameters

- **state_dict** (*dict*) – {*recorder_name: state*} dictionary from *get_state()*
- **time_ms** (*float*) –

Returns {‘stim_name’: *ctrl_signal*} dictionary and output time in milliseconds.

Return type Tuple[dict, float]

put_state(*state_dict*: dict, *sample_time_ms*)

Deliver network state to the IOProcessor.

Parameters

- **state_dict** (dict) – A dictionary of recorder measurements, as returned by [get_state\(\)](#)
- **time** (*brian2 temporal Unit*) – The current simulation timestep. Essential for simulating control latency and for time-varying control.

t_samp_ms: list[float]

Record of sampling times—each time [put_state\(\)](#) is called.

class cleo.ioproc.PIDcontroller(*ref_signal*: callable, *Kp*: float, *Ki*: float = 0, *sample_period_ms*: float = 0, ***kwargs*: Any)

Bases: [cleo.ioproc.base.ProcessingBlock](#)

Simple PI controller.

[compute_output\(\)](#) requires a *sample_time_ms* keyword argument. Only tested on controlling scalar values, but could be easily adapted to controlling a multi-dimensional state.

Parameters

- **ref_signal** (callable) – Must return the target as a function of time in ms
- **Kp** (float) – Gain on the proportional error
- **Ki** (float, optional) – Gain on the integral error, by default 0
- **sample_period_ms** (float, optional) – Rate at which processor takes samples, by default 0. Only used to compute integrated error on first sample

compute_output(*input*: float, ***kwargs*) → float

Compute control input to the system using previously specified gains.

Parameters **input** (Any) – Current system state

Returns Control signal

Return type float

ref_signal: callable[[float], Any]

Callable returning the target as a function of time in ms

class cleo.ioproc.ProcessingBlock(***kwargs*)

Bases: abc.ABC

Abstract signal processing stage or control block.

It's important to use *super().__init__*(***kwargs*) in the base class to use the parent-class logic here.

Keyword Arguments **delay** ([Delay](#)) – Delay object which adds to the compute time

Raises **TypeError** – When *delay* is not a *Delay* object.

abstract compute_output(*input*: Any, ***kwargs*) → Any

Computes output for given input.

This is where the user will implement the desired functionality of the *ProcessingBlock* without regard for latency.

Parameters

- **input** (*Any*) – Data to be processed. Passed in from `process()`.
- ****kwargs** (*Any*) – optional key-value argument pairs passed from `process()`. Could be used to pass in such values as the IO processor’s `walltime` or the measurement time for time- dependent functions.

Returns output

Return type *Any*

delay: `cleo.ioproc.delays.Delay`

The delay object determining compute latency for the block

process(*input: Any, t_in_ms: float, **kwargs*) → *Tuple[Any, float]*

Computes and saves output and output time given input and input time.

The user should implement `compute_output()` for their child classes, which performs the computation itself without regards for timing or saving variables.

Parameters

- **input** (*Any*) – Data to be processed
- **t_in_ms** (*float*) – Time the block receives the input data
- ****kwargs** (*Any*) – Key-value list of arguments passed to `compute_output()`

Returns output, out time in milliseconds

Return type *Tuple[Any, float]*

save_history: *bool*

Whether to record `t_in_ms`, `t_out_ms`, and `values` with every timestep

t_in_ms: *list[float]*

The walltime the block received each input. Only recorded if `save_history`

t_out_ms: *list[float]*

The walltime of each of the block’s outputs. Only recorded if `save_history`

values: *list[Any]*

Each of the block’s outputs. Only recorded if `save_history`

class `cleo.ioproc.RecordOnlyProcessor`(*sample_period_ms, **kwargs*)

Bases: `cleo.ioproc.base.LatencyIOProcessor`

Take samples without performing any control.

Use this if all you are doing is recording.

Parameters `sample_period_ms` (*float*) – Determines how frequently samples are taken from the network.

Keyword Arguments

- **sampling** (*str*) – “fixed” or “when idle”; “fixed” by default
 “fixed” sampling means samples are taken on a fixed schedule, with no exceptions.
 “when idle” sampling means no samples are taken before the previous sample’s output has been delivered. A sample is taken ASAP after an over-period computation: otherwise remains on schedule.

- **processing** (*str*) – “parallel” or “serial”; “parallel” by default

”parallel” computes the output time by adding the delay for a sample onto the sample time, so if the delay is 2 ms, for example, while the sample period is only 1 ms, some of the processing is happening in parallel. Output order matches input order even if the computed output time for a sample is sooner than that for a previous sample.

”serial” computes the output time by adding the delay for a sample onto the output time of the previous sample, rather than the sampling time. Note this may be of limited utility because it essentially means the *entire* round trip cannot be in parallel at all. More realistic is that simply each block or phase of computation must be serial. If anyone cares enough about this, it will have to be implemented in the future.

Note: It doesn’t make much sense to combine parallel computation with “when idle” sampling, because “when idle” sampling only produces one sample at a time to process.

Raises **ValueError** – For invalid *sampling* or *processing* kwargs

process(*state_dict: dict*, *sample_time_ms: float*) → Tuple[dict, float]

Process network state to generate output to update stimulators.

This is the function the user must implement to define the signal processing pipeline.

Parameters

- **state_dict** (*dict*) – {*recorder_name: state*} dictionary from `get_state()`
- **time_ms** (*float*) –

Returns {‘stim_name’: *ctrl_signal*} dictionary and output time in milliseconds.

Return type Tuple[dict, float]

sample_period_ms: float

Determines how frequently the processor takes samples

t_samp_ms: list[float]

Record of sampling times—each time `put_state()` is called.

6.3.8 cleo.viz module

Tools for visualizing models and simulations

```
class cleo.viz.VideoVisualizer(devices: collections.abc.Iterable[Union[cleo.base.InterfaceDevice,
    Tuple[cleo.base.InterfaceDevice, dict]]] = NOTHING, dt:
    brian2.units.fundamentalunits.Quantity = 1. * msecond, *, name: str =
    NOTHING, save_history: bool = True)
```

Bases: `cleo.base.InterfaceDevice`

Device for visualizing a simulation.

Must be injected after all other devices and before the simulation is run.

Method generated by attrs for class VideoVisualizer.

ax: `matplotlib.axes._axes.Axes`

connect_to_neuron_group(*neuron_group*: *brian2.groups.neurongroup.NeuronGroup*, ***kwargs*) → None

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a *NeuronGroup*, *Synapses*, or *Monitor*, make sure to add these to *self.brian_objects*.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwargs** (optional, passed from *inject* or) – *inject*

devices: *collections.abc.Iterable[Union[cleo.base.InterfaceDevice, Tuple[cleo.base.InterfaceDevice, dict]]]*

list of devices or (device, vis_kwargs) tuples to include in the plot, just as in the *plot()* function, by default “all”, which will include all recorders and stimulators currently injected when this visualizer is injected into the simulator.

dt: *brian2.units.fundamentalunits.Quantity*

length of each frame—that is, every *dt* the visualizer takes a snapshot of the network, by default 1*ms

fig: *matplotlib.figure.Figure*

generate_Animation(*plotargs*: *dict*, *slowdown_factor*: *float = 10*, ***figargs*: *Any*) → *matplotlib.animation.Animation*

Create a matplotlib Animation object from the recorded simulation

Parameters

- **plotargs** (*dict*) – dictionary of arguments as taken by *plot()*. can include *xlim*, *ylim*, *zlim*, *colors*, *axis_scale_unit*, *invert_z*, and/or *scatterargs*. neuron groups and devices are automatically added and ***figargs* are specified separately.
- **slowdown_factor** (*float*, *optional*) – how much slower the animation will be rendered, as a multiple of real-time, by default 10
- ****figargs** (*Any*, *optional*) – keyword arguments passed to *plt.figure()*, such as *figsize*

Returns An Animation object capturing the desired visualization. See matplotlib’s docs for saving and rendering options.

Return type *matplotlib.animation.Animation*

init_for_simulator(*simulator*: *cleo.base.CLSimulator*)

Initialize device for simulator on initial injection

This function is called only the first time a device is injected into a simulator and performs any operations that are independent of the individual neuron groups it is connected to.

Parameters **simulator** (*CLSimulator*) – simulator being injected into

neuron_groups: *list[brian2.groups.neurongroup.NeuronGroup]*

cleo.viz.plot(**neuron_groups*: *brian2.groups.neurongroup.NeuronGroup*, *xlim*: *Optional[Tuple[float, float]] = None*, *ylim*: *Optional[Tuple[float, float]] = None*, *zlim*: *Optional[Tuple[float, float]] = None*, *colors*: *Optional[collections.abc.Iterable] = None*, *axis_scale_unit*: *brian2.units.fundamentalunits.Unit = umetre*, *devices*: *collections.abc.Iterable[Union[cleo.base.InterfaceDevice, Tuple[cleo.base.InterfaceDevice, dict]]] = []*, *invert_z*: *bool = True*, *scatterargs*: *dict = {}*, *sim*: *Optional[cleo.base.CLSimulator] = None*, ***figargs*: *Any*) → None

Visualize neurons and interface devices

Parameters

- **xlim** (*Tuple[float, float], optional*) – xlim for plot, determined automatically by default
- **ylim** (*Tuple[float, float], optional*) – ylim for plot, determined automatically by default
- **zlim** (*Tuple[float, float], optional*) – zlim for plot, determined automatically by default
- **colors** (*Iterable, optional*) – colors, one for each neuron group, automatically determined by default
- **axis_scale_unit** (*Unit, optional*) – Brian unit to scale lim params, by default mm
- **devices** (*Iterable[Union[InterfaceDevice, Tuple[InterfaceDevice, dict]]], optional*) – devices to add to the plot or (device, kwargs) tuples. `add_self_to_plot()` is called for each, using the kwargs dict if given. By default []
- **invert_z** (*bool, optional*) – whether to invert z-axis, by default True to reflect the convention that +z represents depth from cortex surface
- **scatterargs** (*dict, optional*) – arguments passed to `plt.scatter()` for each neuron group, such as marker
- **sim** (*CLSimulator, optional*) – Optional shortcut to include all neuron groups and devices
- ****figargs** (*Any, optional*) – keyword arguments passed to `plt.figure()`, such as `figsize`

Raises **ValueError** – When neuron group doesn't have x, y, and z already defined

6.3.9 cleo.recorders module

Contains basic recorders.

class `cleo.recorders.GroundTruthSpikeRecorder`(**, name: str = NOTHING, save_history: bool = True*)

Bases: `cleo.base.Recorder`

Reports the number of spikes seen since last queried for each neuron.

This amounts effectively to the number of spikes per control period. Note: this will only work for one neuron group at the moment.

Method generated by attrs for class `GroundTruthSpikeRecorder`.

connect_to_neuron_group(*neuron_group*)

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a `NeuronGroup`, `Synapses`, or `Monitor`, make sure to add these to `self.brian_objects`.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwargs** (optional, passed from *inject* or) – *inject*

get_state() → `nptyping.types._ndarray.NDArray[Any, nptyping.types._number.UInt]`

Returns `n_neurons`-length array with spike counts over the latest control period.

Return type `NDArray[(n_neurons,), np.uint]`

neuron_group: `brian2.groups.neurongroup.NeuronGroup`

class `cleo.recorders.RateRecorder`(*i: int, *, name: str = NOTHING, save_history: bool = True*)

Bases: `cleo.base.Recorder`

Records firing rate from a single neuron.

Firing rate comes from Brian's `PopulationRateMonitor`

Method generated by attrs for class `RateRecorder`.

connect_to_neuron_group(*neuron_group*)

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a `NeuronGroup`, `Synapses`, or `Monitor`, make sure to add these to *self.brian_objects*.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwargs** (optional, passed from *inject* or) – *inject*

get_state()

Return current measurement.

i: int

index of neuron to record

mon: `brian2.monitors.ratemonitor.PopulationRateMonitor`

class `cleo.recorders.VoltageRecorder`(*voltage_var_name: str = 'v', *, name: str = NOTHING, save_history: bool = True*)

Bases: `cleo.base.Recorder`

Records the voltage of a single neuron group.

Method generated by attrs for class `VoltageRecorder`.

connect_to_neuron_group(*neuron_group*)

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a `NeuronGroup`, `Synapses`, or `Monitor`, make sure to add these to *self.brian_objects*.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwargs** (optional, passed from *inject* or) – *inject*

get_state() → `brian2.units.fundamentalunits.Quantity`

Returns Current voltage of target neuron group

Return type `Quantity`

mon: `brian2.monitors.statemonitor.StateMonitor`

voltage_var_name: str

Name of variable representing membrane voltage

6.3.10 cleo.stimulators module

Contains basic stimulators.

class cleo.stimulators.StateVariableSetter(*default_value: Any = 0, *, name: str = NOTHING, save_history: bool = True, variable_to_ctrl: str, unit: Unit*)

Bases: *cleo.base.Stimulator*

Sets the given state variable of target neuron groups.

Method generated by attrs for class StateVariableSetter.

connect_to_neuron_group(*neuron_group*)

Connect device to given *neuron_group*.

If your device introduces any objects which Brian must keep track of, such as a NeuronGroup, Synapses, or Monitor, make sure to add these to *self.brian_objects*.

Parameters

- **neuron_group** (*NeuronGroup*) –
- ****kwargs** (optional, passed from *inject* or) – *inject*

neuron_groups: list[*NeuronGroup*]

unit: Unit

will be used in *update()*

Type Unit of controlled variable

update(*ctrl_signal: float*) → None

Set state variable of target neuron groups

Parameters **ctrl_signal** (*float*) – Value to update variable to, without unit. The unit provided on initialization is automatically multiplied.

variable_to_ctrl: str

Name of state variable to control

6.3.11 cleo.utilities module

Assorted utilities for developers.

cleo.utilities.add_to_neo_segment(*segment: neo.core.segment.Segment, *objects: neo.core.dataobject.DataObject*)

Taken from *neo.core.group.Group*.

cleo.utilities.analog_signal(*t_ms, values_no_unit, units*) → *neo.core.basesignal.BaseSignal*

cleo.utilities.get_orth_vectors_for_V(*V*)

For nx3 block of row vectors *V*, return nx3 *W1*, *W2* orthogonal vector blocks

cleo.utilities.modify_model_with_eqs(*neuron_group, eqs_to_add*)

Adapted from *_create_variables()* from *neurongroup.py* from Brian2 source code v2.3.0.2

`cleo.utilities.normalize_coords`(*coords*: *brian2.units.fundamentalunits.Quantity*) → *brian2.units.fundamentalunits.Quantity*

Normalize coordinates to unit vectors.

`cleo.utilities.style_plots_for_docs`(*dark*=*True*)

`cleo.utilities.times_are_regular`(*times*)

`cleo.utilities.uniform_cylinder_rz`(*n*, *rmax*, *zmax*)

`cleo.utilities.wavelength_to_rgb`(*wavelength_nm*, *gamma*=0.8)

taken from http://www.noah.org/wiki/Wavelength_to_RGB_in_Python This converts a given wavelength of light to an approximate RGB color value. The wavelength must be given in nanometers in the range from 380 nm through 750 nm (789 THz through 400 THz).

Based on code by Dan Bruton <http://www.physics.sfasu.edu/astro/color/spectra.html>

`cleo.utilities.xyz_from_rz`(*rs*, *thetas*, *zs*, *xyz_start*, *xyz_end*)

Convert from cylindrical to Cartesian coordinates.

6.3.12 cleo.registry module

class `cleo.registry.DeviceInteractionRegistry`(*sim*: *CLSimulator*)

Bases: object

Facilitates the creation and maintenance of ‘neurons’ and ‘synapses’ implementing many-to-many light-opsin/indicator optogenetics

Method generated by attrs for class DeviceInteractionRegistry.

connect_light_to_ldd_for_ng(*light*: *Light*, *ldd*: *LightDependent*, *ng*: *NeuronGroup*)

connections: set[Tuple['Light', 'LightDependent', NeuronGroup]]

init_register_light(*light*: *Light*)

ldds_for_ng: dict[NeuronGroup, set['LightDependent']]

light_prop_model = '\n T : 1\n epsilon : 1\n E_photon : joule\n Irr_post = epsilon * T * Irr0_pre : watt/meter**2 (summed)\n phi_post = Irr_post / E_photon : 1/second/meter**2 (summed)\n '

light_prop_syns: dict[Tuple['LightDependent', NeuronGroup], Synapses]

light_source_ng: NeuronGroup

Represents ALL light sources (multiple devices)

lights_for_ng: dict[NeuronGroup, set['Light']]

register(*device*: *InterfaceDevice*, *ng*: *NeuronGroup*)

Registers a device with the registry

register_ldd(*ldd*: *LightDependent*, *ng*: *NeuronGroup*)

Connects lights previously injected into this neuron group to this opsin

register_light(*light*: *Light*, *ng*: *NeuronGroup*)

Connects light to opsins and indicators already injected into this neuron group

sim: *CLSimulator*

source_for_light(*light*: *Light*) → Subgroup

subgroup_idx_for_light: dict['Light', slice]

cleo.registry.registry_for_sim(*sim*: *CLSimulator*) → *DeviceInteractionRegistry*

Returns the registry for the given simulator

INDICES AND TABLES

- `genindex`
- `modindex`
- `search`

PYTHON MODULE INDEX

C

- `cleo`, [74](#)
- `cleo.coords`, [80](#)
- `cleo.ephys`, [82](#)
- `cleo.imaging`, [98](#)
- `cleo.ioproc`, [105](#)
- `cleo.light`, [89](#)
- `cleo.opto`, [93](#)
- `cleo.recorders`, [112](#)
- `cleo.registry`, [115](#)
- `cleo.stimulators`, [114](#)
- `cleo.utilities`, [114](#)
- `cleo.viz`, [110](#)

A

A (*cleo.imaging.DoubExpCalBindingActivation* attribute), 98
 action_spectrum (*cleo.opto.Opsin* property), 96
 add_self_to_plot() (*cleo.ephys.Probe* method), 82
 add_self_to_plot() (*cleo.imaging.Scope* method), 101
 add_self_to_plot() (*cleo.InterfaceDevice* method), 76
 add_self_to_plot() (*cleo.light.Light* method), 90
 add_signals() (*cleo.ephys.Probe* method), 83
 add_to_neo_segment() (in module *cleo.utilities*), 114
 analog_signal() (in module *cleo.utilities*), 114
 assign_coords() (in module *cleo.coords*), 80
 assign_coords_grid_rect_prism() (in module *cleo.coords*), 80
 assign_coords_rand_cylinder() (in module *cleo.coords*), 80
 assign_coords_rand_rect_prism() (in module *cleo.coords*), 80
 assign_coords_uniform_cylinder() (in module *cleo.coords*), 81
 assign_xyz() (in module *cleo.coords*), 81
 ax (*cleo.viz.VideoVisualizer* attribute), 110

B

B_T (*cleo.imaging.DynamicCalcium* attribute), 99
 BansalFourStateOpsin (class in *cleo.opto*), 93
 bind_act_model (*cleo.imaging.GECI* attribute), 100
 brian_objects (*cleo.ephys.Signal* attribute), 84
 brian_objects (*cleo.InterfaceDevice* attribute), 76

C

Ca_rest (*cleo.imaging.DoubExpCalBindingActivation* attribute), 98
 Ca_rest (*cleo.imaging.DynamicCalcium* attribute), 99
 cal_model (*cleo.imaging.GECI* attribute), 100
 CalBindingActivationModel (class in *cleo.imaging*), 98
 CalciumModel (class in *cleo.imaging*), 98
 chr2_4s() (in module *cleo.opto*), 97
 chr2_b4s() (in module *cleo.opto*), 97

chrimson_4s() (in module *cleo.opto*), 97
 cleo
 module, 74
 cleo.coords
 module, 80
 cleo.ephys
 module, 82
 cleo.imaging
 module, 98
 cleo.ioproc
 module, 105
 cleo.light
 module, 89
 cleo.opto
 module, 93
 cleo.recorders
 module, 112
 cleo.registry
 module, 115
 cleo.stimulators
 module, 114
 cleo.utilities
 module, 114
 cleo.viz
 module, 110
 CLSimulator (class in *cleo*), 74
 compute() (*cleo.ioproc.ConstantDelay* method), 105
 compute() (*cleo.ioproc.Delay* method), 105
 compute() (*cleo.ioproc.GaussianDelay* method), 106
 compute_output() (*cleo.ioproc.FiringRateEstimator* method), 106
 compute_output() (*cleo.ioproc.PIDController* method), 108
 compute_output() (*cleo.ioproc.ProcessingBlock* method), 108
 concat_coords() (in module *cleo.coords*), 81
 connect_light_to_ldd_for_ng()
 (*cleo.registry.DeviceInteractionRegistry* method), 115
 connect_to_neuron_group()
 (*cleo.ephys.MultiUnitSpiking* method), 82
 connect_to_neuron_group() (*cleo.ephys.Probe*

- method), 83
- connect_to_neuron_group() (cleo.ephys.Signal method), 84
- connect_to_neuron_group() (cleo.ephys.SortedSpiking method), 85
- connect_to_neuron_group() (cleo.ephys.Spiking method), 85
- connect_to_neuron_group() (cleo.ephys.TKLFPSignal method), 86
- connect_to_neuron_group() (cleo.imaging.Scope method), 101
- connect_to_neuron_group() (cleo.InterfaceDevice method), 76
- connect_to_neuron_group() (cleo.light.Light method), 90
- connect_to_neuron_group() (cleo.recorders.GroundTruthSpikeRecorder method), 112
- connect_to_neuron_group() (cleo.recorders.RateRecorder method), 113
- connect_to_neuron_group() (cleo.recorders.VoltageRecorder method), 113
- connect_to_neuron_group() (cleo.stimulators.StateVariableSetter method), 114
- connect_to_neuron_group() (cleo.SynapseDevice method), 78
- connect_to_neuron_group() (cleo.viz.VideoVisualizer method), 110
- connections (cleo.registry.DeviceInteractionRegistry attribute), 115
- ConstantDelay (class in cleo.ioproc), 105
- coords (cleo.ephys.Probe attribute), 83
- coords (cleo.light.Light attribute), 90
- coords_from_ng() (in module cleo.coords), 81
- coords_from_xyz() (in module cleo.coords), 81
- cubic_interpolator() (in module cleo.light), 93
- cutoff_probability (cleo.ephys.Spiking attribute), 85
- ## D
- dCa_T (cleo.imaging.DynamicCalcium attribute), 99
- default_value (cleo.light.Light attribute), 90
- default_value (cleo.Stimulator attribute), 77
- Delay (class in cleo.ioproc), 105
- delay (cleo.ioproc.FiringRateEstimator attribute), 106
- delay (cleo.ioproc.ProcessingBlock attribute), 109
- DeviceInteractionRegistry (class in cleo.registry), 115
- devices (cleo.CLSimulator attribute), 74
- devices (cleo.viz.VideoVisualizer attribute), 111
- dFF (cleo.imaging.Scope attribute), 101
- dFF_1AP (cleo.imaging.Scope property), 102
- dFF_1AP (cleo.imaging.Sensor attribute), 103
- dFF_max (cleo.imaging.GECI attribute), 100
- direction (cleo.imaging.Scope attribute), 102
- direction (cleo.light.Light attribute), 91
- DoubExpCalBindingActivation (class in cleo.imaging), 98
- dt (cleo.viz.VideoVisualizer attribute), 111
- DynamicCalcium (class in cleo.imaging), 99
- ## E
- E (cleo.opto.BansalFourStateOpsin attribute), 94
- E (cleo.opto.FourStateOpsin attribute), 95
- epsilon() (cleo.light.LightDependent method), 92
- equal_photon_flux_spectrum() (in module cleo.light), 93
- exc_model (cleo.imaging.GECI attribute), 100
- exc_spectrum (cleo.imaging.Sensor property), 103
- ExcitationModel (class in cleo.imaging), 99
- extra_namespace (cleo.opto.FourStateOpsin attribute), 95
- extra_namespace (cleo.SynapseDevice attribute), 78
- ## F
- fiber473nm() (in module cleo.light), 93
- fig (cleo.viz.VideoVisualizer attribute), 111
- FiringRateEstimator (class in cleo.ioproc), 105
- fluor_model (cleo.imaging.GECI attribute), 100
- focus_coords_per_injct (cleo.imaging.Scope attribute), 102
- focus_depth (cleo.imaging.Scope attribute), 102
- FourStateOpsin (class in cleo.opto), 95
- ## G
- g0 (cleo.opto.BansalFourStateOpsin attribute), 94
- g0 (cleo.opto.FourStateOpsin attribute), 96
- gamma (cleo.imaging.DynamicCalcium attribute), 99
- gamma (cleo.opto.BansalFourStateOpsin attribute), 94
- gamma (cleo.opto.FourStateOpsin attribute), 96
- GaussianDelay (class in cleo.ioproc), 106
- GaussianEllipsoid (class in cleo.light), 89
- Gb0 (cleo.opto.BansalFourStateOpsin attribute), 94
- Gb0 (cleo.opto.FourStateOpsin attribute), 95
- gcamp3() (in module cleo.imaging), 103
- gcamp6f() (in module cleo.imaging), 103
- gcamp6rs06() (in module cleo.imaging), 103
- gcamp6rs09() (in module cleo.imaging), 104
- gcamp6s() (in module cleo.imaging), 104
- Gd1 (cleo.opto.BansalFourStateOpsin attribute), 94
- Gd1 (cleo.opto.FourStateOpsin attribute), 95
- Gd2 (cleo.opto.BansalFourStateOpsin attribute), 94
- Gd2 (cleo.opto.FourStateOpsin attribute), 95
- GECI (class in cleo.imaging), 99
- geci() (in module cleo.imaging), 104
- generate_Animation() (cleo.viz.VideoVisualizer method), 111

- [get_ctrl_signal\(\)](#) (*cleo.ioproc.LatencyIOProcessor method*), 107
[get_ctrl_signal\(\)](#) (*cleo.IOProcessor method*), 75
[get_orth_vectors_for_V\(\)](#) (in module *cleo.utilities*), 114
[get_state\(\)](#) (*cleo.CLSimulator method*), 74
[get_state\(\)](#) (*cleo.ephys.MultiUnitSpiking method*), 82
[get_state\(\)](#) (*cleo.ephys.Probe method*), 83
[get_state\(\)](#) (*cleo.ephys.Signal method*), 84
[get_state\(\)](#) (*cleo.ephys.SortedSpiking method*), 85
[get_state\(\)](#) (*cleo.ephys.Spiking method*), 85
[get_state\(\)](#) (*cleo.ephys.TKLFPSignal method*), 86
[get_state\(\)](#) (*cleo.imaging.GECI method*), 100
[get_state\(\)](#) (*cleo.imaging.Scope method*), 102
[get_state\(\)](#) (*cleo.imaging.Sensor method*), 103
[get_state\(\)](#) (*cleo.Recorder method*), 77
[get_state\(\)](#) (*cleo.recorders.GroundTruthSpikeRecorder method*), 112
[get_state\(\)](#) (*cleo.recorders.RateRecorder method*), 113
[get_state\(\)](#) (*cleo.recorders.VoltageRecorder method*), 113
[Gf0](#) (*cleo.opto.BansalFourStateOpsin attribute*), 94
[Gf0](#) (*cleo.opto.FourStateOpsin attribute*), 95
[Gr0](#) (*cleo.opto.BansalFourStateOpsin attribute*), 94
[Gr0](#) (*cleo.opto.FourStateOpsin attribute*), 95
[GroundTruthSpikeRecorder](#) (class in *cleo.recorders*), 112
[gtacr2_4s\(\)](#) (in module *cleo.opto*), 97
- I**
- [i](#) (*cleo.ephys.Spiking attribute*), 85
[i](#) (*cleo.recorders.RateRecorder attribute*), 113
[I_per_Irr](#) (*cleo.opto.ProportionalCurrentOpsin attribute*), 97
[i_probe_by_i_ng](#) (*cleo.ephys.Spiking attribute*), 85
[i_targets_for_neuron_group\(\)](#) (*cleo.imaging.Scope method*), 102
[i_targets_per_injct](#) (*cleo.imaging.Scope attribute*), 102
[img_width](#) (*cleo.imaging.Scope attribute*), 102
[init_for_probe\(\)](#) (*cleo.ephys.Signal method*), 84
[init_for_probe\(\)](#) (*cleo.ephys.TKLFPSignal method*), 86
[init_for_simulator\(\)](#) (*cleo.InterfaceDevice method*), 76
[init_for_simulator\(\)](#) (*cleo.light.Light method*), 91
[init_for_simulator\(\)](#) (*cleo.viz.VideoVisualizer method*), 111
[init_register_light\(\)](#) (*cleo.registry.DeviceInteractionRegistry method*), 115
[init_syn_vars\(\)](#) (*cleo.imaging.DoubExpCalBindingActivation method*), 98
[init_syn_vars\(\)](#) (*cleo.imaging.DynamicCalcium method*), 99
[init_syn_vars\(\)](#) (*cleo.imaging.GECI method*), 100
[init_syn_vars\(\)](#) (*cleo.opto.BansalFourStateOpsin method*), 94
[init_syn_vars\(\)](#) (*cleo.opto.FourStateOpsin method*), 96
[init_syn_vars\(\)](#) (*cleo.SynapseDevice method*), 78
[inject\(\)](#) (*cleo.CLSimulator method*), 74
[inject_sensor_for_targets\(\)](#) (*cleo.imaging.Scope method*), 102
[InterfaceDevice](#) (class in *cleo*), 76
[io_processor](#) (*cleo.CLSimulator attribute*), 74
[IOProcessor](#) (class in *cleo*), 75
[is_sampling_now\(\)](#) (*cleo.ioproc.LatencyIOProcessor method*), 107
[is_sampling_now\(\)](#) (*cleo.IOProcessor method*), 75
- J**
- [jgcamp7b\(\)](#) (in module *cleo.imaging*), 104
[jgcamp7c\(\)](#) (in module *cleo.imaging*), 104
[jgcamp7f\(\)](#) (in module *cleo.imaging*), 104
[jgcamp7s\(\)](#) (in module *cleo.imaging*), 105
- K**
- [K](#) (*cleo.light.OpticFiber attribute*), 92
[k1](#) (*cleo.opto.BansalFourStateOpsin attribute*), 94
[k1](#) (*cleo.opto.FourStateOpsin attribute*), 96
[k2](#) (*cleo.opto.BansalFourStateOpsin attribute*), 94
[k2](#) (*cleo.opto.FourStateOpsin attribute*), 96
[K_d](#) (*cleo.imaging.GECI attribute*), 100
[kappa_S](#) (*cleo.imaging.DynamicCalcium attribute*), 99
[kb](#) (*cleo.opto.BansalFourStateOpsin attribute*), 94
[kb](#) (*cleo.opto.FourStateOpsin attribute*), 96
[kf](#) (*cleo.opto.BansalFourStateOpsin attribute*), 94
[kf](#) (*cleo.opto.FourStateOpsin attribute*), 96
- L**
- [LatencyIOProcessor](#) (class in *cleo.ioproc*), 106
[ldds_for_ng](#) (*cleo.registry.DeviceInteractionRegistry attribute*), 115
[lfp_uV](#) (*cleo.ephys.TKLFPSignal attribute*), 86
[Light](#) (class in *cleo.light*), 89
[light_agg_ngs](#) (*cleo.light.LightDependent property*), 92
[light_model](#) (*cleo.light.Light attribute*), 91
[light_prop_model](#) (*cleo.registry.DeviceInteractionRegistry attribute*), 115
[light_prop_syns](#) (*cleo.registry.DeviceInteractionRegistry attribute*), 115
[light_source_ng](#) (*cleo.registry.DeviceInteractionRegistry attribute*), 115
[LightDependent](#) (class in *cleo.light*), 92
[LightExcitation](#) (class in *cleo.imaging*), 100

LightModel (class in *cleo.light*), 92
 lights_for_ng (*cleo.registry.DeviceInteractionRegistry* attribute), 115
 linear_interpolator() (in module *cleo.light*), 93
 linear_shank_coords() (in module *cleo.ephys*), 87
 location (*cleo.imaging.GECI* attribute), 100
 location (*cleo.imaging.Scope* attribute), 102
 location (*cleo.imaging.Sensor* attribute), 103

M

max_Irr0_mW_per_mm2 (*cleo.light.Light* attribute), 91
 max_Irr0_mW_per_mm2_viz (*cleo.light.Light* attribute), 91
 model (*cleo.imaging.CalBindingActivationModel* attribute), 98
 model (*cleo.imaging.CalciumModel* attribute), 98
 model (*cleo.imaging.DoubExpCalBindingActivation* attribute), 98
 model (*cleo.imaging.DynamicCalcium* attribute), 99
 model (*cleo.imaging.ExcitationModel* attribute), 99
 model (*cleo.imaging.LightExcitation* attribute), 100
 model (*cleo.imaging.NullBindingActivation* attribute), 100
 model (*cleo.imaging.NullExcitation* attribute), 101
 model (*cleo.imaging.PreexistingCalcium* attribute), 101
 model (*cleo.opto.BansalFourStateOpsin* attribute), 95
 model (*cleo.opto.FourStateOpsin* attribute), 96
 model (*cleo.opto.ProportionalCurrentOpsin* attribute), 97
 model (*cleo.SynapseDevice* attribute), 78
 modify_model_and_params_for_ng() (*cleo.SynapseDevice* method), 79
 modify_model_with_eqs() (in module *cleo.utilities*), 114
 module
 cleo, 74
 cleo.coords, 80
 cleo.ephys, 82
 cleo.imaging, 98
 cleo.ioproc, 105
 cleo.light, 89
 cleo.opto, 93
 cleo.recorders, 112
 cleo.registry, 115
 cleo.stimulators, 114
 cleo.utilities, 114
 cleo.viz, 110
 mon (*cleo.recorders.RateRecorder* attribute), 113
 mon (*cleo.recorders.VoltageRecorder* attribute), 113
 MultiUnitSpiking (class in *cleo.ephys*), 82

N

n (*cleo.ephys.Probe* property), 83
 n (*cleo.imaging.Scope* property), 102

n (*cleo.light.Light* property), 91
 n_H (*cleo.imaging.GECI* attribute), 100
 NAFib (*cleo.light.OpticFiber* attribute), 93
 name (*cleo.ephys.Signal* attribute), 84
 name (*cleo.InterfaceDevice* attribute), 77
 network (*cleo.CLSimulator* attribute), 74
 neuron_group (*cleo.recorders.GroundTruthSpikeRecorder* attribute), 113
 neuron_groups (*cleo.imaging.Scope* attribute), 102
 neuron_groups (*cleo.stimulators.StateVariableSetter* attribute), 114
 neuron_groups (*cleo.viz.VideoVisualizer* attribute), 111
 normalize_coords() (in module *cleo.utilities*), 114
 ntis (*cleo.light.OpticFiber* attribute), 93
 NullBindingActivation (class in *cleo.imaging*), 100
 NullExcitation (class in *cleo.imaging*), 100

O

ogb1() (in module *cleo.imaging*), 105
 on_pre (*cleo.imaging.CalciumModel* attribute), 98
 on_pre (*cleo.imaging.DynamicCalcium* attribute), 99
 on_pre (*cleo.SynapseDevice* attribute), 79
 Opsin (class in *cleo.opto*), 96
 OpticFiber (class in *cleo.light*), 92

P

p (*cleo.opto.BansalFourStateOpsin* attribute), 95
 p (*cleo.opto.FourStateOpsin* attribute), 96
 params (*cleo.imaging.GECI* property), 100
 params (*cleo.SynapseDevice* property), 79
 per_ng_unit_replacements (*cleo.SynapseDevice* attribute), 79
 phim (*cleo.opto.BansalFourStateOpsin* attribute), 95
 phim (*cleo.opto.FourStateOpsin* attribute), 96
 PIController (class in *cleo.ioproc*), 108
 plot() (in module *cleo.viz*), 111
 plot_spectra() (in module *cleo.light*), 93
 poly2_shank_coords() (in module *cleo.ephys*), 87
 poly3_shank_coords() (in module *cleo.ephys*), 87
 PreexistingCalcium (class in *cleo.imaging*), 101
 Probe (class in *cleo.ephys*), 82
 probe (*cleo.ephys.Probe* attribute), 83
 probe (*cleo.ephys.Signal* attribute), 84
 process() (*cleo.ioproc.LatencyIOProcessor* method), 107
 process() (*cleo.ioproc.ProcessingBlock* method), 109
 process() (*cleo.ioproc.RecordOnlyProcessor* method), 110
 ProcessingBlock (class in *cleo.ioproc*), 108
 ProportionalCurrentOpsin (class in *cleo.opto*), 97
 put_state() (*cleo.ioproc.LatencyIOProcessor* method), 107
 put_state() (*cleo.IOProcessor* method), 76

Q

q (*cleo.opto.BansalFourStateOpsin* attribute), 95

q (*cleo.opto.FourStateOpsin* attribute), 96

R

R0 (*cleo.light.OpticFiber* attribute), 93

r_half_detection (*cleo.ephys.Spiking* attribute), 85

r_perfect_detection (*cleo.ephys.Spiking* attribute), 85

rand_seed (*cleo.imaging.Scope* attribute), 102

RateRecorder (class in *cleo.recorders*), 113

Recorder (class in *cleo*), 77

recorders (*cleo.CLSimulator* attribute), 74

RecordOnlyProcessor (class in *cleo.ioproc*), 109

ref_signal (*cleo.ioproc.PIController* attribute), 108

register() (*cleo.registry.DeviceInteractionRegistry* method), 115

register_ldd() (*cleo.registry.DeviceInteractionRegistry* method), 115

register_light() (*cleo.registry.DeviceInteractionRegistry* method), 115

registry_for_sim() (in module *cleo.registry*), 116

required_vars (*cleo.opto.ProportionalCurrentOpsin* attribute), 97

required_vars (*cleo.SynapseDevice* attribute), 79

reset() (*cleo.CLSimulator* method), 74

reset() (*cleo.ephys.Probe* method), 83

reset() (*cleo.ephys.Signal* method), 84

reset() (*cleo.ephys.Spiking* method), 86

reset() (*cleo.ephys.TKLFPSignal* method), 86

reset() (*cleo.imaging.Scope* method), 102

reset() (*cleo.InterfaceDevice* method), 77

reset() (*cleo.IOProcessor* method), 76

reset() (*cleo.Stimulator* method), 77

reset() (*cleo.SynapseDevice* method), 79

rho_rel_per_injct (*cleo.imaging.Scope* attribute), 102

run() (*cleo.CLSimulator* method), 75

S

S (*cleo.light.OpticFiber* attribute), 93

sample_period_ms (*cleo.ioproc.RecordOnlyProcessor* attribute), 110

sample_period_ms (*cleo.IOProcessor* attribute), 76

save_history (*cleo.InterfaceDevice* attribute), 77

save_history (*cleo.ioproc.FiringRateEstimator* attribute), 106

save_history (*cleo.ioproc.ProcessingBlock* attribute), 109

Scope (class in *cleo.imaging*), 101

Sensor (class in *cleo.imaging*), 102

sensor (*cleo.imaging.Scope* attribute), 102

set_io_processor() (*cleo.CLSimulator* method), 75

sigma_axial (*cleo.light.GaussianEllipsoid* attribute), 89

sigma_lateral (*cleo.light.GaussianEllipsoid* attribute), 89

sigma_noise (*cleo.imaging.Scope* property), 102

sigma_noise (*cleo.imaging.Sensor* attribute), 103

sigma_per_injct (*cleo.imaging.Scope* attribute), 102

Signal (class in *cleo.ephys*), 84

signals (*cleo.ephys.Probe* attribute), 83

sim (*cleo.InterfaceDevice* attribute), 77

sim (*cleo.registry.DeviceInteractionRegistry* attribute), 115

snr (*cleo.imaging.Sensor* property), 103

snr_cutoff (*cleo.imaging.Scope* attribute), 102

soma_radius (*cleo.imaging.Scope* attribute), 102

SortedSpiking (class in *cleo.ephys*), 84

source (*cleo.light.Light* property), 91

source_for_light() (*cleo.registry.DeviceInteractionRegistry* method), 116

source_ngs (*cleo.SynapseDevice* attribute), 79

spectrum (*cleo.light.LightDependent* attribute), 92

spectrum (*cleo.opto.Opsin* attribute), 96

spectrum_interpolator (*cleo.light.LightDependent* attribute), 92

spectrum_interpolator (*cleo.opto.Opsin* attribute), 96

Spiking (class in *cleo.ephys*), 85

StateVariableSetter (class in *cleo.stimulators*), 114

Stimulator (class in *cleo*), 77

stimulators (*cleo.CLSimulator* attribute), 75

style_plots_for_docs() (in module *cleo.utilities*), 115

subgroup_idx_for_light
(*cleo.registry.DeviceInteractionRegistry* attribute), 116

SynapseDevice (class in *cleo*), 78

synapses (*cleo.SynapseDevice* attribute), 79

T

t_in_ms (*cleo.ioproc.FiringRateEstimator* attribute), 106

t_in_ms (*cleo.ioproc.ProcessingBlock* attribute), 109

t_ms (*cleo.ephys.Spiking* attribute), 86

t_ms (*cleo.ephys.TKLFPSignal* attribute), 86

t_ms (*cleo.imaging.Scope* attribute), 102

t_ms (*cleo.Stimulator* attribute), 77

t_out_ms (*cleo.ioproc.FiringRateEstimator* attribute), 106

t_out_ms (*cleo.ioproc.ProcessingBlock* attribute), 109

t_samp_ms (*cleo.ephys.Spiking* attribute), 86

t_samp_ms (*cleo.ioproc.LatencyIOProcessor* attribute), 108

t_samp_ms (*cleo.ioproc.RecordOnlyProcessor* attribute), 110

- target_neurons_in_plane() (*cleo.imaging.Scope* method), 102
- target_neurons_in_plane() (in module *cleo.imaging*), 105
- tau_off (*cleo.imaging.DoubExpCalBindingActivation* attribute), 98
- tau_on (*cleo.imaging.DoubExpCalBindingActivation* attribute), 99
- tetrode_shank_coords() (in module *cleo.ephys*), 88
- tile_coords() (in module *cleo.ephys*), 88
- times_are_regular() (in module *cleo.utilities*), 115
- TKLFPSignal (class in *cleo.ephys*), 86
- to_neo() (*cleo.CLSimulator* method), 75
- to_neo() (*cleo.ephys.MultiUnitSpiking* method), 82
- to_neo() (*cleo.ephys.Probe* method), 83
- to_neo() (*cleo.ephys.Spiking* method), 86
- to_neo() (*cleo.ephys.TKLFPSignal* method), 86
- to_neo() (*cleo.light.Light* method), 91
- to_neo() (*cleo.Stimulator* method), 77
- tp_light_from_scope() (in module *cleo.light*), 93
- transmittance() (*cleo.light.GaussianEllipsoid* method), 89
- transmittance() (*cleo.light.Light* method), 91
- transmittance() (*cleo.light.LightModel* method), 92
- transmittance() (*cleo.light.OpticFiber* method), 93
- ## U
- uLFP_threshold_uV (*cleo.ephys.TKLFPSignal* attribute), 87
- uniform_cylinder_rz() (in module *cleo.utilities*), 115
- unit (*cleo.stimulators.StateVariableSetter* attribute), 114
- update() (*cleo.light.Light* method), 91
- update() (*cleo.Stimulator* method), 78
- update() (*cleo.stimulators.StateVariableSetter* method), 114
- update_artists() (*cleo.InterfaceDevice* method), 77
- update_artists() (*cleo.light.Light* method), 91
- update_stimulators() (*cleo.CLSimulator* method), 75
- ## V
- v0 (*cleo.opto.FourStateOpsin* attribute), 96
- v1 (*cleo.opto.FourStateOpsin* attribute), 96
- value (*cleo.Stimulator* attribute), 78
- values (*cleo.ioproc.FiringRateEstimator* attribute), 106
- values (*cleo.ioproc.ProcessingBlock* attribute), 109
- values (*cleo.Stimulator* attribute), 78
- variable_to_ctrl (*cleo.stimulators.StateVariableSetter* attribute), 114
- vfchromson_4s() (in module *cleo.opto*), 98
- VideoVisualizer (class in *cleo.viz*), 110
- viz_params() (*cleo.light.GaussianEllipsoid* method), 89
- viz_params() (*cleo.light.LightModel* method), 92
- viz_params() (*cleo.light.OpticFiber* method), 93
- voltage_var_name (*cleo.recorders.VoltageRecorder* attribute), 113
- VoltageRecorder (class in *cleo.recorders*), 113
- ## W
- wavelength (*cleo.light.Light* attribute), 91
- wavelength_to_rgb() (in module *cleo.utilities*), 115
- ## X
- xs (*cleo.ephys.Probe* property), 83
- xyz_from_rz() (in module *cleo.utilities*), 115
- ## Y
- ys (*cleo.ephys.Probe* property), 83
- ## Z
- zs (*cleo.ephys.Probe* property), 84